

Published paper


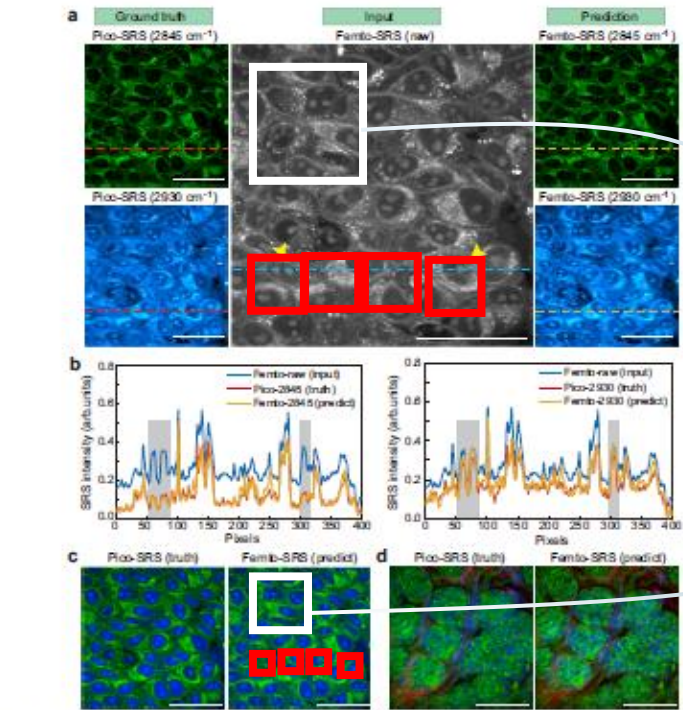


Fig. 2 Deep U-Net based femto-SRS imaging with restored chemical resolution. **a**, A single raw femto-SRS image (middle) was converted by U-Net into the spectrally resolved images at 2845 cm^{-1} and 2930 cm^{-1} (right), compared with the ground truth pico-SRS images (left). **b**, Intensity profiles (dashed lines in **a**) of the predicted and ground-truth data, showing chemical contrast in the cell nucleus regions (shadows, and arrows in **a**). **c**, **d**, Femto-SRS imaging of live cells and gastric tissues, yielding almost identical results to pico-SRS. Green: lipids; blue: protein; red: collagen fibers. Scale bars: 50 μm . Source data are provided in the Source data file.

or nest-like distribution without forming tubular structures (Fig. 4d). In addition, SRS was even able to differentiate signaling cell components similar to that of H&E staining, showing abundant intracellular mucus pushing the nucleus to one side (Fig. 4e).

Despite the overall high agreement between SRS and H&E, a few discrepancies remain and are noteworthy. First of all, SRS maps out the distributions of lipids and protein, whereas H&E mainly stains nucleic acid and protein. Secondly, we demonstrated SRS images on fresh tissues in comparison with H&E staining on thin sections. Although SRS images create thin optical sections of $\sim 1 \mu\text{m}$ thickness due to the nonlinear property of optical signals³⁴, they showed slightly different structural features from that of H&E. For instance, the core of the glandular ducts usually appeared hollow in H&E but filled in SRS (Fig. 3a, b and 4a), which is likely due to the loss and dissolution of contents during the sectioning and deparaffinization processes of H&E. And the dehydration process might result in denser cellular patterns in H&E. Therefore, SRS might reveal more intact and better preserved histoarchitectures of fresh gastric tissues (Fig. S3). Furthermore, label-free SRS is free from the large staining variation of H&E as seen in the different colorings among tissues (Figs. 3 and 4). The more consistent appearance of SRS images may help reduce the uncertainty and improve the accuracy of machine learning based diagnosis.

Extract single training example

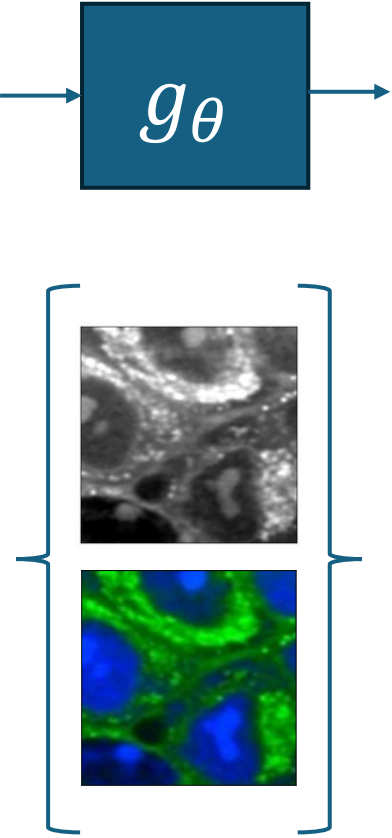

Testing

Input

Original model

Our model

