

(a) Interactive Medical-SAM2 - Initial Setup

Interactive Medical-SAM2 GUI Setup

Mode

☐ Automatic ☒ Manual

Method (auto only)

☐ Detection ☒ Segmentation


Options

☐ Perform preprocessing

☐ Show initial mask (auto segmentation)

Data Path

C:/Users/user/test_data

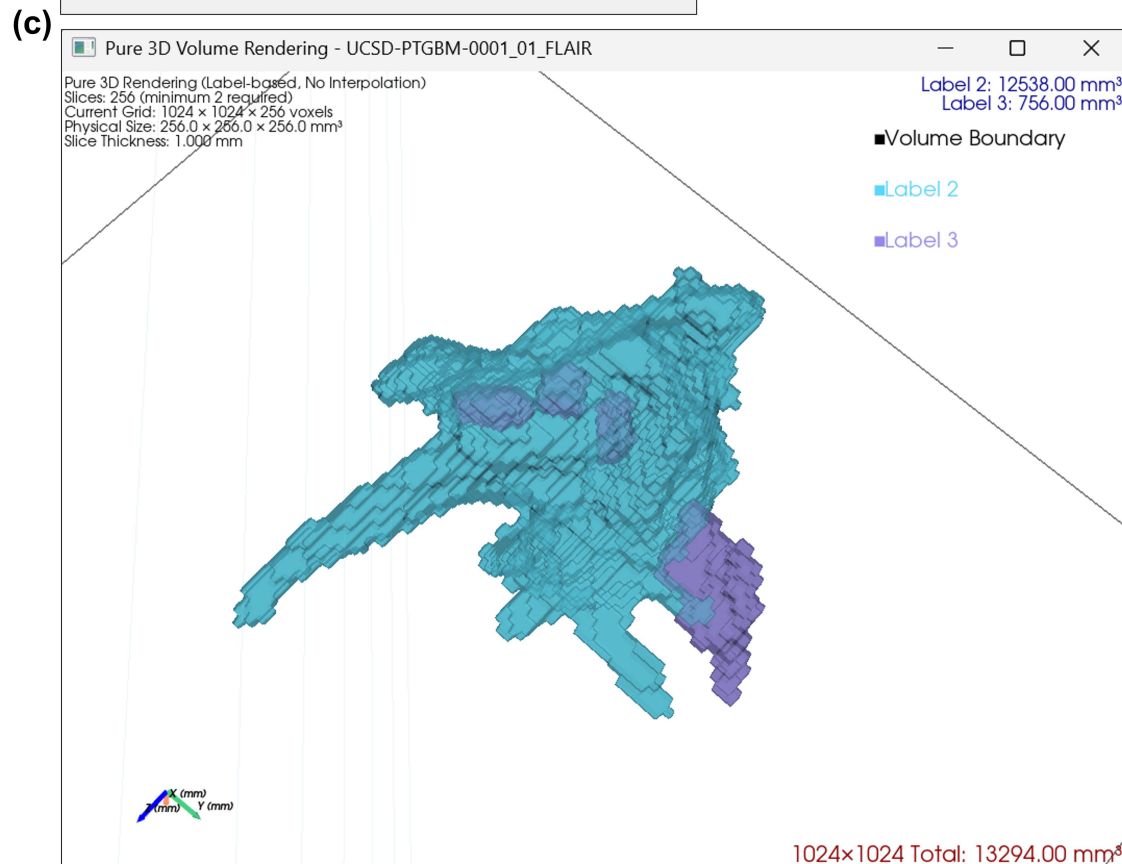

(b) Patient 0: UCSD-PTGBM-0001_01_FLAIR.nii.gz - Selec...

Patient 1: UCSD-PTGBM-0001_01_FLAIR.nii.gz

Select patient:
 UCSD-PTGBM-0001_01_FLAIR.nii.gz

Preprocessing Setup:

☐ Perform preprocessing (N4 bias correction + intensity normalization)
 Correct MRI inhomogeneity with N4 bias field correction
 Apply intensity clipping and normalization
 Preprocessed files are saved in preprocessed folder

Mode Selection:
 manual

Method Selection:
 Segmentation

Detection Model:
 sam2_det
 Enter custom detection model name

Segmentation Model:
 sam2_seg

nnUNetv2 model folder: /path/to/nnunetv2_model

Double Viewer Setup:

☐ Use Double Viewer

Path: Select file/folder path for double viewer

Mode Description:
 - auto: Perform automatic segmentation
 - manual: Manual prompt input

Method Description (for auto mode):
 - Detection: Detection
 - Segmentation: Segmentation with prompts

Preprocessing:
 - N4 bias field correction to correct MRI signal inhomogeneity
 - Intensity clipping (0.5%~99.5%) and Z-score normalization
 - Preprocessing results are automatically saved in preprocessed folder

Double Viewer:
 - Display another image simultaneously in an additional viewer
 - Can select DICOM folder or NIFTI file