


1 Data Source




SAR-RARP50

50 public RARP suturing videos
real robot-assisted
clinical videos

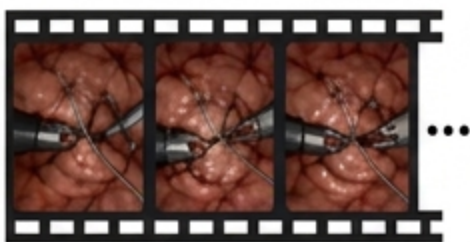

2 Motion Segmentation

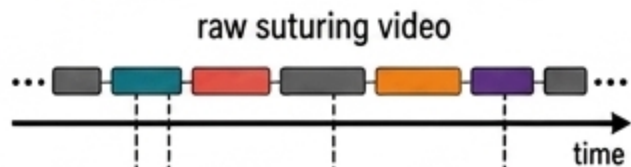

KEPT (valid motions)

effective
insertion



effective
extraction



FILTERED OUT (irrelevant / invalid)


failed
attempt



position
adjustment



needle pose
adjustment

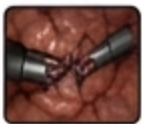

non-suturing
action



Keep only effective needle insertion /
extraction motions

3 20-frame Uniform Sampling

from one valid motion segment



I_1 I_2 I_3 ... I_{18} I_{19} I_{20}

Uniformly sample each valid
motion segment into 20 frames

Segments shorter than
20 frames are discarded

4 2D Anchor-point Annotation

One anchor point per frame



p_1 p_2 p_3 ... p_{20}

20-point instrument trajectory



$$P = \{p_1, \dots, p_{20}\}, \quad p_t \in \mathbb{R}^2$$

5 Benchmark Sample

20 RGB frames



20 annotated
2D anchor points

p_1 p_2 p_3 ... p_{20}

1 instrument
trajectory



Dataset Statistics

1,637 valid motion segments

32,740 annotated frames

Avg. raw segment length: 79.42 frames

Min / Max: 21 / 585 frames

Motion Quality Distribution
(# segments)



High: 1,429 (87.3%)

Medium: 102 (6.2%)

Low: 106 (6.5%)