

IKrNet: A Neural Network for Detecting Specific Drug-Induced Patterns in Electrocardiograms Amidst Physiological Variability

Ahmad Fall¹ Federica Granese¹ Alex Lence¹

¹IRD, Sorbonne Université, Unité de Modélisation Mathématique et Informatique
des Systèmes Complexes, UMMISCO F-93143 Bondy, France

{ahmad.fall, federica.granese, alex.lence}@ird.fr

Dominique Fourer² Blaise Hanczar²

²University of Evry - Paris-Saclay, IBISC Laboratory Evry, France

{dominique.fourer, blaise.hanczar}@univ-evry.fr

Joe-Elie Salem^{3,4}

³Clinical Investigation Center Paris-Est, CIC-1901, INSERM, Department of Pharmacology
Pitié-Salpêtrière University Hospital, Sorbonne Université, France

⁴Department of Medicine, Vanderbilt University Medical Center
Nashville, TN, USA

joe-elie.salem@aphp.fr

Jean-Daniel Zucker^{1,5} Edi Prifti^{1,5}

⁵Sorbonne Université, INSERM, Nutrition et Obésités; systemic approaches,
NutriOmique, AP-HP Hôpital Pitié-Salpêtrière, France

{jean-daniel.zucker, edi.prifti}@ird.fr

Abstract

Monitoring and analyzing electrocardiogram (ECG) signals, even under varying physiological conditions, including those influenced by physical activity, drugs and stress, is crucial to accurately assess cardiac health. However, current AI-based methods often fail to account for how these factors interact and alter ECG patterns, ultimately limiting their applicability in real-world settings. This study introduces IKrNet, a novel neural network model, which identifies drug-specific patterns in ECGs amidst certain physiological conditions. IKrNet's architecture incorporates spatial and temporal dynamics by using a convolutional backbone with varying receptive field size to capture spatial features. A bi-directional Long Short-Term Memory module is also employed to model temporal dependencies. By treating heart rate variability as a surrogate for physiological fluctuations, we evaluated IKrNet's performance across diverse scenarios, including conditions with physical stress, drug intake alone, and a baseline without drug presence. Our assessment follows a clinical protocol in which 990 healthy volunteers were administered 80mg of Sotalol, a drug which

is known to be a precursor to Torsades-de-Pointes, a life-threatening arrhythmia. We show that IKrNet outperforms state-of-the-art models' accuracy and stability in varying physiological conditions, underscoring its clinical viability.

1. Introduction

An Electrocardiogram (ECG) is a diagnostic tool that records the heart's electrical activity through electrodes placed on the body, each producing an ECG lead. ECG leads provide multiple spatial and temporal perspectives on cardiac activities, aiding in diagnosing heart conditions.

In traditional clinical settings, ECGs are recorded over short, 10-second intervals while the patient rests. These recordings are commonly taken during in-clinic medical examinations and effectively detect persistent cardiac conditions. However, they may fail to capture transient or intermittent abnormalities outside controlled environments. For instance, long-term ECG monitoring is crucial for detecting life-threatening arrhythmias, such as Torsades-de-Pointes (TdP) [34, 37] — a specific form of polymorphic ventric-

ular tachycardia marked by rapid, irregular heartbeats that can lead to sudden cardiac arrest.

To overcome this limitation, wearable devices such as Holter monitors facilitate continuous ECG monitoring over 24-48 hours or longer, allowing one to observe real-world cardiac activity as patients go about their daily lives. However, continuous monitoring introduces additional challenges. Holter devices typically operate at a lower frequency range (180-250 Hz) [35], use fewer leads than the usual 12 used in clinical settings, and are susceptible to interference from factors like physical activity, stress, drugs and dietary influences [6]. These factors can affect the cardiovascular system, introducing distinct but complex ECG patterns that are further complicated by high noise levels and physiological variability over extended periods.

Manual analyses of Holter recordings are time-intensive, often requiring disentangling overlapping signal patterns. This challenge is exacerbated by heart rate variability (HRV), which is sensitive to physiological and pharmacological influences [22]. Our study assesses TdP risk by using Sotalol [15] as a pharmacological surrogate to induce QT prolongation, creating conditions that mimic the arrhythmogenic environment leading to TdP. To bridge the clinical challenges of detecting TdP in real-world settings with advancements in deep learning (DL), we focus on using DL models to improve the reliability of ECG analysis under varying physiological conditions.

DL has shown promise in ECG applications, such as noise reduction [8, 12], waveform segmentation [24], feature extraction [2], and abnormality detection [3]. Most architectures, including ResNet [13, 36], Long Short-Term Memory (LSTM) [14, 28], and Transformer [1, 38], have shown potential for automated ECG analyses and arrhythmia detection. With their residual connections, ResNet models are adept at capturing spatial patterns in ECG signals, while LSTMs are well-suited for capturing temporal dependencies in sequential data like ECG waveforms. Transformer models, on the other side, leverage self-attention mechanisms [33] to model long-range dependencies in sequential data, which is beneficial for ECG analysis.

However, these models face significant limitations when applied in real-world scenarios, particularly in the presence of substantial HRV. ResNet’s focus on spatial features limits its effectiveness in continuous temporal tracking, which is crucial for detecting subtle QT prolongations and arrhythmogenic patterns influenced by HRV. For instance, LSTMs struggle to adapt to complex interactions between physiological HRV and drug-induced effects, making them less robust under abrupt changes in ECG signals, as they do not have internal mechanisms to focus on specific parts of the waveform. Transformers, while capable of modeling global dependencies, require extensive data and are prone to noise and variability in ambulatory ECGs — es-

pecially at lower sampling rates typical of Holter monitors. Additionally, the computational demands of Transformer models often preclude their real-time or on-device application due to limited resources [38]. Previous studies, such as [26], have developed CNN-based models for predicting the presence of Sotalol as a surrogate for assessing drug-induced long QT syndromes (diLQTS) and congenital long QT syndromes (cLQTS), both of which are risk factors for Torsades-de-Pointes. Their most effective model, a DenseNet [17] with six blocks of dense convolutional layers, improved diLQTS and cLQTS detection but in expert hand-picked clean ECGs. However, the authors did not evaluate the impact of HRV on the model’s performance, leaving a critical gap in robustness for real-world applications. To address these limitations, our study advances the field by making the following contributions:

- 1. We formally describe our proposed framework**, designed to be flexible enough for application across various contexts. Additionally, we establish a robustness definition tailored to our needs, introducing a fairness-inspired metric that quantitatively assesses the model’s robustness.
- 2. We present a new architecture, IKrNet, specifically designed to detect drug-induced ECG patterns.** Our model leverages both spatial and temporal components of ECG signals through a multi-branch structure, using varying receptive field size to capture features across multiple scales. Spatial features are extracted, followed by bi-directional LSTM layers capturing temporal dependencies.
- 3. We evaluate both state-of-the-art models and our proposed model on a clinician-curated dataset and a real-world dataset obtained from Holter monitors. This comparison reveals the limitations of traditional models in handling data with physiological variability and different sampling rates. Additionally, we rigorously assess model performance by considering typically overlooked ECG features**, such as heart rate variability, which can significantly influence drug-related ECG patterns.

2. Framework

We introduce a new classifier, named IKrNet, designed to detect *drug footprints* from ECGs robustly. By **detecting drug footprint**, we refer to the classifier’s ability to identify unique patterns in ECG signals induced explicitly by the drug. By **robustly**, we mean that the classifier accurately classifies ECGs despite variations in interactions from other footprints, such as stress, physical activity, or natural physiological fluctuations.

We focused on the use case described in [29], where data from 990 healthy participants from the GENEREPOL study received 80mg of Sotalol, an IKr blocker. This drug inhibits the rapid component of the delayed rectifier potassium current (IKr), which is crucial for cardiac repolarization and leads to prolonged QT intervals on the ECG.

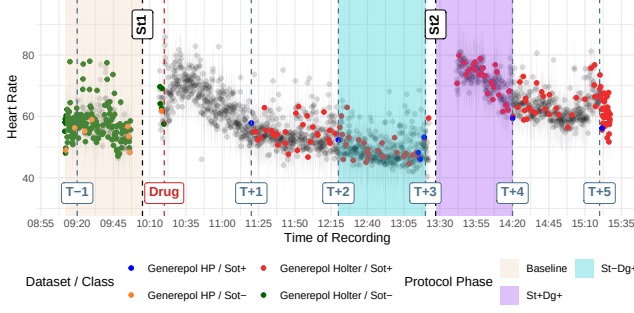


Figure 1. **Clinical protocol setup (the GENEREPOL study):** The heart rate is depicted for each of the 10s ECGs, recorded with a Holter device for one patient. Each ECG is depicted by a dot, which is coloured in grey when not used during the training, in green and blue respectively before drug intake (Sot-) and after drug intake (Sot+) for the Holter dataset and in orange and blue for the expert hand-picked dataset. Because of the unbalanced ECGs before and after drug intake, we selected a proportionate amount of ECGs for the training process (cf. Sec. 3.1)

This prolonged repolarization phase increases the risk of Torsade-de-Pointes. The IKr footprint refers in part to this QT prolongation, although other more subtle alterations are suspected and serves as a surrogate marker for TdP risk. Indeed, actual TdP events are rare and it is difficult to capture the risk before such events in a study setting, especially when vast datasets are needed to train DL models. Study participants were fitted with portable Holter monitors to record their hearts' electrical activity in real time before, during, and after Sotalol administration. As illustrated in Fig. 1, ECGs were recorded continuously over six hours, starting one or two hours before (here denoted as baseline) and continuing for up to five hours after Sotalol intake. Before taking the drug, participants underwent a stress test; they completed a second stress test three hours after taking the drug.

This section formalizes our framework, for broad applicability. We define key concepts and introduce the drug-footprint task as a supervised learning problem, along with the metrics used to evaluate the model's robustness.

2.1. Preliminaries

We consider the case of a single-lead ECG signal x (each observation of the heart activity is referred to as a lead). Here, x can be viewed as a sequence of K time-varying impulsive patterns x_i delimited by consecutive time instants t_{i-1} , t_i expressed in seconds, with $t_{i-1} < t_i$, such as:

$$x(t) = \begin{cases} x_1(t) & \text{if } t \in [t_0, t_1[\\ x_2(t) & \text{if } t \in [t_1, t_2[\\ \vdots & \\ x_K(t) & \text{if } t \in [t_{K-1}, t_K] \end{cases} \quad (1)$$

Each component x_i corresponds to a heart beat.

Discretization/Sampling process. Given a continuous-time signal $x(t)$, we can obtain its discretized version at sampling frequency F_s expressed in Hz, assuming a uniform sampling process as follows:

$$x[n] = \int_{\mathbb{R}} x(t) \delta(t - \frac{n}{F_s}) dt = x(\frac{n}{F_s}) \quad (2)$$

where δ is the Dirac distribution, $n \in \llbracket 0, N-1 \rrbracket$ denotes the time index, where $N = \lfloor t_K F_s \rfloor$, is the fine-length of the input vector $\mathbf{x} = [x[0], x[1], \dots, x[N-1]]$. According to the Nyquist-Shannon theorem, such signal enables a perfect reconstruction if $F_s \geq 2f_{\max}$, where $f_{\max}$ is the maximal frequency present in the signal.

2.2. Drug-footprint detection task

Let X be the random variable (r.v.) for the ECG recordings, with realizations $\mathbf{x} \in \mathcal{X} \subseteq \mathbb{R}^N$, where $\mathcal{X}$ is the space of ECG as defined in Sec. 2.1. Let Y be the r.v. representing the class label, taking values in $\mathcal{Y} = \{0, 1\}$, where 1 indicates the presence of the drug footprint, and 0 indicates its absence. The samples $(\mathbf{x}, y)$ are i.i.d. from P_{XY} .

The classifier $f_\theta : \mathcal{X} \rightarrow \mathcal{Y}$, parameterized by θ , maps an input ECG to a label, estimating whether the recording shows signs of the drug effect. We define this classifier as $f_\theta(\mathbf{x}) = \arg \max_{y \in \mathcal{Y}} p_{\hat{Y}|\mathcal{X}}(y|\mathbf{x}, \theta)$, where $p_{\hat{Y}|\mathcal{X}}(\cdot|\mathbf{x}, \theta)$ is the model's posterior distribution. In Sec. 3.2, we describe our proposed architecture tailored for drug-footprint detection and the rationale behind its design choices.

2.3. Evaluating model robustness

Finally, assuming that f_θ is trained without awareness of subgroup distinctions, we define a fairness-inspired notion of robustness to evaluate model stability across diverse recording conditions. Specifically, we frame this as an accuracy parity problem [4], where f_θ is considered robust across the subgroups $\mathcal{X}^{(1)}, \dots, \mathcal{X}^{(m)}$, which are subsets of $\mathcal{X}$, if for any subgroups $\mathcal{X}^{(i)}$ and $\mathcal{X}^{(j)}$:

$$|\mathbb{P}(Y = f_\theta(X^i) | X^i) - \mathbb{P}(Y = f_\theta(X^j) | X^j)| = 0,$$

where $X^\bullet$ is the r.v. whose realizations are in $\mathcal{X}^\bullet$ and Y is the true label.

Practical implementation. In our clinical use case, we consider the protocol time for identifying the subgroups we would like the model to achieve similar performances. By considering the protocol in Fig. 1, we recognise the following subgroups (or zones)¹, (i) Baseline – from the start

¹The Stress zone, which includes ECGs with a high heart rate and no drug administration, is omitted due to the lack of available data.

of the recording to the first stress test (StI); (ii) No Stress and Drug Present ($St-Dg_+$) – from $T+2h$ to $T+3h$, where no stress test is conducted but the drug is active; (iii) Stress and Drug Present ($St+Dg_+$) – immediately after the second stress test, where both the stress test and drug are present. In this context, we will call such subgroups protocol zones. We say f_θ is robust w.r.t. the protocol zones if²

$$|acc_{\max} - acc_{\min}| = 0,$$

where $acc_{\text{func}} \stackrel{\text{def}}{=} \text{func}\{acc(\text{Baseline}), acc(St-Dg_+), acc(St+Dg_+)\}$ with $\text{func} \in \{\min, \max\}$, and $acc(\mathcal{S}) \stackrel{\text{def}}{=} \frac{1}{|\mathcal{S}|} \sum_{(x,y) \in \mathcal{S}} \mathbb{1}[f_\theta(x) = y]$.

3. Experimental Setting

In the following section, we provide details of the experimental setup. In particular, Sec. 3.2 contains a thorough description of the proposed architecture.

3.1. Datasets and preprocessing

We used the Generepol-Holter and Generepol-Holter_{HP} datasets for our experiments. The Generepol-Holter dataset contains all the 10s consecutive ECGs from 990 participants recorded at a sampling rate of 500 Hz, following the clinical protocol presented in Fig. 1 and in [29]. The Generepol-Holter_{HP} dataset comprises a total of 15119 ECGs (15/participant) manually selected by expert cardiologists from the Generepol-Holter dataset. These hand picked (HP) ECGs are relatively clean (*in the clinical sense*) and contain a low noise level. However, they do not accurately reflect real-life scenarios, as most recorded ECGs are subject to different perturbations, including noise. Both datasets label ECGs as Sot- (before Sotalol intake) or Sot+ (after) and are divided into training (70%), validation (10%), evaluation (10%, to pick the best hyper-parameters configuration), and holdout (10%) partitions. To prevent class imbalance in Generepol-Holter, an equal number of Sot+ and Sot- ECGs were randomly selected per study participant. Indeed, Sot- ECGs are only recorded for a maximum of 2 hours, whereas Sot+ ECGs are recorded for 5 hours. Each participant had approximately the same amount of selected ECGs in the different protocol times ($T+1$, $T+2$, $T+3$, $T+4$, $T+5$ in Fig. 1) for label Sot+. The Generepol-Holter_{HP} training partition contained 6315 ECGs, and 61544 ECGs for the Generepol-Holter dataset.

Data augmentation techniques. To enhance data variability and address the varying sampling rates of recording devices, which typically operate at $F_s \in$

180, 250, 500, 1000 Hz, we applied a sampling rate augmentation technique [11, 20]. This process involves downsampling an original ECG signal (sampled at 500 Hz) to a lower rate, $F'_s \in 180, 250$ Hz, and then upsampling it back to 500 Hz for the training, validation, and evaluation sets. This simulates ECG recordings initially captured at lower sampling rates and then upsampled, reflecting real-world variability in ECG data [5]. We used spline interpolation for both downsampling and upsampling, in line with standard ECG processing techniques [32]. After data augmentation, the Generepol-Holter_{HP} training partition contained 18945 ECGs, and the Generepol-Holter dataset 184631 ECGs. Detailed number of ECGs per partition is provided in Appendix Tab. 3.

Data preprocessing. We standardized the ECG signals before inputting them into the network. This process aligns input feature scales, which is essential to prevent the model from applying disproportionate weight to ECGs with higher amplitudes [31]. By doing so, we mitigate the influence of amplitude variations that may arise from factors unrelated to the underlying cardiac conditions, such as differences in recording equipment, or individual patient characteristics.

3.2. IKrNet: our novel proposed architecture

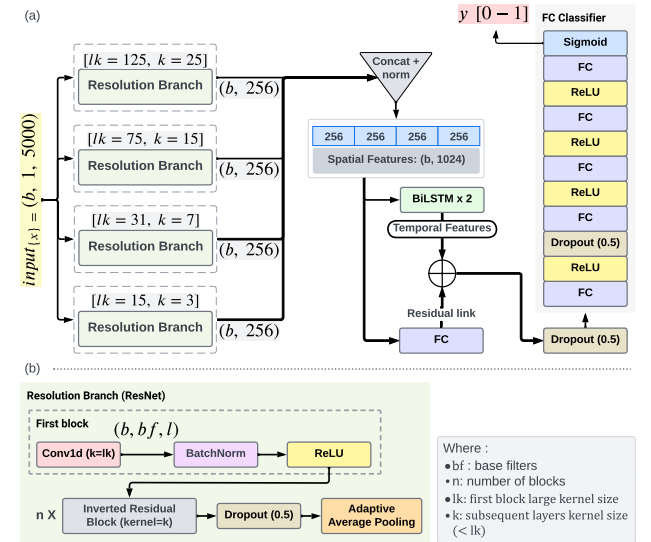


Figure 2. (a) Global architecture of IKrNet, featuring a CNN backbone with multiple branches, each capturing low- and high-frequency features through different receptive fields. (b) Detailed branch structure with inverted residual blocks, optimizing feature extraction and reducing complexity.

In Fig. 2, we illustrate the IKrNet architecture. A CNN backbone [21] first extracts spatial patterns from the raw ECG signal, using multiple branches with varying receptive field widths to capture features across different scales.

²The optimal scenario is identical accuracy across protocol zones. In practice, we consider a model more robust if this difference is *reasonably small*.

After spatial feature extraction, bi-directional Long Short-Term Memory (BiLSTM) layers [7, 18] analyze temporal dependencies in the ECG signal. Both spatial and temporal features are summed to produce a consolidated feature vector of the ECG, then fed to a fully connected classifier, which leverages the combined feature set to predict the occurrence of the drug footprint.

Spatial feature extractor. It comprises multiple resolution branches that process the same signal input, each with a distinct receptive field, enabling a multi-resolution analysis that adapts to varying ECG sampling rates and captures both low- and high-frequency variations. This design enhances the model’s ability to detect subtle and pronounced drug-induced changes, such as QT prolongation and heart rate changes associated with the drug. Within each branch, we employ a ResNet [13, 36] module known for its ability to facilitate the training of deeper network structures through identity shortcut connections [13]. These connections help mitigate the vanishing gradient problem by allowing gradients to flow directly through the network layers.

Each ResNet begins with a convolution layer with a large initial kernel size (lk), designed to capture macro features. This is followed by a series of residual blocks that employ a smaller, consistent kernel size (k), extracting more profound and detailed patterns. The distinct kernel pair (lk, k) varies across each branch, enabling the architecture to achieve multi-resolution analysis effectively. A larger kernel (lk) captures information over longer distances, while a smaller kernel (k) focuses on nearby information. Instead of basic blocks [13] as in the original ResNet architecture, we used inverted residual structures [30] as described in Appendix Fig. 6, which enhance model efficiency. They significantly reduce the trainable parameters, thus reducing computational costs compared to traditional basic blocks by limiting the scaling of parameters with the kernel size and number of channels. Squeeze-and-excite layers [16] (see Appendix Fig. 7), are integrated to recalibrate channels dynamically, compressing feature maps via global average pooling and modulating them with a fully connected layer and sigmoid activation. This process enhances sensitivity to clinically relevant ECG patterns. Features from each branch are then pooled with adaptive average pooling to reduce dimensionality and consolidate key information for temporal analysis. Initial inputs undergo progressive downsampling through convolution strides.

Temporal feature extractor. After extracting spatial features, the model processes temporal features to learn the dynamic nature of ECGs. This requires recognizing static patterns and capturing their temporal evolution to predict the drug intake footprint. Two BiLSTM layers process the data bidirectionally, capturing dependencies across the se-

quence to characterize cardiac events. Each BiLSTM layer integrates the concatenated spatial features extracted by the CNNs, building a complete temporal profile of the ECG.

Drug footprint detection. At this stage, we combine the temporal and the spatial features, which involves a skip connection strategy that directly adds the output from the fully connected layer to the output of the last BiLSTM layer. This approach preserves both high-resolution spatial details and deep temporal relationships. The final integration of features is then fed into a fully connected classifier that produces a score between 0 and 1.

Training process We trained our model using the binary cross-entropy loss function alongside the AdamW optimizer [23], an improved version of the traditional Adam optimizer, introducing decoupled weight decay regularization allowing better regularization, and reduced overfitting by penalizing large weights. The training was conducted over ten epochs, using 4 NVIDIA A100 GPUs, each with 80 GB of memory, in approximately 4 hours. The model exhibited rapid convergence throughout the training, with both training and validation losses decreasing steadily and remaining closely aligned, suggesting the model does not overfit the training data. A hyperparameter optimization process and an ablation study (cf. Appendix C.2) were conducted to identify the optimal hyperparameter set. The best-performing model, selected based on evaluation metrics, had the following hyperparameters: batch size of 64 per GPU (256 globally), a learning rate of 0.001, 4 branches in the CNN backbone with kernel pairs (lk, k) (125, 25), (75, 15), (31, 7), (15, 3). Each CNN branch had the following strides [1, 5, 1, 5, 1, 4, 1, 4, 1, 3], an initial convolution filter of 64, and 10 inverted residual blocks. The convolution filters were increased by 2 after every 4 block. The final architecture has ~ 24 million parameters.

4. Experimental Results and Analysis

This section describes and evaluates eight models, summarized in Appendix Tab. 2. Specifically, we assess the performance of the state-of-the-art (SOTA) model, DenseNet₅₀₀^{HP} for Sotalol footprint detection, as described in [26], concerning its behavior on ECGs downsampled at different sampling frequencies and its robustness across the protocol zone analyzed in Sec. 2.2. In their paper, [26] trained and tested on Generepol-Holter_{HP}. Their model was trained for 100 epochs using the Adam optimizer, with a learning rate of 0.001 and a dropout rate of 0.2. We considered three variations of the model based on the DenseNet architecture. These models are trained on Generepol-Holter_{HP} with data augmentation (DenseNet₊^{HP}), Generepol-Holter alone (DenseNet₅₀₀), and Generepol-Holter with data

Table 1. **Cross-dataset model accuracies and Accuracy Parity Difference (APD) by protocol zone across sampling rates:** The results are presented as mean $\pm$ standard deviation for ECGs sampled at 150, 180, 215, 250, 300, 350, 425, and 500 Hz.

	Generepol-Holter _{HP}				Generepol-Holter			
	Baseline	St-Dg ₊	St ₊ Dg ₊	-	Baseline	St-Dg ₊	St ₊ Dg ₊	-
	Accuracy ($\uparrow$)	Accuracy ($\uparrow$)	Accuracy ($\uparrow$)	APD ($\downarrow_0$)	Accuracy ($\uparrow$)	Accuracy ($\uparrow$)	Accuracy ($\uparrow$)	APD ($\downarrow_0$)
Models trained on Generepol-Holter_{HP}								
DenseNet ₅₀₀ ^{HP}	88.82 $\pm$ 8.98	95.08 $\pm$ 2.68	98.61 $\pm$ 0.41	11.71 $\pm$ 10.18	49.44 $\pm$ 14.27	99.64 $\pm$ 0.19	96.96 $\pm$ 1.41	59.49 $\pm$ 19.74
DenseNet ₊ ^{HP}	91.64 $\pm$ 5.27	94.15 $\pm$ 2.01	98.92 $\pm$ 0.7	8.7 $\pm$ 5.96	65.13 $\pm$ 12.83	99.59 $\pm$ 0.15	93.98 $\pm$ 1.49	39.37 $\pm$ 16.84
IKrNet ₅₀₀ ^{HP}	92.87 $\pm$ 3.86	98.0 $\pm$ 0.37	97.65 $\pm$ 0.66	5.68 $\pm$ 4.38	78.43 $\pm$ 6.43	99.2 $\pm$ 0.12	92.97 $\pm$ 0.89	22.75 $\pm$ 7.58
IKrNet ₊ ^{HP}	94.92 $\pm$ 0.92	94.21 $\pm$ 1.54	98.37 $\pm$ 0.21	4.53 $\pm$ 1.59	79.89 $\pm$ 2.08	99.08 $\pm$ 0.13	92.31 $\pm$ 0.59	20.91 $\pm$ 2.33
Models trained on Generepol-Holter								
DenseNet ₅₀₀	94.73 $\pm$ 1.43	90.23 $\pm$ 1.36	94.9 $\pm$ 0.78	5.52 $\pm$ 1.31	85.68 $\pm$ 2.05	96.39 $\pm$ 0.63	86.23 $\pm$ 0.56	12.67 $\pm$ 1.28
DenseNet ₊	96.28 $\pm$ 0.69	89.9 $\pm$ 1.1	95.42 $\pm$ 0.54	6.97 $\pm$ 1.21	89.74 $\pm$ 1.85	97.17 $\pm$ 0.38	89.37 $\pm$ 0.42	9.31 $\pm$ 0.73
IKrNet ₅₀₀	97.57 $\pm$ 1.6	91.89 $\pm$ 0.25	97.21 $\pm$ 0.65	6.52 $\pm$ 0.53	89.69 $\pm$ 4.07	99.1 $\pm$ 0.12	96.08 $\pm$ 0.61	9.88 $\pm$ 4.45
IKrNet ₊	97.72 $\pm$ 1.28	95.48 $\pm$ 0.4	98.65 $\pm$ 0.41	3.39 $\pm$ 0.7	92.22 $\pm$ 1.66	98.92 $\pm$ 0.19	95.35 $\pm$ 0.27	6.98 $\pm$ 1.67

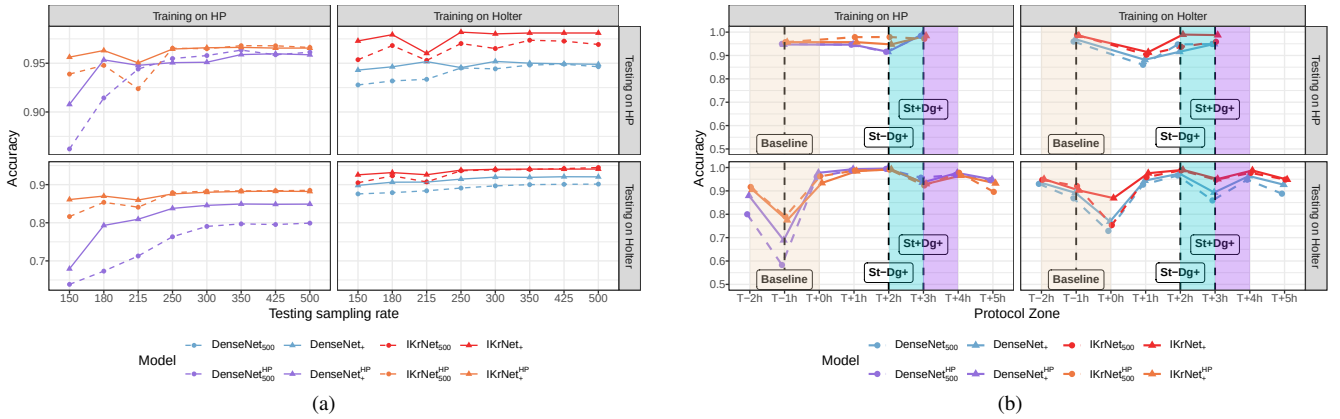


Figure 3. **Cross-dataset evaluation of IKrNet and DenseNet** on different sampling rates (a), and across full protocol at 500 Hz (b)

augmentation (DenseNet₊). Finally, we evaluated the proposed IKrNet architecture, considering models trained on Generepol-Holter_{HP} (IKrNet₅₀₀^{HP}), Generepol-Holter_{HP} with data augmentation (IKrNet₊^{HP}), Generepol-Holter alone (IKrNet₅₀₀), and Generepol-Holter with data augmentation (IKrNet₊).

4.1. Performance analysis of models

Analysis on Generepol-Holter_{HP} and Generepol-Holter.

The results in Fig. 3a show that when tested as in the original setup on the Generepol-Holter_{HP} dataset at 500 Hz, DenseNet₅₀₀^{HP} performed very closely (within a one percentage point difference) to the results claimed in the original paper, achieving an accuracy of 95.77%, a precision of 96.03%, a recall of 96.41%, and an F1 score of 96.22%.

Those results also reveal a drop in performance when DenseNet₅₀₀^{HP} is tested on the Generepol-Holter dataset, regardless of the sampling frequency (~ 2 percentage points difference between the accuracy of DenseNet₅₀₀^{HP} on ECGs at 500 Hz from Generepol-Holter_{HP} and the same setting in Generepol-Holter). However, DenseNet₅₀₀

is more stable across multiple sampling rates and shows higher accuracy than DenseNet₅₀₀^{HP} when tested on Generepol-Holter. From Fig. 3a, we also observe that IKrNet₅₀₀^{HP} and IKrNet₅₀₀ display higher accuracies compared to previous models and have similar performance on multiple sampling rates. When trained and tested on the Generepol-Holter_{HP} dataset, IKrNet₅₀₀^{HP} achieved an accuracy of 96.66%, surpassing DenseNet₅₀₀^{HP}'s 95.77%. A similar improvement was observed on the Generepol-Holter dataset, where IKrNet₅₀₀ reached 94.37% in accuracy (higher than the 89.98% accuracy of DenseNet₅₀₀). In Appendix Tabs. 9 and 11, we provide performance results on cross-dataset evaluations.

Analysis on data-augmented ECGs.

DenseNet models do not show consistency when tested in settings with variability in sampling frequency (cf. Fig. 3a and Appendix Tabs. 5 and 7). Specifically, we observed a decrease in accuracy performance by up to approximately 9% when resampling the signal from 500 Hz to 150 Hz then to 500 Hz. Similarly, the F1 score decreased up to ap-

proximately 7 percentage points. This performance drop could be expected due to the inevitable loss of information during downsampling. To assess whether the performance drop of $\text{DenseNet}_{500}^{\text{HP}}$ was due to the architecture itself or some epistemic uncertainty caused by the lack of data in the training set, we retrained the models by considering ECGs resampled at different sampling frequencies, as explained in Sec. 3.1. The numerical results for this second set of experiments are in Appendix Tabs. 10 and 12.

The experiments confirmed that augmenting the training dataset significantly improves the model’s generalization capability across different sampling rates. On ECGs at 150 Hz, accuracy increased from 86.08% with $\text{DenseNet}_{500}^{\text{HP}}$ to 90.31% with $\text{DenseNet}_{+}^{\text{HP}}$ when testing on Generepol-Holter_{HP}. However, the results on Generepol-Holter remain unsatisfactory, with $\text{DenseNet}_{+}^{\text{HP}}$ failing to achieve an overall accuracy greater than 70% (cf. Fig. 5). Conversely, DenseNet_{500} and DenseNet_{+} , benefiting from the diversity in the Generepol-Holter training set, achieved accuracies of 87.41% and 89.58% on Generepol-Holter at 150 Hz, respectively. The accuracy of IKrNet improved by considering more sampling rates when training. For instance, $\text{IKrNet}_{500}^{\text{HP}}$ accuracy increased from 93.43% to 95.10% with $\text{IKrNet}_{+}^{\text{HP}}$ when testing on Generepol-Holter_{HP}. On the Generepol-Holter dataset, IKrNet_{500} accuracy improved from 90.30% to 92.46% at 150 Hz with IKrNet_{+} model. At 500 Hz on the same dataset, IKrNet_{+} has an accuracy of **94.00%** outperforming DenseNet_{+} ’s 91.90%, and when tested on the Generepol-Holter_{HP} dataset, they achieved accuracies of **97.88%** and **94.77%** respectively. IKrNet demonstrated greater performance compared to DenseNet , consistently outperforming it across different datasets and conditions.

Robustness to the protocol zones. Tab. 1 shows the numerical results across the protocol zones (cf. Sec. 2.2). The Accuracy Parity Difference (APD) notably revealed model instability across protocol zones. On Generepol-Holter_{HP}, all models performed similarly, with an APD difference of less than 7% between the worst-performing model ($\text{DenseNet}_{500}^{\text{HP}}$) and the best (IKrNet_{+}). In contrast, on Generepol-Holter, models, particularly those trained on Generepol-Holter_{HP}, show a significant performance decline. $\text{DenseNet}_{500}^{\text{HP}}$ exhibits a nearly 60% accuracy gap between the Baseline and St-Dg₊ results. Interestingly, IKrNet_{+} consistently outperformed the other models in this protocol time zone. This can be explained by several factors: (i) Although the training set from Generepol-Holter contains more perturbations compared to Generepol-Holter_{HP}, it helps the model to better generalize and reduces overfitting to the specific dataset. This is further supported by the fact that *all* models trained on this dataset generally show better performance. (ii) The data augmentation

technique also seems to play a key role, particularly as the models were tested on various sampling frequencies that were not present during training. (iii) The multi-resolution CNN in IKrNet which can adapt to different waveform segments helps it in detecting baseline variations, including subtle QT prolongation changes. In contrast, the fixed kernel size of DenseNet limits its ability to focus on shorter or longer sections of the waveform.

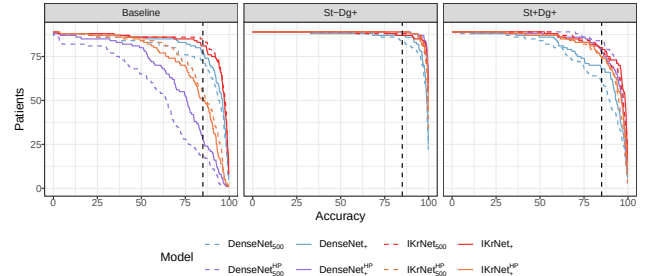


Figure 4. **Results on Generepol-Holter (ECGs at 500 Hz) over protocol zones and patients.** Number of patients (vertical axis) achieving accuracy above the threshold (horizontal axis).


In Fig. 4, to further investigate the model performances from a patient perspective, we considered all the ECGs recorded during the specified protocol periods for each model we measured the number of patients achieving accuracy above the examined threshold accuracy value (horizontal axis). The IKrNet models trained on Generepol-Holter exhibit consistent performance across study participants, with more than 75 out of 90 (83%), in the holdout set, achieving an accuracy higher than 85% (red lines). In contrast, the other models showed more variability. Notably, $\text{DenseNet}_{500}^{\text{HP}}$ performed poorly in the Baseline zone, with fewer than 25 out of 90 (28%) participants with an accuracy above 85% (violet dashed line).

$\text{DenseNet}_{500}^{\text{HP}}$ mistakes the most when considering the ECGs recorded one hour before Sotalol intake, as seen in Fig. 5b. During this period, the model’s accuracy was less than 85% for fewer than 25 out of 90 participants (27%).


Interestingly, some of the $\text{DenseNet}_{500}^{\text{HP}}$ limitations were also reported in Granese et al. [10], where it was demonstrated that the model lacked robustness when ECGs are processed through denoising methods, even after re-training.

4.2. Heart rate effect on drug footprint inference


Heart rate variability (HRV) [27] is an important factor that can affect model performance, particularly in real-world scenarios where physiological states fluctuate over time [19]. Sotalol is known to lower heart rate [15] (cf. Fig. 1, zone St-Dg₊), an effect that models can typically detect. However, sudden increases in heart rate due to external factors among which stress efforts (cf. Fig. 1,

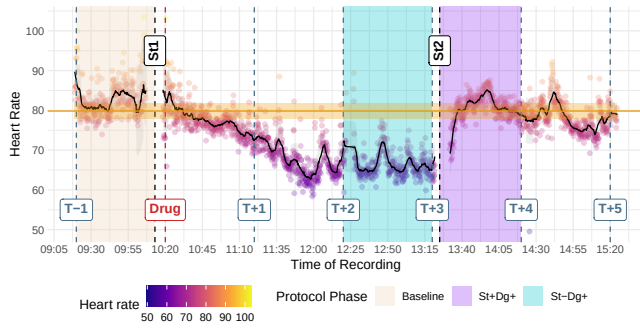


(a) Heart rate variability.

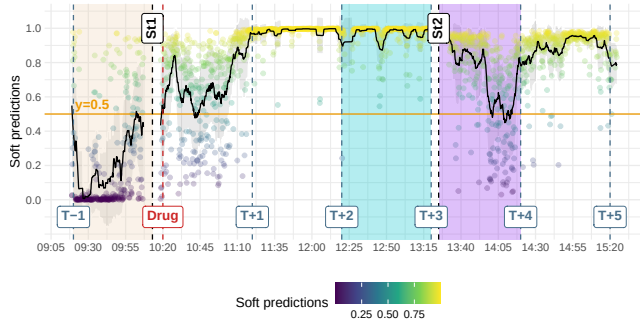


(b) DenseNet₅₀₀^{HP} predictions.

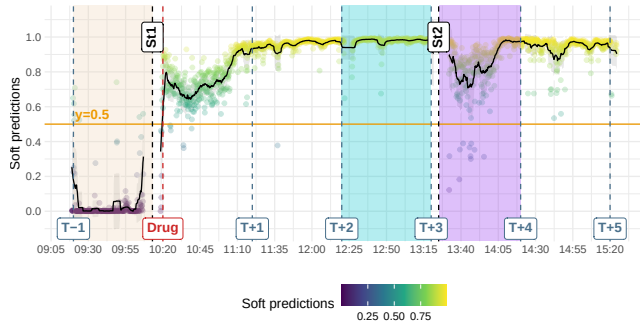


(c) IKrNet₊ predictions.

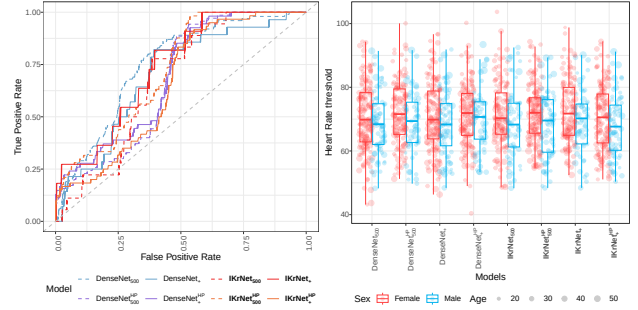

(d) ROC curves with heart rate as a predictor of misclassifications (e) Heart Rate threshold per model on all patients

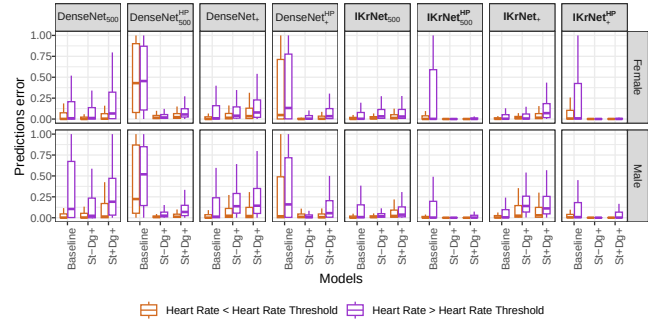

(f) Model's prediction errors by heart rate groups on all sampling rates

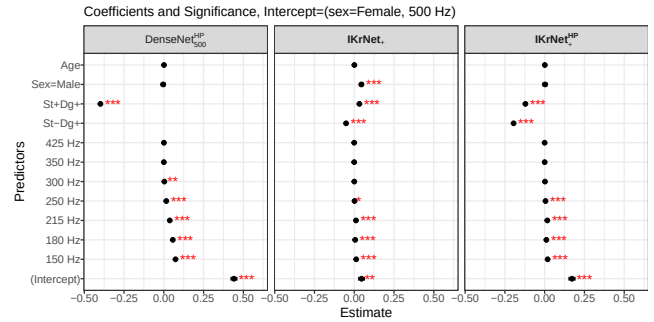

(g) LME models [25]: $\text{predictions.error} \sim \text{sampling.rate} + \text{sex} + \text{age} + (1|\text{patient})$

Figure 5. **DenseNet₅₀₀^{HP} and IKrNet₊ on single patient and influence of heart rate variability across time:** DenseNet₅₀₀^{HP} is negatively affected by sampling rate changes and stress zones where heart rate increases, whereas IKrNet₊ is more robust.

zone St+Dg+), can obscure the drug's effect, potentially confusing the models and impacting their predictive accuracy. In Fig. 5, we illustrate a comparison between DenseNet₅₀₀^{HP} and IKrNet₊, focusing on data from a single study participant across the entire protocol duration. We assessed how HRV impacted both models' performance and how effectively each model can adapt to varying heart rate conditions.

Identifying optimal heart rate (OHR) thresholds We identified heart rate thresholds by analyzing the ROC curves for each model, as shown in Fig. 5d, using heart rate as a predictor of the models' misclassification errors. By calcu-

lating the Youden Index [9] on the ROC curves, which maximize the difference between sensitivity and specificity, we determined the optimal heart rate threshold for each model and each study participant, beyond which models start to struggle. Interestingly, we observed a high variability of OHR throughout the study participants, as illustrated by boxplots in Fig. 5e. This was partially explained by gender and age.

Predictions errors analysis Additionally, we computed the prediction errors for the ECGs above the OHR and below and observed that the predicted error was systematically higher above the OHR. Overall this error was shown to be

lower for the IKrNet models compared with DenseNet Fig. 5f. Furthermore, we computed Linear Mixed Effects models (LME) [25] to further assess the impact of sampling rates and stress on the models. IKrNet models show relatively stable performance across low- and high-heart rate conditions in all zones, particularly at $St+Dg+$ where drug influence and physiological stress are present. In contrast, DenseNet models exhibit increased variability, with more prediction errors at Baseline but also at higher heart rates during the $St-Dg+$ and $St+Dg+$ zones.

5. Conclusion

In this work, we presented IKrNet, a novel neural network architecture, to detect drug-induced patterns in ECGs. This architecture is tailored explicitly for ECG analysis, leveraging a multi-resolution CNN backbone with varying receptive fields to capture low- and high-frequency features. By integrating spatial feature extraction through the CNN backbone with temporal dependency analysis via BiLSTM layers, IKrNet effectively adapts to the complex, time-evolving patterns, which are typical of ambulatory ECG Holter recordings. These are affected by varying heart rate and an integration of different footprints induced by external stressors. We considered a clinical use case where study participants were administered Sotalol, a pharmacological agent used to induce QT prolongation and mimic arrhythmogenic conditions that can lead to Torsades-de-Pointes. Our contributions include establishing a formal framework for model robustness to evaluate consistency quantitatively under heart rate variability and stress conditions. This framework offers a standardized approach to assess model reliability across diverse and fluctuating physiological states, which is critical for real-world applications.

Acknowledgments

This study was supported by the ANR-20-CE17-0022 DeepECG4U funding from the French National Research Agency (ANR).

References

- [1] Taymaz Akan, Sait Alp, and Mohammad Alfrad Nobel Bhuiyan. Ecgformer: Leveraging transformer for ecg heart-beat arrhythmia classification. In *2023 International Conference on Computational Science and Computational Intelligence (CSCI)*, pages 1412–1417. IEEE, 2023.
- [2] Ridwan Alam, Aaron D Aguirre, and Collin M Stultz. Qtnet: Deep learning for estimating qt intervals using a single lead ecg. In *2023 45th Annual International Conference of the IEEE Engineering in Medicine & Biology Society (EMBC)*, pages 1–4. IEEE, 2023.
- [3] Yaqoob Ansari, Omar Mourad, Khalid Qaraqe, and Erchin Serpedin. Deep learning for ecg arrhythmia detection and classification: an overview of progress for period 2017–2023. *Frontiers in Physiology*, 14:1246746, 2023.
- [4] Solon Barocas, Moritz Hardt, and Arvind Narayanan. *Fairness and machine learning: Limitations and opportunities*. MIT press, 2023.
- [5] Mathias Baumert, Martin Schmidt, Sebastian Zaunseder, and Alberto Porta. Effects of ecg sampling rate on qt interval variability measurement. *Biomedical Signal Processing and Control*, 25:159–164, 2016.
- [6] Andrea Bravi, Geoffrey Green, Christophe Herry, Heather E Wright, André Longtin, Glen P Kenny, and Andrew JE Seely. Do physiological and pathological stresses produce different changes in heart rate variability? *Frontiers in physiology*, 4: 197, 2013.
- [7] Jinyong Cheng, Qingxu Zou, and Yunxiang Zhao. ECG signal classification based on deep CNN and BiLSTM. *BMC Med Inform Decis Mak*, 21(1):365, 2021.
- [8] Hsin-Tien Chiang, Yi-Yen Hsieh, Szu-Wei Fu, Kuo-Hsuan Hung, Yu Tsao, and Shao-Yi Chien. Noise reduction in ecg signals using fully convolutional denoising autoencoders. *IEEE Access*, 7:60806–60813, 2019.
- [9] Ronen Fluss, David Faraggi, and Benjamin Reiser. Estimation of the youden index and its associated cutoff point. *Biometrical Journal*, 47(4):458–472, 2005.
- [10] Federica Granese, Ahmad Fall, Alex Lence, Joe-Elie Salem, Jean-Daniel Zucker, and Edi Prifti. The negative impact of denoising on automated classification of electrocardiograms. In *Deep Generative Models for Health Workshop NeurIPS 2023*, 2023.
- [11] Marjan Gusev, Stojancho Tudjarski, and Ana Anagelevska. Sampling rate impact on heart rate variability. In *2022 30th Telecommunications Forum (TELFOR)*, pages 1–4, 2022.
- [12] C. Haritha, M. Ganesan, and E. P. Sumesh. A survey on modern trends in ecg noise removal techniques. In *2016 International Conference on Circuit, Power and Computing Technologies (ICCPCT)*, pages 1–7, 2016.
- [13] Kaiming He, Xiangyu Zhang, Shaoqing Ren, and Jian Sun. Deep residual learning for image recognition, 2015.
- [14] Sepp Hochreiter and Jürgen Schmidhuber. Long short-term memory. *Neural Comput.*, 9(8):1735–1780, 1997.
- [15] Stefan H. Hohnloser, Thomas Klenghagen, Markus Zabel, and Hanjörg Just. Effect of sotalol on heart rate variability assessed by holter monitoring in patients with ventricular arrhythmias. *The American Journal of Cardiology*, 72(4): A67–A71, 1993. Sotalol in Life-Threatening Ventricular Arrhythmias: A Unique Class III Antiarrhythmic.
- [16] Jie Hu, Li Shen, and Gang Sun. Squeeze-and-excitation networks. In *2018 IEEE/CVF Conference on Computer Vision and Pattern Recognition*, pages 7132–7141, 2018.
- [17] Gao Huang, Zhuang Liu, Laurens van der Maaten, and Kilian Q. Weinberger. Densely connected convolutional networks, 2018.
- [18] Zhiheng Huang, Wei Xu, and Kai Yu. Bidirectional lstm-crf models for sequence tagging, 2015.
- [19] Yamini Kotriwar, Sneha Kachhara, K. P. Harikrishnan, and G. Ambika. Higher order spectral analysis of ecg signals, 2018.

- [20] Ohhwan Kwon, Jinwoo Jeong, Hyung Bin Kim, In Ho Kwon, Song Yi Park, Ji Eun Kim, and Yuri Choi. Electrocardiogram sampling frequency range acceptable for heart rate variability analysis. *Health Inform Res*, 24(3):198–206, 2018.
- [21] Yann LeCun, Yoshua Bengio, and Geoffrey Hinton. Deep learning. *Nature*, 521(7553):436–444, 2015.
- [22] Xiangchen Li, Yuting Song, Huang Wang, Xinyu Su, Mengyao Wang, Jing Li, Zhiqiang Ren, Daidi Zhong, and Zhiyong Huang. Evaluation of measurement accuracy of wearable devices for heart rate variability. *iScience*, 26(11):108128, 2023.
- [23] Ilya Loshchilov and Frank Hutter. Decoupled weight decay regularization, 2019.
- [24] Viktor Moskalenko, Nikolai Zolotykh, and Grigory Osipov. Deep learning for ecg segmentation. In *international conference on Neuroinformatics*, pages 246–254. Springer, 2019.
- [25] José C Pinheiro and Douglas M Bates. Linear mixed-effects models: basic concepts and examples. *Mixed-effects models in S and S-Plus*, pages 3–56, 2000.
- [26] Edi Prifti, Ahmad Fall, Giovanni Davogustto, Alfredo Pulini, Isabelle Denjoy, Christian Funck-Brentano, Yasmin Khan, Alexandre Durand-Salmon, Fabio Badilini, Quinn S Wells, et al. Deep learning analysis of electrocardiogram for risk prediction of drug-induced arrhythmias and diagnosis of long qt syndrome. *European Heart Journal*, 42(38):3948–3961, 2021.
- [27] U Rajendra Acharya, K Paul Joseph, N Kannathal, Choo Min Lim, and Jasjit S Suri. Heart rate variability: a review. *Medical and Biological Engineering and Computing*, 44(12):1031–1051, 2006.
- [28] Moumita Roy, Sukanta Majumder, Anindya Halder, and Utpal Biswas. Ecg-net: A deep lstm autoencoder for detecting anomalous ecg. *Engineering Applications of Artificial Intelligence*, 124:106484, 2023.
- [29] Joe-Elie Salem, Marine Germain, Jean-Sébastien Hulot, Pascal Voiriot, Bruno Lebourgeois, Jean Waldura, David-Alexandre Tregouet, Beny Charbit, and Christian Funck-Brentano. Genome wide analysis of sotalol-induced ikr inhibition during ventricular repolarization, “generepol study”: Lack of common variants with large effect sizes. *PLoS One*, 12(8):e0181875, 2017.
- [30] Mark Sandler, Andrew Howard, Menglong Zhu, Andrey Zhmoginov, and Liang-Chieh Chen. Mobilenetv2: Inverted residuals and linear bottlenecks, 2019.
- [31] M. Shanker, M.Y. Hu, and M.S. Hung. Effect of data standardization on neural network training. *Omega*, 24(4):385–397, 1996.
- [32] Khairul Azami Sidek and Ibrahim Khalil. Enhancement of low sampling frequency recordings for ecg biometric matching using interpolation. *Computer Methods and Programs in Biomedicine*, 109(1):13–25, 2013.
- [33] Ashish Vaswani, Noam Shazeer, Niki Parmar, Jakob Uszkoreit, Llion Jones, Aidan N. Gomez, Lukasz Kaiser, and Illia Polosukhin. Attention is all you need, 2023.
- [34] Sami Viskin. Long qt syndromes and torsade de pointes. *The Lancet*, 354(9190):1625–1633, 1999.
- [35] Eric S Winokur, Maggie K Delano, and Charles G Sadini. A wearable cardiac monitor for long-term data acquisition and analysis. *IEEE Trans. Biomed. Eng.*, 60(1):189–192, 2013.
- [36] Jianxin Xie, Stavros Stavrakis, and Bing Yao. Automated identification of atrial fibrillation from single-lead ecgs using multi-branching resnet. *Frontiers in Physiology*, 15, 2024.
- [37] Yee Guan Yap and A John Camm. Drug induced qt prolongation and torsades de pointes. *Heart*, 89(11):1363–1372, 2003.
- [38] Zibin Zhao. Transforming ecg diagnosis: an in-depth review of transformer-based deep learning models in cardiovascular disease detection, 2023.

IKrNet: A Neural Network for Detecting Specific Drug-Induced Patterns in Electrocardiograms Amidst Physiological Variability

Supplementary Material

A. Heart rate estimation

The heart rate refers to the number of times the heart beats per minute (bpm). It can be estimated from the time difference between two successive heartbeats as follows:

$$\text{HR}(\mathbf{x}_i, \mathbf{x}_{i+1}) = \frac{60F_s}{n_{R_{i+1}} - n_{R_i}} \quad (3)$$

with, $n_{R_i} = \arg \max_{n \in [n_{i-1}, n_i]} |\mathbf{x}_i[n]|$

n_{R_i} being the sample index of the R impulse local peak position in the signal pattern x_i . This value can be averaged over the whole signal made of K observed heartbeats as:

$$\text{HR}_{\text{avg}}(\mathbf{x}) = \frac{1}{K-1} \sum_{i=1}^{K-1} \text{HR}(\mathbf{x}_i, \mathbf{x}_{i+1}). \quad (4)$$

B. Additional material to Sec. 3

B.1. Notation

Table 2 provides an overview of the model names referenced throughout the paper.

Table 2. Table of notations for the experiments.

Symbol	Description
$\text{DenseNet}_{\{\cdot\}}^{\text{HP}}$	SOTA DenseNet trained on Generepol-Holter _{HP} with ECGs at 500 Hz as in [26] ($\{\cdot\} = 500$), or [180, 215 and 500] Hz ($\{\cdot\} = +$).
$\text{DenseNet}_{\{\cdot\}}$	SOTA DenseNet trained on Generepol-Holter with ECGs at 500 Hz ($\{\cdot\} = 500$), or [180, 215 and 500] Hz ($\{\cdot\} = +$).
$\text{IKrNet}_{\{\cdot\}}^{\text{HP}}$	IKrNet trained on Generepol-Holter _{HP} with ECGs at 500 Hz ($\{\cdot\} = 500$), or [180, 215 and 500] Hz ($\{\cdot\} = +$).
$\text{IKrNet}_{\{\cdot\}}$	IKrNet trained on Generepol-Holter with ECGs at 500 Hz ($\{\cdot\} = 500$), or [180, 215 and 500] Hz ($\{\cdot\} = +$).

B.2. Datasets

Table 3 describes the dataset per partition. The number of samples in the table refers to the total amount after data augmentation as described in Sec. 3.1.

B.3. IKrNet architecture

Figures 6 and 7 provide a detailed view of the Inverted Residual Block and Squeeze-and-Excite Block modules. We refer to Sec. 3.2 for the general architecture description.

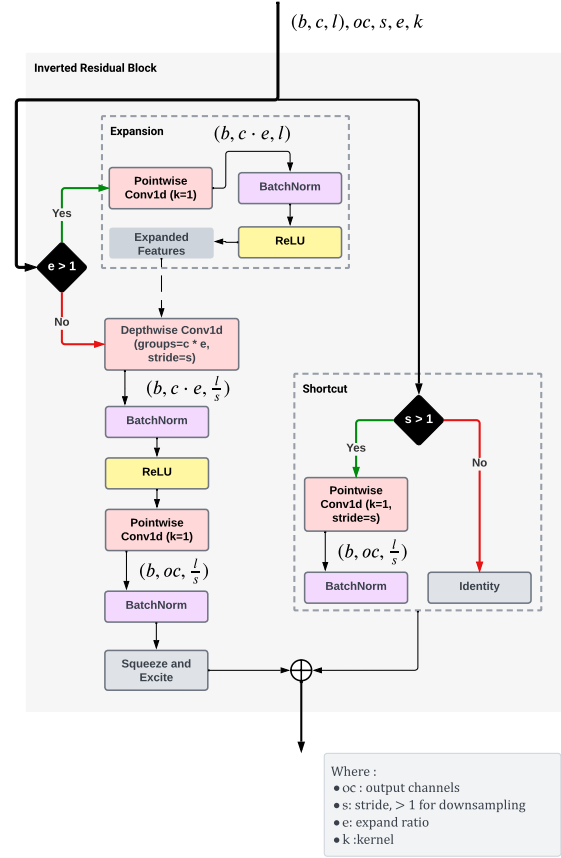


Figure 6. **Inverted residual block** structure enables optimal feature extraction while reducing the complexity of traditional basic blocks.

C. Additional material to Sec. 4

C.1. LME models to evaluate models per zone

In Fig. 8 we build 4 Linear Mixed Effects models to evaluate the prediction errors of all models across protocol zones and on all sampling rates. Results confirm observations from Tab. 1 where DenseNet models perform poorly on the Baseline zone compared to IKrNet models.

Table 3. **Generepol-Holter and Generepol-Holter_{HP} datasets.** The amount of raw signal (at 500 Hz) is obtained by dividing each partition value by 3, except in Holdout, where it is divided by 8. F_s represents the sampling rates to augment the datasets.

Partition details			Generepol-Holter _{HP}		Generepol-Holter	
Name	N°Patients	F_s (Hz)	Sot- Raw / Augmented	Sot+ Raw / Augmented	Sot- Raw / Augmented	Sot+ Raw / Augmented
Training	629	180, 250, 500	2739/8217	3576/10728	30736/92207	30808/92424
Validation	90	180, 250, 500	413/1239	539/1617	4698/14092	4739/14216
Evaluation	90	180, 250, 500	339/1017	451/1353	4898/14692	3287/9860
Holdout	90	150, 180, 215, 250, 300, 350, 425, 500	396/3168	502/4016	4673/37384	4710/37680

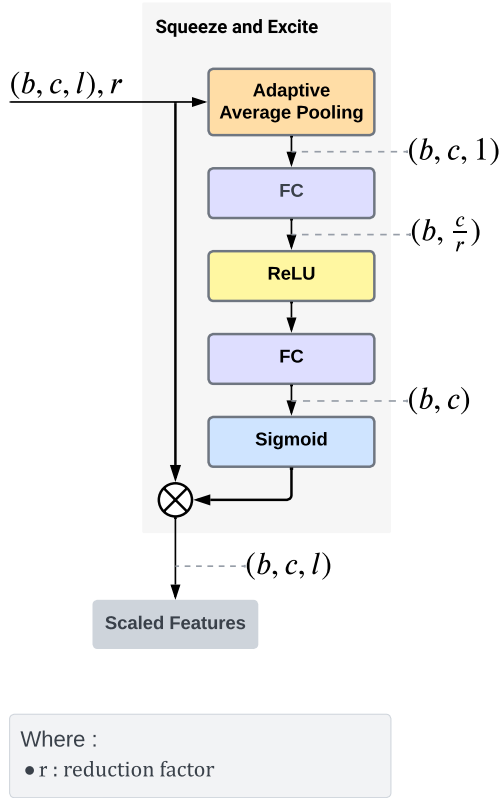


Figure 7. **Squeeze-and-Excite Block.** This structure is designed to enhance model performance by dynamically re-calibrating feature channels and optimizing feature extraction and prioritization within the network.

IKrNet₊ is generally more robust when considering all zones.

C.2. Ablation study

We conducted an ablation study in which we evaluated different configurations of the IKrNet architecture in the pro-

tol zones. Specifically, we evaluated the number of resolution branches in the CNN backbone, the type of convolution block in each branch (Inverted Residual Block Fig. 6, or basic block), the output size of each branch, the number of BiLSTM layers and the use of the residual skip connection leading to the combination of spatial and temporal features.

Each configuration is trained on the Generepol-Holter dataset with data augmentation. Results of the best-performing configurations are presented in Table 4; the first row corresponds to the final architecture structure that yields the best performances.

We observed the model’s performance drop when reducing the number of resolution branches in the CNN backbone to capture spatial features. This suggests that a less diverse set of receptive fields diminishes the model’s ability to detect abrupt signal changes, as it limits the range of spatial patterns the model can effectively analyze. Moreover, we observed that increasing the output size of each branch reduces the accuracy and leads to overfitting. BiLSTM layers also play a critical role; when removed from the architecture, we notice a significant drop in accuracy, particularly in the Sot + class in the St-Dg₊ and St+Dg₊ zones. This highlights the importance of capturing temporal dependencies, as BiLSTM layers enable the model to analyze the progression of ECG signals over time. Without these layers, the architecture struggles between transient heart rate induced changes and sustained drug-induced patterns, leading to a higher misclassification rate. The skip connection used to combine spatial and temporal features is also crucial for maintaining robust model performance, we noticed a drop of the accuracy by 7% in the St+Dg₊ zone, where stress coincides with the drug’s effects. The skip connection ensures that critical spatial features, such as subtle morphological changes in the ECG waveform, are retained and integrated with temporal patterns captured by the BiLSTM layers. Without this connection, the model struggles to effectively link the localized spatial features with their

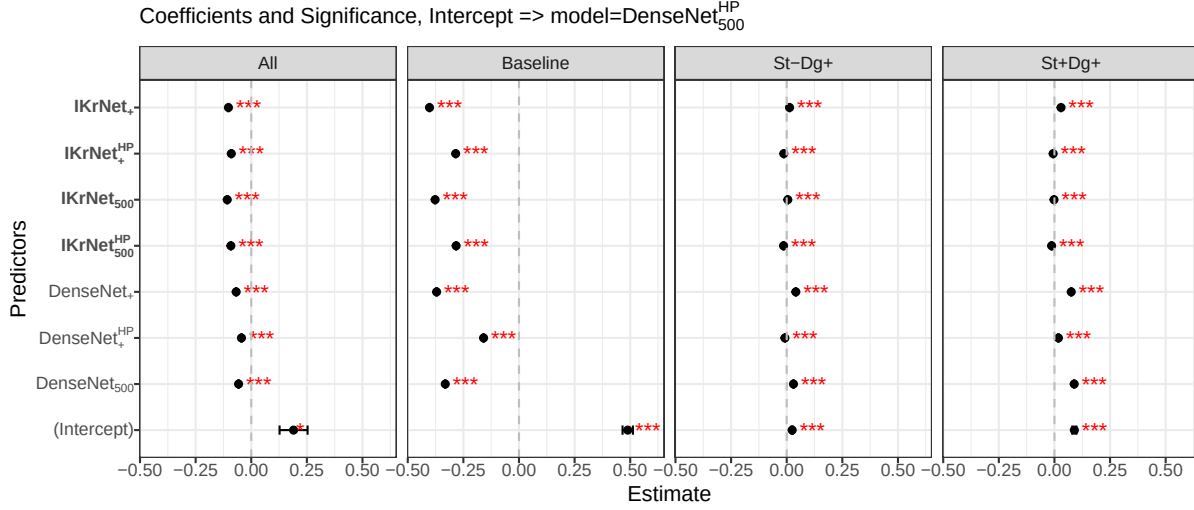


Figure 8. Linear Mixed Effects models per zone on all sampling rates to evaluate model predictions errors. In *all zones* plot (first facet) we used model: $predictions_error \text{ model} + (1|sampling_rate) + (1|zone) + (1|patient_id) + (1|ecg_id)$. For each zone individually, we used the model: $predictions_error \text{ model} + (1|sampling_rate) + (1|patient_id) + (1|ecg_id)$.

temporal progression.

C.3. Per sampling-rates results

Tables 5 and 6 present the numerical results of the eight models trained on the Generepol-Holter_{HP} dataset; Tables 7 and 8 summarize the results of the eight models trained on the Generepol-Holter datasets. In Fig. 9 we tested all eight models on a single patient similar to Fig. 5. Finally, Tabs. 9 to 12 present the results per patient and averaged over all the patients.

In this set of experiments, we evaluated models' performance under varying ECG sampling rates without distinguishing results by protocol zones as in Sec. 4. Overall, the IKrNet architecture consistently outperforms the DenseNet architecture, regardless of the training or testing set.

Similar to the robustness setting of the protocol zones, we assessed the models' robustness by defining subgroups based on ECGs at different sampling rates: $acc_{func} \stackrel{\text{def}}{=} func\{acc(150), acc(180), acc(215), acc(250), acc(300), acc(350), acc(425), acc(500)\}$, with $func \in \{\min, \max\}$. Below, we indicate the accuracy parity difference (APD) results for this setting. Results indicate that IKrNet+, trained in both datasets, is the most robust model among all frequency sampling rates.

APD_{↓0} of the models tested on Generepol-Holter_{HP}

- DenseNet₅₀₀^{HP}: 9.69 (Max accuracy: 95.99)
- DenseNet₊^{HP}: 5.57 (Max accuracy: 95.99)
- DenseNet₅₀₀: 2.33 (Max accuracy: 94.65)
- DenseNet₊: **0.67** (Max accuracy: 94.88)

- IKrNet₅₀₀^{HP}: 3.45 (Max accuracy: 96.88)
- IKrNet₊^{HP}: 1.45 (Max accuracy: 96.55)
- IKrNet₅₀₀: 1.89 (Max accuracy: 96.99)
- IKrNet₊: **0.9** (Max accuracy: 98)

APD_{↓0} of the models tested on Generepol-Holter


- DenseNet₅₀₀^{HP}: 16.51 (Max accuracy: 79.47)
- DenseNet₊^{HP}: 17.84 (Max accuracy: 84.80)
- DenseNet₅₀₀: 2.57 (Max accuracy: 89.98)
- DenseNet₊: 2.34 (Max accuracy: 91.92)
- IKrNet₅₀₀^{HP}: 7.14 (Max accuracy: 88.51)
- IKrNet₊^{HP}: 2.27 (Max accuracy: 88.19)
- IKrNet₅₀₀: **4.07** (Max accuracy: **94.37**)
- IKrNet₊: **1.54** (Max accuracy: 94)

Table 4. **Ablation study results:** different configurations are evaluated on the protocol zones. Each configuration is trained on Generepol-Holter with data-augmentation


Configurations					Accuracy		
N° resolution branches	Type of conv block	Output size per resolution branch	N° BiLSTM layers	Use residual skip link	Baseline	St-Dg ₊	St ₊ Dg ₊
4	Inverted residual	256	2	✓	92.22	98.92	95.35
4	Inverted residual	256	2	✗	92.45	91.01	88.23
4	Inverted residual	256	1	✓	91.92	96.34	92.67
4	Inverted residual	256	0	✓	97.89	85.90	82.04
2	Inverted residual	256	2	✓	90.65	92.93	88.67
1	Inverted residual	256	2	✓	89.31	89.94	85.34
4	Inverted residual	512	2	✓	92.06	96.68	94.19
4	Basic	256	2	✓	90.25	94.12	92.57

Table 5. **Comparison between DenseNet₅₀₀^{HP} and IKrNet₅₀₀^{HP}.** The models are trained on Generepol-Holter_{HP} with ECGs at 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline.


F_s	DenseNet ₅₀₀ ^{HP}				IKrNet ₅₀₀ ^{HP}			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	86.08	80.75	98.61	88.79	93.43	90.49	98.61	<u>94.38</u>
180	90.98	86.74	99.00	92.47	94.88	93.35	97.81	<u>95.53</u>
215	93.99	91.33	98.61	<u>94.83</u>	93.10	90.74	97.61	94.05
250	94.88	94.36	96.61	95.47	96.55	96.27	97.61	<u>96.93</u>
300	95.32	95.63	96.02	95.83	96.55	96.27	97.61	<u>96.93</u>
350	95.99	96.05	96.81	96.43	96.88	96.65	97.81	<u>97.23</u>
425	95.43	95.64	96.22	95.93	96.88	96.65	97.81	<u>97.23</u>
500	95.77	96.03	96.41	96.22	96.66	96.64	97.41	<u>97.02</u>
Test on Generepol-Holter dataset								
150	62.96	57.58	99.23	72.88	81.37	73.90	97.16	<u>83.95</u>
180	66.53	60.06	99.28	74.84	85.24	78.70	96.75	<u>86.80</u>
215	70.68	63.28	98.88	77.17	84.10	77.35	96.57	<u>85.90</u>
250	75.95	68.01	98.26	80.38	87.85	82.42	96.30	<u>88.82</u>
300	78.80	71.04	97.42	82.17	88.26	83.07	96.19	<u>89.15</u>
350	79.47	71.74	97.42	82.63	88.34	83.16	96.21	<u>89.21</u>
425	79.28	71.50	97.58	82.53	88.38	83.22	96.21	<u>89.25</u>
500	79.66	71.98	97.31	82.75	88.51	83.54	95.99	<u>89.34</u>

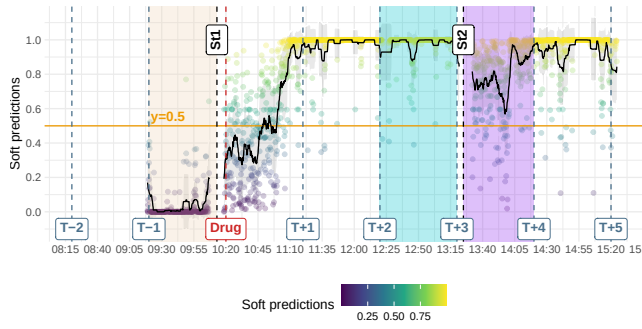


(a) DenseNet₅₀₀ predictions.

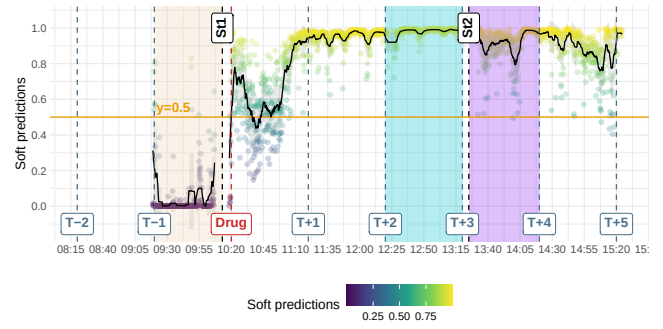


(e) IKrNet₅₀₀ predictions.

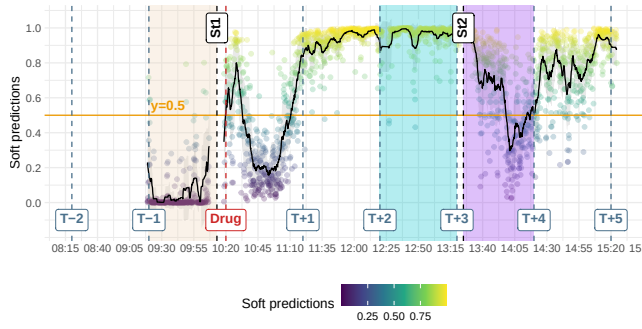


(b) DenseNet₊ predictions.

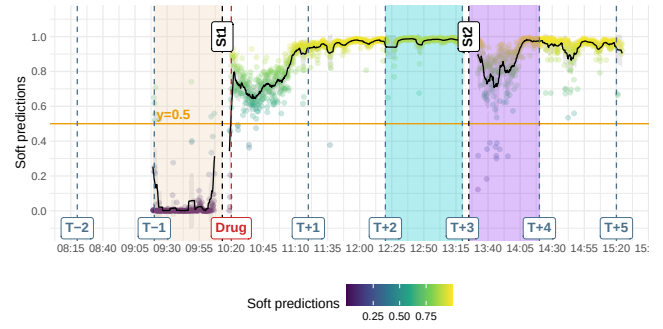


(f) IKrNet₊ predictions.

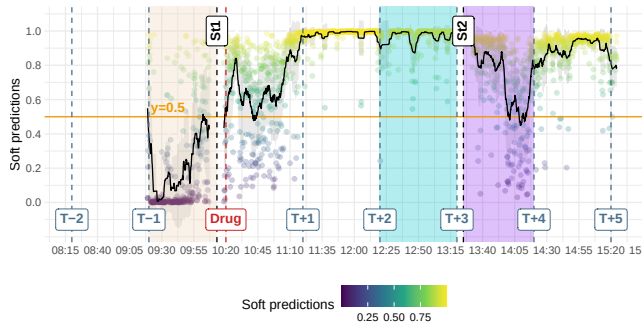


(c) DenseNet₅₀₀^{HP} predictions.

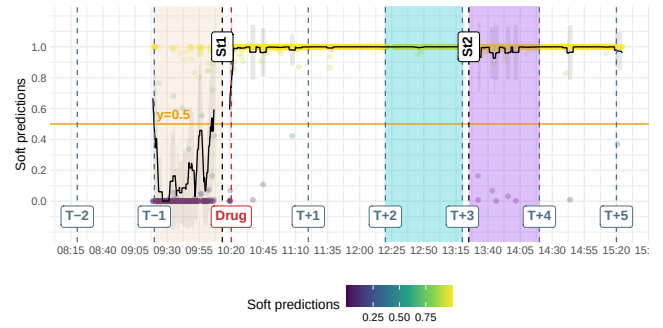

(g) IKrNet₅₀₀^{HP} predictions.

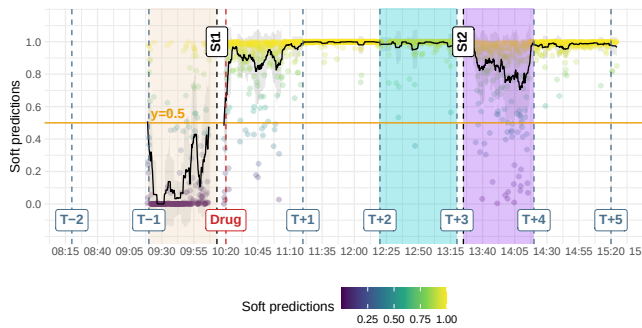

(d) DenseNet₊^{HP} predictions.

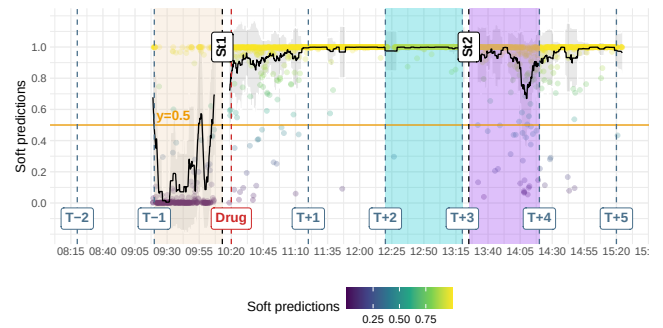

(h) IKrNet₊^{HP} predictions.

Figure 9. All models predictions on a single patient.

Table 6. **Comparison between DenseNet₊^{HP} and IKrNet₊^{HP}**. The models are trained on Generepol-Holter_{HP} augmented with ECGs at 180, 215, and 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline.

F_s	DenseNet ₊ ^{HP}				IKrNet ₊ ^{HP}			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	90.31	85.71	99.20	91.97	95.10	95.80	95.42	<u>95.61</u>
180	95.10	93.87	97.61	95.70	95.77	95.85	96.61	<u>96.23</u>
215	94.77	93.17	97.81	95.43	95.10	94.55	96.81	<u>95.67</u>
250	94.99	94.89	96.22	95.55	95.99	95.87	97.01	<u>96.44</u>
300	94.77	95.23	95.42	95.32	96.21	96.43	96.81	<u>96.62</u>
350	95.88	95.68	97.01	96.34	96.55	96.45	97.41	<u>96.93</u>
425	95.99	95.87	97.01	96.44	96.55	96.63	97.21	<u>96.92</u>
500	95.88	95.86	96.81	96.33	96.55	96.63	97.21	<u>96.92</u>
Test on Generepol-Holter dataset								
150	66.96	60.40	99.05	75.04	85.97	80.36	95.32	<u>87.20</u>
180	78.92	71.05	97.82	82.31	86.88	81.28	95.92	<u>88.00</u>
215	80.77	73.13	97.44	83.56	85.92	79.91	96.08	<u>87.25</u>
250	83.62	76.52	97.15	85.60	87.47	82.02	96.08	<u>88.50</u>
300	84.41	77.59	96.90	86.17	87.91	82.58	96.18	<u>88.86</u>
350	84.80	77.97	97.11	86.50	88.12	82.79	96.35	<u>89.05</u>
425	84.73	77.84	97.21	86.46	88.19	82.85	96.41	<u>89.12</u>
500	84.77	77.87	97.26	86.49	88.13	82.77	96.38	<u>89.06</u>

Table 7. **Comparison between DenseNet₅₀₀ and IKrNet₅₀₀**. The models are trained on Generepol-Holter_{HP} with ECGs at 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline.

F_s	DenseNet ₅₀₀				IKrNet ₅₀₀			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	92.32	94.82	91.24	92.99	95.10	95.80	95.42	<u>95.61</u>
180	92.76	94.14	92.83	93.48	96.33	98.15	95.22	<u>96.66</u>
215	93.54	94.40	94.02	94.21	95.43	96.38	95.42	<u>95.90</u>
250	94.32	95.93	93.82	94.86	96.99	98.77	95.82	<u>97.27</u>
300	94.21	96.30	93.23	94.74	96.21	98.35	94.82	<u>96.55</u>
350	94.54	96.70	93.43	95.04	96.88	98.97	95.42	<u>97.16</u>
425	94.65	96.90	93.43	95.13	96.77	98.96	95.22	<u>97.06</u>
500	94.43	96.69	93.23	94.93	96.33	98.96	94.42	<u>96.64</u>
Test on Generepol-Holter dataset								
150	87.41	85.49	90.21	87.79	90.30	85.41	97.27	<u>90.96</u>
180	87.69	85.05	91.54	88.18	92.17	88.42	97.10	<u>92.56</u>
215	88.18	85.43	92.15	88.66	90.74	86.02	97.34	<u>91.33</u>
250	88.93	86.71	92.01	89.28	93.51	90.65	97.07	<u>93.75</u>
300	89.51	87.71	91.97	89.79	93.76	91.04	97.12	<u>93.98</u>
350	89.84	88.17	92.10	90.09	93.87	91.19	97.16	<u>94.08</u>
425	89.94	88.36	92.06	90.17	94.10	91.71	97.00	<u>94.28</u>
500	89.98	88.68	91.73	90.18	94.37	92.69	96.38	<u>94.50</u>

Table 8. **Comparison between DenseNet₊ and IKrNet₊.** The models are trained on Generepol-Holter_{HP} augmented with ECGs at 180, 215, and 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline.

F_s	DenseNet ₊				IKrNet ₊			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	94.21	96.30	93.23	94.74	97.10	97.98	96.81	<u>97.39</u>
180	94.54	96.13	94.02	95.07	97.66	98.20	97.61	<u>97.90</u>
215	94.88	96.72	94.02	95.35	96.44	95.90	97.81	<u>96.84</u>
250	94.32	96.88	92.83	94.81	98.00	98.59	97.81	<u>98.20</u>
300	94.88	97.50	93.23	95.32	97.77	98.59	97.41	<u>98.00</u>
350	94.77	97.30	93.23	95.22	97.88	98.79	97.41	<u>98.09</u>
425	94.77	97.30	93.23	95.22	97.88	98.79	97.41	<u>98.09</u>
500	94.77	97.49	93.03	95.21	97.88	98.79	97.41	<u>98.09</u>
Test on Generepol-Holter dataset								
150	89.58	87.19	92.87	89.94	92.46	89.13	96.77	<u>92.79</u>
180	90.40	87.59	94.21	90.78	93.03	89.68	97.29	<u>93.33</u>
215	90.51	87.99	93.88	90.84	92.55	88.73	97.53	<u>92.92</u>
250	91.28	89.40	93.71	91.50	93.69	90.80	97.27	<u>93.93</u>
300	91.81	90.19	93.87	91.99	93.91	91.23	97.20	<u>94.12</u>
350	91.77	90.18	93.79	91.95	93.98	91.31	97.24	<u>94.18</u>
425	91.92	90.28	94.00	92.10	93.97	91.29	97.24	<u>94.17</u>
500	91.90	90.37	93.84	92.07	94.00	91.42	97.14	<u>94.19</u>

Table 9. **Per patient comparison between DenseNet₅₀₀^{HP} and IKrNet₅₀₀^{HP}.** The models are trained on Generepol-Holter_{HP} with ECGs at 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline. Borderline cases are omitted, i.e., when the models predict only negative samples for a patient or when the patient contains only true negative samples.

F_s	DenseNet ₅₀₀ ^{HP}				IKrNet ₅₀₀ ^{HP}			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	86.2 ± 20.9	81.5 ± 26.4	98.8 ± 4.9	91.5 ± 10.5	93.9 ± 13.5	91.7 ± 17.3	98.8 ± 4.5	<u>95.6</u> ± 7.8
180	91.4 ± 14.6	86.9 ± 23.3	99.1 ± 4.5	94.3 ± 9.0	94.8 ± 13.1	93.7 ± 16.4	98.0 ± 6.5	<u>96.4</u> ± 7.2
215	94.4 ± 12.7	92.1 ± 18.4	98.8 ± 5.1	<u>95.5</u> ± 9.6	92.4 ± 17.9	90.8 ± 21.9	97.9 ± 6.5	95.3 ± 9.5
250	95.5 ± 10.1	95.3 ± 13.4	97.0 ± 8.5	96.1 ± 7.5	96.5 ± 11.7	96.9 ± 12.4	97.9 ± 6.5	<u>97.5</u> ± 6.3
300	95.8 ± 10.4	96.3 ± 12.9	96.4 ± 10.0	96.3 ± 8.4	96.5 ± 11.8	96.9 ± 12.4	97.8 ± 6.7	<u>97.5</u> ± 6.3
350	96.3 ± 9.7	96.5 ± 12.8	97.2 ± 8.2	96.9 ± 7.3	96.8 ± 11.7	97.2 ± 12.3	98.1 ± 6.3	<u>97.8</u> ± 6.2
425	95.9 ± 10.4	96.2 ± 13.0	96.6 ± 9.5	96.4 ± 8.4	96.8 ± 11.7	97.2 ± 12.3	98.1 ± 6.4	<u>97.8</u> ± 6.2
500	96.1 ± 10.2	96.6 ± 12.8	96.8 ± 9.3	96.7 ± 8.1	96.6 ± 11.7	97.2 ± 12.3	97.7 ± 6.8	<u>97.6</u> ± 6.3
Test on Generepol-Holter dataset								
150	63.8 ± 12.4	59.6 ± 10.8	99.2 ± 1.7	73.9 ± 8.3	81.6 ± 9.0	75.6 ± 10.4	97.1 ± 3.3	<u>84.6</u> ± 6.8
180	67.3 ± 12.6	62.3 ± 11.4	99.3 ± 1.4	75.9 ± 8.5	85.4 ± 8.4	80.2 ± 9.7	96.8 ± 4.3	<u>87.3</u> ± 6.3
215	71.3 ± 12.3	65.5 ± 11.8	98.8 ± 2.1	78.2 ± 8.5	84.1 ± 9.5	79.1 ± 10.8	96.5 ± 4.3	<u>86.5</u> ± 6.8
250	76.3 ± 11.3	70.1 ± 11.6	98.2 ± 2.8	81.2 ± 7.9	87.9 ± 8.1	83.9 ± 9.7	96.3 ± 4.8	<u>89.2</u> ± 6.3
300	79.0 ± 11.1	73.1 ± 11.7	97.4 ± 3.6	82.9 ± 7.8	88.3 ± 7.9	84.5 ± 9.4	96.1 ± 4.9	<u>89.5</u> ± 6.1
350	79.7 ± 11.2	73.8 ± 11.8	97.4 ± 3.7	83.4 ± 7.9	88.4 ± 7.9	84.6 ± 9.5	96.2 ± 4.7	<u>89.6</u> ± 6.1
425	79.5 ± 10.9	73.6 ± 11.6	97.5 ± 3.5	83.3 ± 7.8	88.4 ± 7.9	84.6 ± 9.5	96.2 ± 4.8	<u>89.6</u> ± 6.1
500	79.9 ± 10.8	74.0 ± 11.5	97.3 ± 3.9	83.5 ± 7.7	88.5 ± 7.9	84.9 ± 9.4	95.9 ± 4.9	<u>89.7</u> ± 6.1

Table 10. **Per patient comparison between DenseNet₊^{HP} and IKrNet₊^{HP}**. The models are trained on Generepol-Holter_{HP} augmented with ECGs at 180, 215 and 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline. Borderline cases are omitted, i.e., when the models predict only negative samples for a patient or when the patient contains only true negative samples.

F_s	DenseNet ₊ ^{HP}				IKrNet ₊ ^{HP}			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	90.8 ± 17.0	87.7 ± 20.4	99.3 ± 3.7	93.8 ± 8.8	95.6 ± 9.1	97.2 ± 8.4	95.7 ± 9.6	<u>95.9 ± 7.2</u>
180	95.3 ± 11.3	94.2 ± 16.6	97.8 ± 6.5	96.5 ± 7.5	96.3 ± 8.4	97.2 ± 8.4	97.0 ± 8.1	<u>96.6 ± 6.6</u>
215	94.8 ± 13.3	93.5 ± 17.5	98.1 ± 5.7	<u>96.1 ± 8.5</u>	95.0 ± 11.2	94.4 ± 16.9	97.1 ± 8.1	96.0 ± 8.3
250	95.1 ± 12.5	95.8 ± 13.1	95.7 ± 12.2	96.3 ± 7.0	96.5 ± 8.5	97.3 ± 8.3	97.3 ± 7.7	<u>96.9 ± 6.4</u>
300	95.1 ± 11.2	96.0 ± 13.1	95.6 ± 10.7	95.7 ± 8.8	96.6 ± 8.6	97.7 ± 8.1	97.0 ± 8.2	<u>96.9 ± 6.6</u>
350	95.9 ± 12.4	96.5 ± 12.8	97.3 ± 7.0	97.0 ± 6.7	96.6 ± 9.4	96.9 ± 12.3	97.7 ± 7.6	<u>97.3 ± 6.4</u>
425	96.0 ± 12.4	96.6 ± 12.8	97.3 ± 7.0	97.1 ± 6.7	96.6 ± 9.5	97.0 ± 12.3	97.5 ± 7.8	<u>97.3 ± 6.4</u>
500	95.9 ± 12.4	96.6 ± 12.8	97.1 ± 7.3	97.0 ± 6.8	96.6 ± 9.5	97.0 ± 12.3	97.5 ± 7.8	<u>97.3 ± 6.4</u>
Test on Generepol-Holter dataset								
150	67.9 ± 12.1	62.6 ± 10.8	99.0 ± 2.2	76.1 ± 8.3	86.1 ± 8.2	82.0 ± 10.2	95.2 ± 6.6	<u>87.6 ± 6.7</u>
180	79.3 ± 10.3	73.0 ± 11.1	97.8 ± 3.4	83.1 ± 7.5	87.0 ± 8.0	82.9 ± 10.0	95.8 ± 6.2	<u>88.4 ± 6.5</u>
215	80.9 ± 10.1	74.9 ± 11.0	97.5 ± 4.0	84.2 ± 7.3	86.0 ± 8.7	81.6 ± 10.3	96.0 ± 6.3	<u>87.6 ± 6.8</u>
250	83.8 ± 9.2	78.1 ± 10.4	97.2 ± 4.2	86.1 ± 6.9	87.6 ± 8.0	83.5 ± 9.9	96.0 ± 6.4	<u>88.8 ± 6.5</u>
300	84.6 ± 9.2	79.2 ± 10.4	96.9 ± 4.0	86.7 ± 6.9	88.0 ± 7.8	84.1 ± 9.7	96.1 ± 6.2	<u>89.2 ± 6.4</u>
350	84.9 ± 9.3	79.5 ± 10.4	97.1 ± 4.0	87.0 ± 6.9	88.2 ± 7.9	84.3 ± 9.7	96.3 ± 6.0	<u>89.4 ± 6.4</u>
425	84.8 ± 9.3	79.4 ± 10.4	97.2 ± 3.9	87.0 ± 6.9	88.3 ± 7.8	84.3 ± 9.7	96.3 ± 5.9	<u>89.4 ± 6.4</u>
500	84.9 ± 9.3	79.4 ± 10.3	97.3 ± 3.7	87.0 ± 6.9	88.2 ± 7.9	84.3 ± 9.8	96.3 ± 5.9	<u>89.4 ± 6.4</u>

Table 11. **Per patient comparison between DenseNet₅₀₀ and IKrNet₅₀₀**. The models are trained on Generepol-Holter with ECGs at 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline. Borderline cases are omitted, i.e., when the models predict only negative samples for a patient or when the patient contains only true negative samples.

F_s	DenseNet ₅₀₀				IKrNet ₅₀₀			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	92.8 ± 13.4	95.8 ± 13.5	91.5 ± 16.9	92.7 ± 12.9	95.4 ± 11.1	96.2 ± 12.4	95.9 ± 13.1	<u>95.7 ± 10.4</u>
180	93.2 ± 13.5	95.3 ± 13.5	93.3 ± 15.8	93.4 ± 12.5	96.8 ± 8.1	98.6 ± 5.3	95.8 ± 12.0	<u>96.6 ± 8.4</u>
215	93.4 ± 14.7	94.0 ± 18.5	94.7 ± 13.0	94.9 ± 10.1	95.3 ± 13.1	96.2 ± 15.0	96.1 ± 11.8	<u>96.4 ± 8.8</u>
250	94.5 ± 11.5	95.7 ± 15.4	94.3 ± 13.1	95.1 ± 9.7	97.0 ± 8.9	98.2 ± 10.5	96.4 ± 11.4	<u>97.3 ± 8.1</u>
300	94.4 ± 12.6	96.8 ± 12.2	93.7 ± 14.8	94.7 ± 11.2	96.5 ± 9.7	97.9 ± 10.8	95.5 ± 13.6	<u>96.4 ± 10.4</u>
350	94.8 ± 12.1	96.9 ± 12.6	94.0 ± 13.0	95.2 ± 10.0	97.4 ± 8.1	99.2 ± 4.3	96.0 ± 12.8	<u>96.9 ± 9.6</u>
425	94.9 ± 12.1	97.1 ± 12.5	94.0 ± 13.0	95.3 ± 10.0	97.3 ± 8.4	99.1 ± 4.6	95.8 ± 12.9	<u>96.8 ± 9.8</u>
500	94.6 ± 12.2	96.9 ± 12.6	93.7 ± 13.5	95.0 ± 10.4	96.9 ± 8.9	99.1 ± 4.6	95.1 ± 14.2	<u>96.3 ± 10.6</u>
Test on Generepol-Holter dataset								
150	87.6 ± 9.4	88.2 ± 12.0	90.0 ± 13.5	87.7 ± 10.7	90.5 ± 8.0	87.1 ± 10.4	97.2 ± 5.3	<u>91.4 ± 6.6</u>
180	87.9 ± 9.2	87.8 ± 12.1	91.4 ± 12.1	88.3 ± 9.6	92.3 ± 7.3	89.8 ± 9.3	97.0 ± 5.4	<u>92.9 ± 6.2</u>
215	88.4 ± 9.2	88.1 ± 12.0	92.0 ± 11.7	88.8 ± 9.3	90.7 ± 8.9	87.8 ± 10.9	97.3 ± 5.2	<u>91.8 ± 6.9</u>
250	89.1 ± 8.9	89.2 ± 11.6	91.8 ± 11.6	89.4 ± 9.0	93.6 ± 7.1	91.9 ± 8.8	97.0 ± 5.4	<u>94.0 ± 6.0</u>
300	89.7 ± 8.6	90.0 ± 11.2	91.8 ± 11.5	89.8 ± 8.9	93.9 ± 7.0	92.3 ± 8.7	97.0 ± 5.4	<u>94.2 ± 5.9</u>
350	90.0 ± 8.5	90.4 ± 10.9	92.0 ± 11.5	90.1 ± 8.8	94.0 ± 7.0	92.4 ± 8.8	97.1 ± 5.3	<u>94.3 ± 6.0</u>
425	90.1 ± 8.4	90.6 ± 10.9	91.9 ± 11.4	90.2 ± 8.7	94.2 ± 7.0	92.9 ± 8.6	96.9 ± 5.6	<u>94.5 ± 6.0</u>
500	90.1 ± 8.4	90.9 ± 10.9	91.6 ± 11.5	90.2 ± 8.7	94.5 ± 6.9	93.9 ± 8.4	96.3 ± 5.9	<u>94.7 ± 6.0</u>

Table 12. **Per patient comparison between DenseNet₊ and IKrNet₊.** The models are trained on Generepol-Holter augmented with ECGs at 180, 215 and 500 Hz. The best accuracy is highlighted in **bold**, and the best F1 score is highlighted with underline. Borderline cases are omitted, i.e., when the models predict only negative samples for a patient or when the patient contains only true negative samples.

F_s	DenseNet ₊				IKrNet ₊			
	Accuracy	Precision	Recall	F1 Score	Accuracy	Precision	Recall	F1 Score
Test on Generepol-Holter_{HP} dataset								
150	94.3 ± 11.7	96.4 ± 13.1	92.8 ± 15.9	94.9 ± 10.3	97.3 ± 7.1	98.6 ± 5.7	97.0 ± 8.7	<u>97.4</u> ± 6.0
180	94.6 ± 11.0	96.2 ± 13.1	93.6 ± 15.4	95.3 ± 9.8	97.9 ± 6.2	98.8 ± 5.6	97.8 ± 6.8	<u>98.1</u> ± 5.1
215	95.2 ± 11.5	96.9 ± 12.8	93.5 ± 16.3	95.4 ± 10.7	96.0 ± 12.0	95.6 ± 16.4	98.0 ± 5.8	<u>97.3</u> ± 7.9
250	94.6 ± 10.9	97.8 ± 8.0	92.4 ± 16.0	94.8 ± 10.1	98.2 ± 6.1	99.1 ± 5.2	98.0 ± 6.5	<u>98.3</u> ± 5.0
300	95.2 ± 11.5	98.2 ± 8.6	92.8 ± 16.5	95.2 ± 11.1	98.0 ± 6.4	99.0 ± 5.3	97.6 ± 7.0	<u>98.1</u> ± 5.3
350	95.0 ± 11.7	98.1 ± 8.7	92.7 ± 16.2	95.2 ± 10.9	98.1 ± 5.8	99.1 ± 4.8	97.6 ± 7.0	<u>98.2</u> ± 5.1
425	94.9 ± 11.5	98.1 ± 8.4	92.8 ± 16.3	95.2 ± 10.7	98.1 ± 5.8	99.1 ± 4.8	97.6 ± 7.0	<u>98.2</u> ± 5.1
500	94.9 ± 11.6	98.3 ± 8.4	92.6 ± 16.4	95.1 ± 10.8	98.1 ± 5.8	99.1 ± 4.8	97.6 ± 7.0	<u>98.2</u> ± 5.1
Test on Generepol-Holter dataset								
150	89.8 ± 8.8	89.2 ± 10.9	92.6 ± 12.2	89.8 ± 10.2	92.6 ± 7.3	90.5 ± 9.1	96.8 ± 6.5	<u>93.1</u> ± 6.3
180	90.6 ± 8.2	89.5 ± 10.5	94.0 ± 10.5	90.8 ± 8.6	93.2 ± 7.0	91.0 ± 9.0	97.3 ± 5.0	<u>93.7</u> ± 5.8
215	90.7 ± 8.1	89.8 ± 10.1	93.7 ± 10.5	90.8 ± 8.4	92.6 ± 7.7	90.2 ± 9.7	97.5 ± 4.3	<u>93.3</u> ± 6.1
250	91.4 ± 7.9	91.1 ± 9.8	93.5 ± 10.5	91.5 ± 8.2	93.8 ± 6.8	92.0 ± 8.7	97.3 ± 4.8	<u>94.2</u> ± 5.6
300	92.0 ± 7.6	91.8 ± 9.4	93.7 ± 10.1	92.0 ± 8.0	94.0 ± 6.6	92.4 ± 8.6	97.2 ± 4.8	<u>94.4</u> ± 5.6
350	91.9 ± 7.7	91.7 ± 9.4	93.6 ± 10.3	91.9 ± 8.2	94.1 ± 6.7	92.5 ± 8.5	97.2 ± 4.6	<u>94.5</u> ± 5.5
425	92.1 ± 7.6	91.8 ± 9.4	93.8 ± 10.3	92.1 ± 8.3	94.1 ± 6.7	92.5 ± 8.5	97.3 ± 4.7	<u>94.5</u> ± 5.6
500	92.0 ± 7.6	91.9 ± 9.4	93.6 ± 10.2	92.0 ± 8.2	94.1 ± 6.7	92.6 ± 8.5	97.1 ± 4.8	<u>94.5</u> ± 5.6