

CARDIOLOGENT: Multi-Agent Clinical Decision Support for Patient-Level Arrhythmia Assessment, Urgency, and Management

Sukju Oh¹, Moo-Yong Rhee², Jae-Sik Jang², Sukkyu Sun^{1*}

¹Department of Computer Science and Artificial Intelligence, Dongguk University, Seoul, 04620, Republic of Korea

²Department of Cardiology, Dongguk University Ilsan Hospital, Goyang, 10326, Republic of Korea
dhtjrwn119@dgu.ac.kr, mooyong.rhee@gmail.com, jsjang71@gmail.com, sukkyu.sun@dgu.ac.kr

Abstract

The same episode of atrial fibrillation is a minor finding in a healthy adult and grounds for anticoagulation in an elderly patient with hypertension: identical signal, opposite decision. Naming the rhythm is only the start; what determines a patient's outcome is the judgement that follows—what the arrhythmia is across the whole record, what it means for this patient, and what should be done about it. Recent work pairing large language models with the ECG stops short of this, reading one recording without assembling it into a patient-level finding; and agentic systems built around it either receive the arrhythmia a device has already detected or target a different diagnostic task, stopping before the decision this task requires. We formulate patient-level arrhythmia decision support as a task and present CARDIOLOGENT, a multi-agent system that spans it from detection to decision. An agent for each signal—a single ECG lead and the photoplethysmogram a wearable acquires—grounds its window reading in measured features rather than a bare label; the readings are assembled into the patient's rhythm profile and, with the patient's own data, reasoned against clinical guidelines retrieved for the case, with a critic checking each conclusion against the guideline it cites. We evaluate the clinical decision rather than the report, across integrated diagnosis, clinical significance, and urgency and management. CARDIOLOGENT scores highest on every axis, first on every patient-level task under both cardiologists and an at-scale LLM judge—whose agreement with the cardiologists (ICC 0.74, 0.66) matches theirs with each other (0.67). Because each conclusion traces to a cited guideline and is validated against expert cardiologists, it yields decisions a clinician can audit rather than act on blindly—a step toward systems that could be used in continuous monitoring. Code and data are available at <https://github.com/sukjuoh/Cardiologent>.

Introduction

Cardiac arrhythmias are a major cause of cardiovascular morbidity, contributing to stroke, heart failure, and sudden cardiac death (Fu 2015; Zeppenfeld et al. 2022). Atrial fibrillation, the most common arrhythmia, is asymptomatic in roughly 32% of patients (Frykman et al. 2001) and often goes undetected until a severe event occurs. Detecting such arrhythmias reliably requires monitoring a patient's cardiac signals continuously.

Deep learning has driven much of the progress on this problem. End-to-end networks classify arrhythmias from

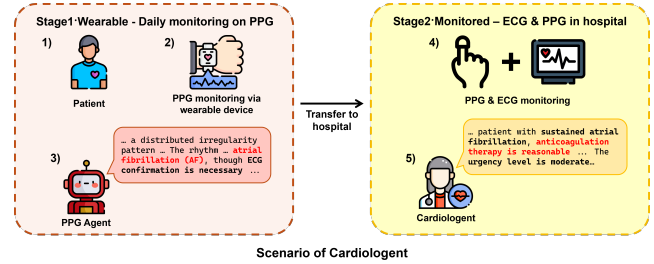


Figure 1: **Deployment scenario.** *Stage 1, wearable:* the pulse is monitored on its own, and the PPG agent reads it as far as the pulse allows—an irregular rhythm it calls atrial fibrillation, flagged as needing an ECG to confirm. *Stage 2, monitored:* the ECG is added, and CARDIOLOGENT reads it together with the pulse and the patient's own data to reach the decision—the therapy indicated, and how urgent. Excerpts are verbatim.

single- and 12-lead electrocardiogram (ECG) at cardiologist-level accuracy (Hannun et al. 2019; Ribeiro et al. 2020; Siontis et al. 2021; Li et al. 2025b), and comparable methods detect atrial fibrillation from wearable photoplethysmogram (PPG) at population scale (Perez et al. 2019; Guo et al. 2019). These models predict a class label but do not explain it: they give no account of which signal features support the prediction or why competing rhythms are excluded (Wang et al. 2026). Many arrhythmias are separated by narrow margins—atrial fibrillation from other irregular rhythms, supraventricular from ventricular ectopy—and the distinction turns on specific, localizable features of the signal. A clinician given only a label cannot see which feature the decision rested on, and so cannot confirm or overturn it; since these distinctions carry directly into management, an unverifiable label is of limited clinical use.

Recent work supplies the grounds a bare label lacks, on the discovery that large language models can be made to operate on time series at all (Zhou et al. 2023; Jin et al. 2024; Al-negheimish et al. 2024; Zhou and Yu 2024; Quinlan, Li, and Zhu 2026). Applied to the ECG, this produces models that pair the waveform with a language model and explain their reading rather than only labeling it, from a textual diagnosis to a full report grounded in the waveform (Liu et al. 2026;

*Corresponding authors.

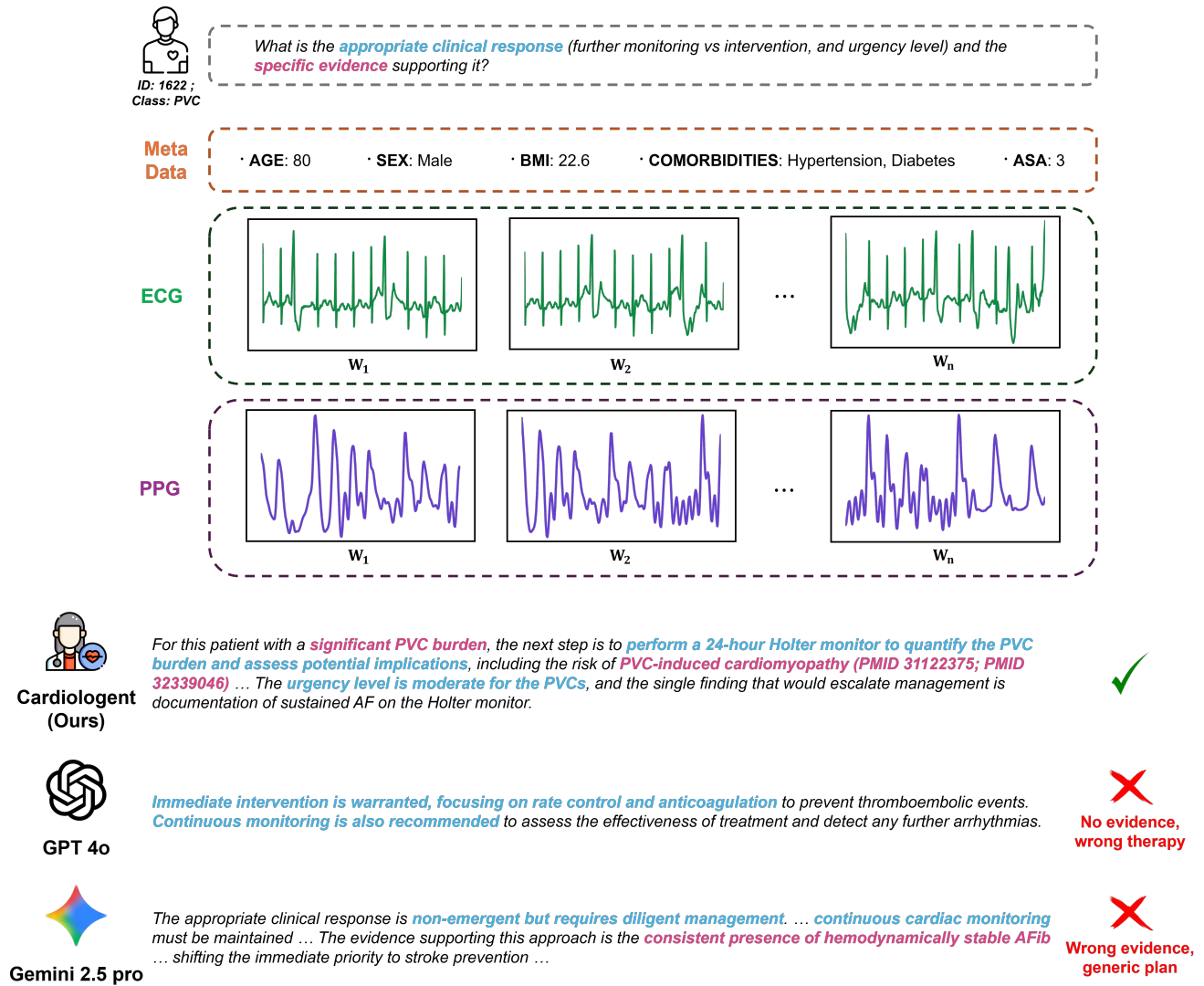


Figure 2: **Patient-level arrhythmia decision support on a monitored record.** An 80-year-old man with hypertension and diabetes, whose record holds a burden of premature ventricular contractions (PVC). Every system gets the same input: his metadata and the windowed lead II ECG and PPG. Blue marks the clinical response, pink the evidence for it. GPT-4o prescribes rate control and anticoagulation with no evidence at all, never naming the ectopy. Gemini 2.5 Pro names a rhythm the record does not show—atrial fibrillation—and builds on it a plan that would fit any patient. CARDIOLOGENT quantifies the finding it has and cites the guideline behind the recommendation. Excerpts are verbatim; ellipses mark omitted text.

Wan et al. 2025; Lan et al. 2026; Jin et al. 2026a,b; Chung et al. 2026; Li et al. 2025a). But they stop at the reading: each reports what the signal shows, and none assembles it into a finding about the patient.

A second line of work reaches further, to the document a physician actually reads. Zodiac assembles a patient-specific report from a monitoring study, synthesizing findings against clinical guidelines under cardiologist validation (Zhou et al. 2024); CardAIC-Agents ingests the raw signal and coordinates tools over it (Zhang et al. 2025). But neither spans the task we pose. Zodiac is handed the findings a recording device already computed and a strip already selected to depict the significant rhythm—detection, and the choice of what to

look at, have happened before the model is invoked—and its interpretation stops at what those findings mean, validated as a report: its accuracy against the patient’s data and its succinctness. CardAIC reads the signal itself, but its task is disease-level diagnosis scored by classification accuracy, not the rhythm-level reading of a monitored record assembled into a patient-level finding; the management plans it forms along the way are never scored.

Between naming a rhythm and knowing what to do about it lies most of the clinical act, and it is not a step that follows from the waveform. The same episode of atrial fibrillation is a minor finding in a healthy 40-year-old and grounds for anticoagulation in an elderly woman with hypertension (Joglar et al.

2024): identical signal, opposite decision. What separates them is the patient—age, sex, comorbidity, medication—and none of it is visible in the trace. Nor is the decision made on one window. Whether an abnormal beat is incidental or part of a sustained arrhythmia, and whether an episode occurred or happened once, are properties of the record rather than of any window in it, and they are what set the urgency. A patient-level judgement therefore needs both: the record assembled, and the patient’s own data alongside it. GEM’s authors observe the second directly: without clinical context, their model overstates the severity of its findings (Lan et al. 2026).

The signals such a judgement can be built on are constrained in turn. The 12-lead ECG demands electrode placement and a stationary patient, and is ill-suited to monitoring a record over hours; what monitors record instead is a small montage in which lead II is the rhythm lead—its axis aligns with that of the P wave, so atrial activity shows itself most clearly there—alongside the PPG, acquired continuously by clinical monitors and by every wearable on the market. The two are not always both in hand (Figure 1): the pulse runs continuously, on a ward monitor and on a wearable alike, while the ECG is attached only when something in the pulse calls for it. A system for this setting must therefore reach a reading from the pulse by itself, enough to say that an ECG is needed, and a fuller one once the ECG arrives. Hence an agent per signal rather than one reader over the pair, which is defined only when both are present.

We present **CARDIOLOGENT**, a multi-agent system that works on a monitored record—lead II ECG and PPG—at two levels. At the window level, an agent for each signal reads its own waveform through tools that measure it—a specialist fine-tuned on that signal, and analysts for the features the reading turns on—and grounds every claim in a measured quantity; the two readings are then reconciled into one verdict for the window. At the patient level, the windows are assembled into the patient’s rhythm profile—which arrhythmias appear and which dominates—and only then does the patient enter: retrieved guidelines, a reasoning stage judging what the profile means for this patient, what to do about it and how urgently, and a critic holding each conclusion to the guideline it cites. Figure 2 shows what this is worth on one patient: given the same record and the same metadata, GPT-4o prescribes anticoagulation for a rhythm the record does not hold, citing nothing; Gemini 2.5 Pro names that rhythm as its evidence and returns a plan that would read the same for any patient; **CARDIOLOGENT** quantifies the finding actually present and cites the guideline its recommendation comes from. Work on the PPG has concentrated almost entirely on atrial fibrillation (Pereira et al. 2020; Liu et al. 2022); our window level reads all seven classes from the pulse alone.

In summary, the main contributions of this work are as follows:

- We formulate patient-level arrhythmia decision support as a task: not the classification of a window alone, nor the interpretation of findings a device has already computed, but the whole span between them—what arrhythmia a record holds, what it means for this patient, and what should be done about it.

- We present **CARDIOLOGENT**, the first system to span this task, as we define it, from detection to decision. Where prior work receives its findings from a recording device, **CARDIOLOGENT** assembles the patient’s profile from its own window-level evidence; at the window level it also runs on the PPG alone, the signal a wearable acquires.
- We evaluate the clinical decision itself, not the fluency of a report: cardiologists score **CARDIOLOGENT** against frontier baselines on what the system concludes and would do, the level at which such a system must ultimately be judged.

Dataset

The task we formulate calls for data of a particular kind: two signals recorded at once, continuously, over a monitoring procedure, with the rhythm in them annotated. VitalDB is an open database of intraoperative monitoring from 6,388 non-cardiac surgical cases, recording ECG and PPG concurrently alongside each patient’s clinical record but saying nothing about the rhythm they contain (Lee et al. 2022); the VitalDB Arrhythmia Database supplies that, annotating 734,528 seconds of that ECG from 482 of the patients with 661,894 heartbeats across four beat types and ten rhythm categories, read by at least two of five anesthesiologists at $\kappa = 0.930$ (Eun et al. 2026).

Of the 482 annotated patients, 426 have both signals, giving 166,026 windows. Because the labels are on the ECG, we transplant them onto the PPG, aligning each R peak with the pulse it drove by pulse arrival time. On the ECG a window carries one of nine classes—*normal*, atrial fibrillation, multifocal atrial tachycardia, supraventricular tachycardia, ventricular tachycardia, PVC, PAC, AV block, and sinus node dysfunction. On the PPG these collapse to a coarser seven, since the pulse cannot resolve some distinctions the ECG can; the merged set and its rationale are given with the PPG agent (Method).

The split is by patient, so no case contributes windows to more than one side of it: 302 patients and 126,959 windows train the specialists, 59 patients and 20,194 windows are held out for validation, and the remaining 65 patients and 18,873 windows are the test set, on which every result reported here is computed.

Problem Formulation

A monitored record consists of two synchronized signals—a lead II ECG and a PPG—acquired continuously through a procedure, together with the patient’s own data: age, sex, comorbidity, and the rest of the case record. The record is divided into short windows, each spanning a fixed number of beats (around ten seconds; Dataset).

Window level. Each window is assigned one arrhythmia class—nine when both the ECG and the PPG are present, and a coarser seven in the PPG-only setting (Dataset)—together with the span of signal that is abnormal, which the class alone does not give.

Patient level. Given every window of a record and the patient’s own data, the task is to say what the record means for that patient, in three parts. Which arrhythmias appear, where, and which is the dominant rhythm—a property of the

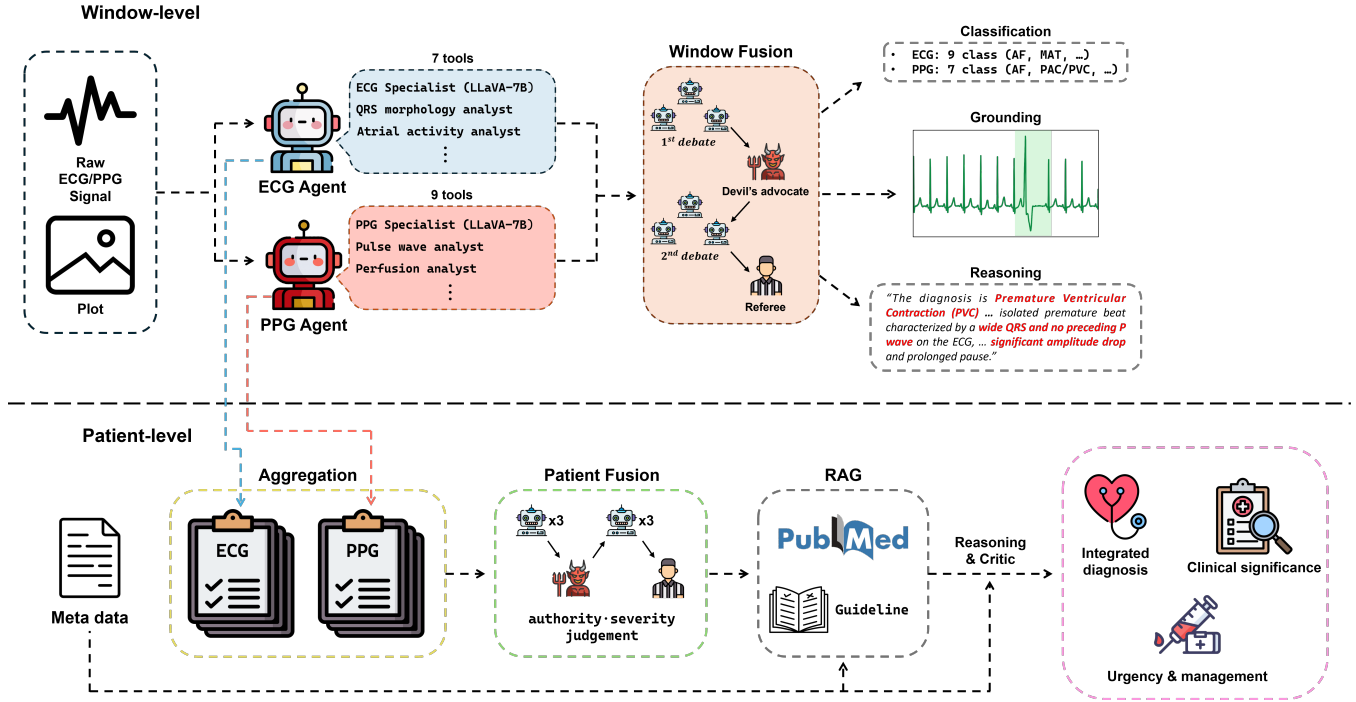


Figure 3: **CARDIOLOGENT**. *Window level*: an agent per signal reads the waveform through its tools—seven on the ECG, nine on the PPG—and returns a grounded class, and window fusion reconciles the ECG and PPG readings by debate into one verdict. *Patient level*: the verdicts are aggregated into a rhythm profile, reconciled across signals by the same debate, and only then read against the patient’s metadata through guideline retrieval, a reasoning stage, and a critic, yielding the integrated diagnosis, its significance, and the urgency and management.

record and of no window in it. What that profile signifies for this patient, and which single finding is the one to act on. And what should be done: further monitoring or intervention, at what urgency, on what evidence.

Method

CARDIOLOGENT works at the two levels defined above, shown in Figure 3. At the window level, an agent per signal reads the ECG or the PPG and the two readings are reconciled into one verdict per window. At the patient level, each agent’s window readings are assembled across the record into a rhythm profile, the two signals’ profiles are reconciled, and only afterwards is the profile read against the patient’s own data and the clinical guidelines. Both levels reconcile the two signals by the same debate, described under window fusion below; we take the pipeline in order.

ECG and PPG Agents

Each signal is read by its own agent, one for the ECG and one for the PPG, identical in construction. At the centre of each is a *specialist*: a LLaVA-7B multimodal model (Liu et al. 2023a) fine-tuned to read the raw waveform and its plot. A language model reads a raw physiological trace only weakly on its own (Zhou et al. 2023; Jin et al. 2024; Alnegheimish et al. 2024; Zhou and Yu 2024; Quinlan, Li, and Zhu 2026), so we adapt it to the signal, as recent ECG language models

do (Lan et al. 2026; Liu et al. 2026; Jin et al. 2026b). Around it sit deterministic tools that measure the quantities a rhythm reading turns on—rate, R–R regularity, and P-wave and QRS descriptors on the ECG, pulse-wave morphology and perfusion on the PPG—so every number the agent reports is a measured value rather than one read off the trace. Each agent reaches its reading through its own internal deliberation (below) and returns one reading: a class, the abnormal span, and a rationale citing the measurements behind it. Because each agent is complete on its own, the system matches whatever a deployment provides: the PPG agent alone in a wearable, both together where the ECG is monitored too. Running alone costs resolution—the pulse cannot separate ventricular from atrial ectopy, nor MAT from VT (Page et al. 2016; Al-Khatib et al. 2018; Han et al. 2020), so PVC merges with PAC and VT with MAT, leaving the coarser seven. Architecture and prompts are in the supplement.

Window Fusion

The two readings are reconciled into one verdict by a *debate*, several instances arguing the reading out across rounds rather than one model refining its own answer, which tends to defend a confident error (Du et al. 2023; Liang et al. 2024). Because the agents use different label sets—nine classes on the ECG, seven on the PPG with ectopy and VT/MAT merged where the pulse cannot separate them—the readings are first checked for compatibility, telling one finding named two

ways from a real conflict, then argued through three roles:

- **Analysts** re-examine both readings, each along one axis—rate, atrial activity, and ventricular morphology—weighing each modality by how well its call is grounded in its own measured evidence rather than by a fixed precedence, since neither signal is uniformly the better witness.
- A **devil’s advocate** then argues against the emerging verdict, pressing its weakest points; this adversarial role is what separates debate from self-refinement, and keeps a confident but wrong reading from passing unchallenged (Kim, Kim, and Yoon 2024). Whereas that work uses the devil’s advocate to *score* generated text, here it contests a diagnostic decision and forces the analysts to re-examine the measured evidence behind it. The analysts answer it in a further round.
- A **referee** weighs the exchange and issues the verdict, so a correct minority reading can prevail on its evidence rather than be outvoted (Du et al. 2023). It grants authority feature by feature to the signal that can see the feature: the ECG’s QRS width and P wave—which the pulse cannot resolve—split a merged PPG class such as VT/MAT into VT or a supraventricular rhythm, while the PPG’s pulse regularity and perfusion carry the haemodynamic significance the ECG does not show.

The verdict is one nine-class label; the abnormal span is the union of the two. In the wearable scenario the PPG agent’s own verdict stands.

Patient-Level Assembly

- **Aggregation.** Each agent’s window readings are gathered, in time order, into a record-level profile per signal—which arrhythmias occur, over how much of the record, and which dominates—that no single window contains; the PPG profile is held to what the pulse can support, flagging rather than asserting an atrial rhythm.
- **Patient fusion.** The ECG and PPG profiles are reconciled by the same window-fusion debate, judged now on authority and severity rather than per-feature evidence: authority weighs how confident each profile is, so the ECG leads where it is sure and yields where it is not, and severity breaks an even exchange toward the more dangerous reading, so the finding a clinician must not miss survives.

Patient Reasoning

Only once the profile is fixed does the patient enter, and the remaining stages turn the finding into a decision.

- **Retrieval (RAG).** The profile and the metadata generate a few clinical queries about the significance, management, and risk of the finding—not its identity, which is already settled—and retrieve matching passages from two sources: a curated corpus of arrhythmia guideline statements and the PubMed literature. Each retrieved passage carries its citation, so the evidence it supplies arrives already referenced; the corpus and retrieval details are in the supplement.

- **Reasoning.** A debate over the metadata, the profile, and the retrieved passages produces the three patient-level outputs: the integrated diagnosis, its clinical significance, and the urgency with the management, each tied to a retrieved passage and citing the PubMed reference that passage carries.
- **Critic.** A final pass checks the assessment against those passages—that each claim is grounded in the profile and supported by a citation actually retrieved, not one invented—and revises it where it is not, before the outputs are returned.

The order is deliberate: the record is read and reconciled before the patient is seen, so that what the arrhythmia *is* rests on the waveform, and only what it *means* draws on who the patient is.

Experiments

We evaluate CARDIOLOGENT at two levels—the window level and the patient level—against one fixed set of baselines: five general-purpose vision–language models (GPT-4o (Hurst et al. 2024), Gemini 2.5 Pro (Comanici et al. 2025), Qwen2.5-VL-72B (?), Llama-3.2-90B-Vision, MedGemma-27B (Sellersgren et al. 2026)) given identical inputs and questions. The baselines are prompted directly at both levels, so the gap at each measures what the construction adds over a capable model given the same inputs. Evaluation is on 3,000 windows sampled from the held-out test patients, covering every test patient and preserving class proportions. Reasoning is scored by an LLM judge (Zheng et al. 2023; Gu et al. 2025)—DeepSeek-V4-Pro, outside the baseline set, so no system grades its own outputs, and one whose scores agree with blinded cardiologists as closely as the cardiologists agree with each other (patient level; Figure 4). The orchestration (analysts, devil’s advocate, referee, and patient-level stages) runs on Qwen2.5-72B, an openweight model, so the whole pipeline can be run on-premises without sending patient data to an external API; since one baseline is the vision variant of that same model at the same scale, the gap over it isolates the agent construction from the backbone. Rubrics and sampling are in the supplement.

Metrics. Classification is reported as macro-F1 and balanced accuracy, not inflated by the dominant *normal* class. Grounding uses IoU and, since a whole-window span scores a high IoU without localizing, a per-sample MCC (Matthews 1975; Chicco, Tötsch, and Jurman 2021) we treat as primary. Reasoning is judged on four axes, each 1–4, on every window probe:

- **Correctness** — is the diagnosis itself right (the rhythm matches ground truth)?
- **Faithfulness** — are the features it cites real, with measured values that match the signal rather than fabricated?
- **Validity** — do those features actually support the diagnosis, or are they irrelevant to it?
- **Completeness** — are the discriminators that separate this rhythm from its neighbors all addressed, or is the reasoning partial?

The window agent is drawn out through three probes:

- **T1** — diagnosis with supporting features.
- **T2** — differential, excluding each alternative.
- **T3** — the one change that would alter the diagnosis.

Each probe is scored by the judge against a per-window gold reference anchored to the ground-truth label and cross-checked by a separate model (supplement). The patient pipeline is scored on three axes, each two sub-axes (0–4, to 8):

- **P1** — integrated diagnosis: *correctness* (the dominant rhythm matches ground truth and secondary arrhythmias are named) and *faithfulness* (every claim traces back to the waveform).
- **P2** — clinical significance: *prioritization* (the single most dangerous finding is foregrounded) and *context* (age, comorbidity, and surgical setting are used to drive the risk).
- **P3** — urgency and management: *management* (a specific, guideline-concordant recommendation tied to the class and patient) and *urgency and safety* (the stated urgency is right and the plan is safe).

Per-grade rubrics are in the supplement.

CARDIOLOGENT at the Window Level

The window level wraps the two specialists in the full system—deterministic tools, the internal debate, and window fusion across the two signals—and reads out a per-window verdict. The specialists are also measured on their own signal in the supplement, where each, at 7B, leads every baseline on both classification and reasoning.

Classification. CARDIOLOGENT raises nine/seven-class macro-F1 to 0.449 (both) and 0.425 (PPG), roughly double the strongest baseline (Table 1). As a binary normal-vs-arrhythmia decision it reaches balanced accuracy above 0.9 in both settings, while most baselines collapse toward chance from one side, calling nearly all windows normal or nearly all abnormal. GPT-4o is the apparent exception at 0.83 on *both*, but its low macro-F1 shows this rests on the two frequent classes: it separates AF from normal, misses the rarer rhythms, and falls to 0.56 once the ECG is removed. Detection is the easier tier; the nine-class macro-F1 above it is where the gap is widest.

Reasoning. The window agent is probed on diagnosis (T1), differential (T2), and the one change that would alter the call (T3); CARDIOLOGENT leads all six task-setting combinations with non-overlapping confidence intervals (Table 2), by a similar margin whether both signals or the PPG alone is read. Absolute scores run higher under PPG, but this reflects its coarser label space and the correspondingly relaxed judging, so we read the gap to the best baseline, stable across settings, not the raw level. T2 is our hardest probe—ruling alternatives out is harder than naming a rhythm—and the smaller baselines fall below 6 there, listing possibilities without excluding them, while we tie each exclusion to a measured feature. T3, counterfactual reasoning, is where the general models are most competitive and our margin is narrowest.

Table 1: CARDIOLOGENT **window-level** classification: nine/seven-class macro-F1 (top) and binary normal-vs-arrhythmia balanced accuracy (bottom). $\pm$ is a 95% bootstrap CI.

	both	PPG
<i>macro-F1</i>		
CARDIOLOGENT (Ours)	0.449$\pm$0.026	0.425$\pm$0.027
GPT-4o	0.255 $\pm$ 0.045	0.147 $\pm$ 0.024
Gemini-2.5-Pro	0.186 $\pm$ 0.042	0.194 $\pm$ 0.020
Qwen2.5-VL-72B	0.160 $\pm$ 0.015	0.198 $\pm$ 0.011
MedGemma-27B	0.084 $\pm$ 0.002	0.131 $\pm$ 0.011
Llama-3.2-90B-Vision	0.021 $\pm$ 0.004	0.169 $\pm$ 0.003
<i>binary accuracy</i>		
CARDIOLOGENT (Ours)	0.94$\pm$0.01	0.92$\pm$0.01
GPT-4o	0.83 $\pm$ 0.01	0.56 $\pm$ 0.01
Gemini-2.5-Pro	0.61 $\pm$ 0.01	0.61 $\pm$ 0.01
Qwen2.5-VL-72B	0.58 $\pm$ 0.01	0.53 $\pm$ 0.01
MedGemma-27B	0.54 $\pm$ 0.01	0.53 $\pm$ 0.01
Llama-3.2-90B-Vision	0.50 $\pm$ 0.00	0.50 $\pm$ 0.00

Table 2: CARDIOLOGENT **window-level** reasoning per probe (SUM/16). $\pm$ is a 95% bootstrap CI.

	T1	T2	T3
<i>both (ECG+PPG)</i>			
CARDIOLOGENT (Ours)	12.51$\pm$0.11	10.50$\pm$0.17	12.88$\pm$0.19
GPT-4o	11.57 $\pm$ 0.16	8.99 $\pm$ 0.20	12.02 $\pm$ 0.19
Gemini-2.5-Pro	8.66 $\pm$ 0.16	9.94 $\pm$ 0.24	11.30 $\pm$ 0.21
Qwen2.5-VL-72B	10.74 $\pm$ 0.18	5.55 $\pm$ 0.14	8.21 $\pm$ 0.22
MedGemma-27B	9.75 $\pm$ 0.16	4.87 $\pm$ 0.09	7.21 $\pm$ 0.19
Llama-3.2-90B-Vision	4.47 $\pm$ 0.06	4.78 $\pm$ 0.09	7.58 $\pm$ 0.20
<i>PPG-only</i>			
CARDIOLOGENT (Ours)	12.90$\pm$0.12	8.26$\pm$0.13	12.19$\pm$0.21
GPT-4o	8.21 $\pm$ 0.18	7.01 $\pm$ 0.09	11.30 $\pm$ 0.16
Gemini-2.5-Pro	10.50 $\pm$ 0.17	6.56 $\pm$ 0.16	9.29 $\pm$ 0.17
Qwen2.5-VL-72B	11.03 $\pm$ 0.18	6.65 $\pm$ 0.11	7.63 $\pm$ 0.14
MedGemma-27B	10.60 $\pm$ 0.18	5.78 $\pm$ 0.09	7.22 $\pm$ 0.13
Llama-3.2-90B-Vision	5.80 $\pm$ 0.11	7.07 $\pm$ 0.11	8.38 $\pm$ 0.18

Grounding. Grounding is scored only on focal rhythms (PVC, PAC, VT, SVT), where the abnormality occupies part of the window; for sustained rhythms any span scores trivially. On overlap alone (focal-IoU) the baselines look competitive, but overlap rewards marking the whole window, so we take the per-sample focal-MCC as primary. There the separation is stark (Table 3): CARDIOLOGENT reaches 0.305 (both) and 0.319 (PPG), while every general model sits at or below 0.18—most within noise of zero, MedGemma negative—despite respectable IoU, shading the window rather than localizing the abnormality.

CARDIOLOGENT at the Patient Level

The patient level assembles the window verdicts, with the patient’s metadata and retrieved guidelines, into the three-part assessment; the baselines are asked the same three questions on the same record. CARDIOLOGENT leads on all three (Table 4), by the widest margin on P1, the integrated diagnosis

Table 3: **CARDIOLOGENT window-level grounding** on focal rhythms. MCC is primary; IoU alone rewards marking the whole window. both $n=230$, PPG $n=209$. $\pm$ is a 95% bootstrap CI.

	both		PPG	
	IoU	MCC	IoU	MCC
CARDIOLOGENT (Ours)	0.584	0.305	0.569	0.319
Llama-3.2-90B-Vision	0.533	0.001	0.544	0.000
GPT-4o	0.489	0.005	0.357	0.029
Gemini-2.5-Pro	0.440	0.178	0.535	0.109
Qwen2.5-VL-72B	0.292	0.056	0.433	0.037
MedGemma-27B	0.230	-0.006	0.182	-0.011

Table 4: **Patient-level** assessment on the three tasks, each scored on two axes (SUM/8). Top: the LLM judge (DeepSeek-V4-Pro) over all 65 test patients; bottom: two cardiologists, blind to system identity, over a 20-patient subset drawn at random but stratified to preserve the arrhythmia-class proportions of the full test set (per-class counts in the supplement). $\pm$ is a 95% cluster bootstrap CI ($B=10,000$).

	P1	P2	P3
<i>LLM judge (n=65)</i>			
CARDIOLOGENT (Ours)	4.66$\pm$0.61	4.37$\pm$0.64	4.55$\pm$0.58
Gemini-2.5-Pro	3.20 $\pm$ 0.50	3.31 $\pm$ 0.59	3.94 $\pm$ 0.54
GPT-4o	2.78 $\pm$ 0.38	2.83 $\pm$ 0.42	2.94 $\pm$ 0.36
Qwen2.5-VL-72B	2.12 $\pm$ 0.26	2.23 $\pm$ 0.21	2.86 $\pm$ 0.22
MedGemma-27B	1.91 $\pm$ 0.07	2.09 $\pm$ 0.23	2.89 $\pm$ 0.21
Llama-3.2-90B-Vision	2.06 $\pm$ 0.20	1.49 $\pm$ 0.29	2.91 $\pm$ 0.32
<i>Cardiologists (n=20, two raters)</i>			
CARDIOLOGENT (Ours)	4.58$\pm$1.32	4.70$\pm$1.08	5.03$\pm$1.01
Gemini-2.5-Pro	3.33 $\pm$ 1.36	3.75 $\pm$ 1.09	4.22 $\pm$ 0.99
GPT-4o	2.98 $\pm$ 1.11	3.48 $\pm$ 0.88	3.85 $\pm$ 0.95
Qwen2.5-VL-72B	1.48 $\pm$ 0.78	2.48 $\pm$ 0.54	2.92 $\pm$ 0.45
MedGemma-27B	1.38 $\pm$ 0.60	2.62 $\pm$ 0.54	2.80 $\pm$ 0.47
Llama-3.2-90B-Vision	1.00 $\pm$ 0.35	1.88 $\pm$ 0.21	2.25 $\pm$ 0.31

the assembly is built to reach. Two cardiologists blind to system identity rank the systems in the judge’s order and place CARDIOLOGENT first throughout; their agreement with the judge (ICC 0.74, 0.66) matches their agreement with each other (0.67, Figure 4), so the at-scale judge stands in for a clinician about as well as one clinician stands in for another.

We also score the diagnosis directly—the dominant rhythm named exactly, within its management family, or anywhere at all (supplement)—where CARDIOLOGENT leads at every grade and the baselines exceed a third only at the loosest.

Limitations

This work has several limitations. Our results come from one institution’s recordings, so the class balance, the patient mix, and the recording conditions are particular to it. Coverage is thinnest on the rare classes—SVT, VT, AV block—which hold the macro-F1 down and are the rhythms it would matter most to catch. The PPG-only setting has never been tested on

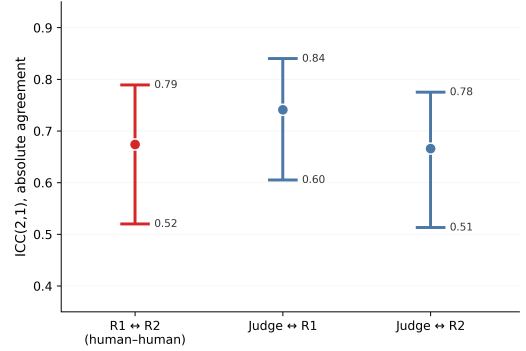


Figure 4: **Rater agreement** (ICC(2,1), absolute agreement, with 95% CI) on the 20-patient subset. The judge’s agreement with each cardiologist (0.74, 0.66) is comparable to the two cardiologists’ agreement with each other (0.67): the judge agrees with a cardiologist about as closely as the two cardiologists agree with each other.

a wearable: our pulse comes from a monitor, on a still patient with a fixed sensor. The evaluation is also retrospective—we score the decision itself, not whether acting on it would have changed the patient’s course—and the two deployment states are evaluated separately rather than as one pathway: every record carries both signals, so the hand-off from pulse-only surveillance to a combined read is never exercised.

Conclusion and Future Work

We formulate patient-level arrhythmia decision support as a task—what arrhythmia a monitored record holds, what it means for this patient, and what should be done about it—and present CARDIOLOGENT, a multi-agent system that spans it from the window reading to the decision, with each conclusion held to the guideline it cites. Across integrated diagnosis, clinical significance, and urgency and management, CARDIOLOGENT scores highest on every axis—under cardiologists and under the LLM judge alike—and its window-level macro-F1 (0.449) is roughly double the strongest baseline’s. We hope a decision a clinician can audit, rather than one more unexplained alert, is a step toward systems usable in continuous monitoring.

The limitations above each point to a next step: external validation on multi-centre data; targeted collection of the high-risk rhythms rather than more of the same recordings; validation on ambulatory recordings; the pathway itself, evaluated end to end—a wearable running the PPG agent alone for continuous, low-cost surveillance, and, on a finding that warrants it, the patient moving to a monitored setting where the ECG is added and the two agents combine, raising the level of scrutiny without changing systems; a retrospective study on a hospital’s own monitoring data; and, above all, placing the system in front of clinicians prospectively, to measure whether it changes what they do and to what end.

References

- Al-Khatib, S. M.; Stevenson, W. G.; Ackerman, M. J.; Bryant, W. J.; Callans, D. J.; Curtis, A. B.; Deal, B. J.; Dickfeld, T.; Field, M. E.; Fonarow, G. C.; et al. 2018. 2017 AHA/ACC/HRS guideline for management of patients with ventricular arrhythmias and the prevention of sudden cardiac death: a report of the American College of Cardiology/American Heart Association Task Force on Clinical Practice Guidelines and the Heart Rhythm Society. *Journal of the American College of Cardiology*, 72(14): e91–e220.
- Alnegheimish, S.; Nguyen, L.; Berti-Equille, L.; and Veeramachaneni, K. 2024. Large language models can be zero-shot anomaly detectors for time series? *arXiv preprint arXiv:2405.14755*.
- Chicco, D.; Tötsch, N.; and Jurman, G. 2021. The Matthews correlation coefficient (MCC) is more reliable than balanced accuracy, bookmaker informedness, and markedness in two-class confusion matrix evaluation. *BioData mining*, 14(1): 13.
- Chung, H.; Oh, J.; Kyung, D.; Kim, J.; Kwon, Y.; Kim, M.-G.; and Choi, E. 2026. ECG-Agent: On-Device Tool-Calling Agent for ECG Multi-Turn Dialogue. *arXiv preprint arXiv:2601.20323*.
- Comanici, G.; Bieber, E.; Schaeckermann, M.; Pasupat, I.; Sachdeva, N.; Dhillon, I.; Blistein, M.; Ram, O.; Zhang, D.; Rosen, E.; et al. 2025. Gemini 2.5: Pushing the frontier with advanced reasoning, multimodality, long context, and next generation agentic capabilities. *arXiv preprint arXiv:2507.06261*.
- Du, Y.; Li, S.; Torralba, A.; Tenenbaum, J. B.; and Mordatch, I. 2023. Improving factuality and reasoning in language models through multiagent debate. *arXiv preprint arXiv:2305.14325*.
- Eun, D.-I.; Shim, K.; Lee, H.; Lim, Y.; Lim, H.; Lee, H.; Lee, J.; and Lee, H.-C. 2026. Vitaldb arrhythmia database: An anesthesiologist-validated large-scale intraoperative arrhythmia dataset with beat and rhythm labels. *Scientific Data*, 13(1): 838.
- Fallahpour, A.; Ma, J.; Munim, A.; Lyu, H.; and Wang, B. 2025. MedRAX: Medical Reasoning Agent for Chest X-ray. In *Proceedings of the 42nd International Conference on Machine Learning*.
- Frykman, V.; Frick, M.; Jensen-Urstad, M.; Östergren, J.; and Rosenqvist, M. 2001. Asymptomatic versus symptomatic persistent atrial fibrillation: clinical and noninvasive characteristics. *Journal of internal medicine*, 250(5): 390–397.
- Fu, D.-g. 2015. Cardiac arrhythmias: diagnosis, symptoms, and treatments. *Cell biochemistry and biophysics*, 73(2): 291–296.
- Gu, J.; Jiang, X.; Shi, Z.; Tan, H.; Zhai, X.; Xu, C.; Li, W.; Shen, Y.; Ma, S.; Liu, H.; Wang, S.; Zhang, K.; Wang, Y.; Gao, W.; Ni, L.; and Guo, J. 2025. A Survey on LLM-as-a-Judge. *arXiv:2411.15594*.
- Guo, Y.; Wang, H.; Zhang, H.; Liu, T.; Liang, Z.; Xia, Y.; Yan, L.; Xing, Y.; Shi, H.; Li, S.; et al. 2019. Mobile photoplethysmographic technology to detect atrial fibrillation. *Journal of the American College of Cardiology*, 74(19): 2365–2375.
- Han, D.; Bashar, S. K.; Mohagheghian, F.; Ding, E.; Whitcomb, C.; McManus, D. D.; and Chon, K. H. 2020. Premature atrial and ventricular contraction detection using photoplethysmographic data from a smartwatch. *Sensors*, 20(19): 5683.
- Hannun, A. Y.; Rajpurkar, P.; Haghpanahi, M.; Tison, G. H.; Bourn, C.; Turakhia, M. P.; and Ng, A. Y. 2019. Cardiologist-level arrhythmia detection and classification in ambulatory electrocardiograms using a deep neural network. *Nature medicine*, 25(1): 65–69.
- He, K.; Chen, X.; Xie, S.; Li, Y.; Dollár, P.; and Girshick, R. 2022. Masked autoencoders are scalable vision learners. In *Proceedings of the IEEE/CVF conference on computer vision and pattern recognition*, 16000–16009.
- Hurst, A.; Lerer, A.; Goucher, A. P.; Perelman, A.; Ramesh, A.; Clark, A.; Ostrow, A.; Welihinda, A.; Hayes, A.; Radford, A.; et al. 2024. Gpt-4o system card. *arXiv preprint arXiv:2410.21276*.
- Jin, J.; Wang, H.; Lan, X.; Li, J.; Li, H.; and Hong, S. 2026a. UniECG: Understanding and Generating ECG in One Unified Model. *arXiv:2509.18588*.
- Jin, J.; Wang, H.; Wu, X.; Fang, X.; Lan, X.; Wang, Z.; Zhang, D.; Liu, B.; Zhang, Y.; Wu, X.; et al. 2026b. ECG-R1: Protocol-Guided and Modality-Agnostic MLLM for Reliable ECG Interpretation. *arXiv preprint arXiv:2602.04279*.
- Jin, M.; Wang, S.; Ma, L.; Chu, Z.; Zhang, J.; Shi, X.; Chen, P.-Y.; Liang, Y.; Li, Y.-F.; Pan, S.; et al. 2024. Time-llm: Time series forecasting by reprogramming large language models. In *International conference on learning representations*, volume 2024, 23857–23880.
- Joglar, J. A.; Chung, M. K.; Armbruster, A. L.; Benjamin, E. J.; Chyou, J. Y.; Cronin, E. M.; Deswal, A.; Eckhardt, L. L.; Goldberger, Z. D.; Gopinathannair, R.; et al. 2024. 2023 ACC/AHA/ACCP/HRS guideline for the diagnosis and management of atrial fibrillation: a report of the American College of Cardiology/American Heart Association Joint Committee on Clinical Practice Guidelines. *Circulation*, 149(1): e1–e156.
- Kim, A.; Kim, K.; and Yoon, S. 2024. DEBATE: Devil’s advocate-based assessment and text evaluation. In *Findings of the Association for Computational Linguistics: ACL 2024*, 1885–1897.
- Kim, Y.; Park, C.; Jeong, H.; Chan, Y. S.; Xu, X.; McDuff, D.; Lee, H.; Ghassemi, M.; Breazeal, C.; and Park, H. W. 2024. MDAgents: An Adaptive Collaboration of LLMs for Medical Decision-Making. *arXiv:2404.15155*.
- Kusumoto, F. M.; Schoenfeld, M. H.; Barrett, C.; Edgerton, J. R.; Ellenbogen, K. A.; Gold, M. R.; Goldschlager, N. F.; Hamilton, R. M.; Joglar, J. A.; Kim, R. J.; et al. 2019. 2018 ACC/AHA/HRS guideline on the evaluation and management of patients with bradycardia and cardiac conduction delay: a report of the American College of Cardiology/American Heart Association Task Force on Clinical Practice Guidelines and the Heart Rhythm Society. *Journal of the American College of Cardiology*, 74(7): e51–e156.
- Lan, X.; Wu, F.; He, K.; Zhao, Q.; Hong, S.; and Feng, M. 2026. Gem: Empowering mllm for grounded ecg under-

- standing with time series and images. *Advances in Neural Information Processing Systems*, 38: 94421–94455.
- Lee, H.-C.; Park, Y.; Yoon, S. B.; Yang, S. M.; Park, D.; and Jung, C.-W. 2022. VitalDB, a high-fidelity multi-parameter vital signs database in surgical patients. *Scientific Data*, 9(1): 279.
- Li, B.; Yan, T.; Pan, Y.; Luo, J.; Ji, R.; Ding, J.; Xu, Z.; Liu, S.; Dong, H.; Lin, Z.; and Wang, Y. 2024. MMedAgent: Learning to Use Medical Tools with Multi-modal Agent. arXiv:2407.02483.
- Li, H.; Li, Z.; Mao, Y.; Liu, Z.; Sun, Z.; and Huang, Z. 2025a. anyECG-chat: A Generalist ECG-MLLM for Flexible ECG Input and Multi-Task Understanding. arXiv:2506.00942.
- Li, J.; Aguirre, A. D.; Junior, V. M.; Jin, J.; Liu, C.; Zhong, L.; Sun, C.; Clifford, G.; Brandon Westover, M.; and Hong, S. 2025b. An electrocardiogram foundation model built on over 10 million recordings. *Nejm ai*, 2(7): AIoa2401033.
- Liang, T.; He, Z.; Jiao, W.; Wang, X.; Wang, Y.; Wang, R.; Yang, Y.; Shi, S.; and Tu, Z. 2024. Encouraging divergent thinking in large language models through multi-agent debate. In *Proceedings of the 2024 conference on empirical methods in natural language processing*, 17889–17904.
- Liu, H.; Li, C.; Wu, Q.; and Lee, Y. J. 2023a. Visual instruction tuning. *Advances in neural information processing systems*, 36: 34892–34916.
- Liu, R.; Bai, Y.; Yue, X.; and Zhang, P. 2026. Teaching multimodal LLMs to comprehend 12-lead electrocardiographic images. *npj Digital Medicine*, 9(1): 349.
- Liu, Y.; Iter, D.; Xu, Y.; Wang, S.; Xu, R.; and Zhu, C. 2023b. G-Eval: NLG Evaluation using GPT-4 with Better Human Alignment. arXiv:2303.16634.
- Liu, Z.; Zhou, B.; Jiang, Z.; Chen, X.; Li, Y.; Tang, M.; and Miao, F. 2022. Multiclass arrhythmia detection and classification from photoplethysmography signals using a deep convolutional neural network. *Journal of the American Heart Association*, 11(7): e023555.
- Matthews, B. W. 1975. Comparison of the predicted and observed secondary structure of T4 phage lysozyme. *Biochimica et Biophysica Acta (BBA)-Protein Structure*, 405(2): 442–451.
- Page, R. L.; Joglar, J. A.; Caldwell, M. A.; Calkins, H.; Conti, J. B.; Deal, B. J.; Estes III, N. M.; Field, M. E.; Goldberger, Z. D.; Hammill, S. C.; et al. 2016. 2015 ACC/AHA/HRS guideline for the management of adult patients with supraventricular tachycardia: a report of the American College of Cardiology/American Heart Association Task Force on Clinical Practice Guidelines and the Heart Rhythm Society. *Journal of the American College of Cardiology*, 67(13): e27–e115.
- Pereira, T.; Tran, N.; Gadhoumi, K.; Pelter, M. M.; Do, D. H.; Lee, R. J.; Colorado, R.; Meisel, K.; and Hu, X. 2020. Photoplethysmography based atrial fibrillation detection: a review. *NPJ digital medicine*, 3(1): 3.
- Perez, M. V.; Mahaffey, K. W.; Hedlin, H.; Rumsfeld, J. S.; Garcia, A.; Ferris, T.; Balasubramanian, V.; Russo, A. M.; Rajmane, A.; Cheung, L.; et al. 2019. Large-scale assessment of a smartwatch to identify atrial fibrillation. *New England Journal of Medicine*, 381(20): 1909–1917.
- Quinlan, P.; Li, Q.; and Zhu, X. 2026. Chat-ts: Enhancing multi-modal reasoning over time-series and natural language data. In *Proceedings of the 19th Conference of the European Chapter of the Association for Computational Linguistics (Volume 1: Long Papers)*, 5621–5647.
- Ribeiro, A. H.; Ribeiro, M. H.; Paixão, G. M.; Oliveira, D. M.; Gomes, P. R.; Canazart, J. A.; Ferreira, M. P.; Andersson, C. R.; Macfarlane, P. W.; Meira Jr, W.; et al. 2020. Automatic diagnosis of the 12-lead ECG using a deep neural network. *Nature communications*, 11(1): 1760.
- Sellergren, A.; Kazemzadeh, S.; Jaroensri, T.; Kiraly, A.; Traverse, M.; Kohlberger, T.; Xu, S.; Jamil, F.; Hughes, C.; Lau, C.; Chen, J.; Mahvar, F.; Yatziv, L.; Chen, T.; Sterling, B.; Baby, S. A.; Baby, S. M.; Lai, J.; Schmidgall, S.; Yang, L.; Chen, K.; Bjornsson, P.; Reddy, S.; Brush, R.; Philbrick, K.; Asiedu, M.; Mezerreg, I.; Hu, H.; Yang, H.; Tiwari, R.; Jansen, S.; Singh, P.; Liu, Y.; Azizi, S.; Kamath, A.; Ferret, J.; Pathak, S.; Vieillard, N.; Merhej, R.; Perrin, S.; Matejovicova, T.; Ramé, A.; Riviere, M.; Rouillard, L.; Mesnard, T.; Cideron, G.; bastien Grill, J.; Ramos, S.; Yvinec, E.; Casbon, M.; Buchatskaya, E.; Alayrac, J.-B.; Lepikhin, D.; Feinberg, V.; Borgeaud, S.; Andreev, A.; Hardin, C.; Dadashi, R.; Hussenot, L.; Joulin, A.; Bachem, O.; Matias, Y.; Chou, K.; Hassidim, A.; Goel, K.; Farabet, C.; Barral, J.; Warkentin, T.; Shlens, J.; Fleet, D.; Cotruta, V.; Sansevero, O.; Martins, G.; Kirk, P.; Rao, A.; Shetty, S.; Steiner, D. F.; Kirmizibayrak, C.; Pilgrim, R.; Golden, D.; and Yang, L. 2026. MedGemma Technical Report. arXiv:2507.05201.
- Siontis, K. C.; Noseworthy, P. A.; Attia, Z. I.; and Friedman, P. A. 2021. Artificial intelligence-enhanced electrocardiography in cardiovascular disease management. *Nature Reviews Cardiology*, 18(7): 465–478.
- Wan, Z.; Liu, C.; Wang, X.; Tao, C.; Shen, H.; Xiong, J.; Arcucci, R.; Yao, H.; and Zhang, M. 2025. MEIT: Multimodal electrocardiogram instruction tuning on large language models for report generation. In *Findings of the association for computational linguistics: ACL 2025*, 14510–14527.
- Wang, X.; Chang, C.; Cao, D.; Han, K.; Sun, F.; Huang, Y.; Wang, M.; Xu, C.; Luo, X.; Yan, R.; Zhang, X.; Hu, X.; Liu, Y.; Sun, Y.; Wang, W.; and Yang, C. 2026. Position: Beyond Prediction: Toward Verifiable Physiological Waveform Reasoning with Foundation Models and Agentic LLMs. *TechRxiv*, 2026(0212).
- Yang, K.; Hong, M.; Zhang, J.; Luo, Y.; Zhao, S.; Zhang, O.; Yu, X.; Zhou, J.; Yang, L.; Zhang, P.; et al. 2025. ECG-LM: understanding electrocardiogram with a large language model. *Health Data Science*, 5: 0221.
- Zeppenfeld, K.; Tfelt-Hansen, J.; De Riva, M.; Winkel, B. G.; Behr, E. R.; Blom, N. A.; Charron, P.; Corrado, D.; Dagres, N.; De Chillou, C.; et al. 2022. 2022 ESC Guidelines for the management of patients with ventricular arrhythmias and the prevention of sudden cardiac death: Developed by the task force for the management of patients with ventricular arrhythmias and the prevention of sudden cardiac death of the

European Society of Cardiology (ESC) Endorsed by the Association for European Paediatric and Congenital Cardiology (AEPC). *European heart journal*, 43(40): 3997–4126.

Zhang, Y.; Bunting, K. V.; Champai, A.; Wang, X.; Lu, W.; Thorley, A.; Hothi, S. S.; Qiu, Z.; Buyukates, B.; Kotecha, D.; et al. 2025. Cardiac-agents: A multimodal framework with hierarchical adaptation for cardiac care support. *arXiv preprint arXiv:2508.13256*.

Zhao, Y.; Kang, J.; Zhang, T.; Han, P.; and Chen, T. 2025. Ecg-chat: A large ecg-language model for cardiac disease diagnosis. In *2025 IEEE International Conference on Multimedia and Expo (ICME)*, 1–6. IEEE.

Zheng, L.; Chiang, W.-L.; Sheng, Y.; Zhuang, S.; Wu, Z.; Zhuang, Y.; Lin, Z.; Li, Z.; Li, D.; Xing, E. P.; Zhang, H.; Gonzalez, J. E.; and Stoica, I. 2023. Judging LLM-as-a-Judge with MT-Bench and Chatbot Arena. *arXiv:2306.05685*.

Zhou, T.; Niu, P.; Sun, L.; Jin, R.; et al. 2023. One fits all: Power general time series analysis by pretrained lm. *Advances in neural information processing systems*, 36: 43322–43355.

Zhou, Y.; Zhang, P.; Song, M.; Zheng, A.; Lu, Y.; Liu, Z.; Chen, Y.; and Xi, Z. 2024. Zodiac: A cardiologist-level llm framework for multi-agent diagnostics. *arXiv preprint arXiv:2410.02026*.

Zhou, Z.; and Yu, R. 2024. Can llms understand time series anomalies? *arXiv preprint arXiv:2410.05440*.

Related Work

Language models for time series. Whether a language model can read a raw numerical signal has been asked in a fairly clear progression. Early work reprogrammed a frozen model for forecasting and classification, showing a general backbone could be adapted to time series at all (Zhou et al. 2023; Jin et al. 2024). Attention then moved to anomaly detection, where a model was prompted zero-shot to flag the abnormal part of a series (Alnegheimish et al. 2024); the limits of that prompting were examined soon after, finding the ability shallow and inconsistent (Zhou and Yu 2024), and more recent systems pair the signal with text for multimodal reasoning rather than raw prompting (Quinlan, Li, and Zhu 2026). The throughline is that a general model reads a raw trace only weakly, which is why our specialists are adapted to the waveform rather than prompted on it; and unlike this line, the signal here is not the end product but a tool the agent calls.

ECG and signal language models. Attaching a language model to the ECG began with instruction tuning for report generation and question answering, which gave a waveform a textual reading (Wan et al. 2025; Zhao et al. 2025; Yang et al. 2025). Attention then turned to making that reading trustworthy rather than fluent: grounding each statement in the time series it refers to (Lan et al. 2026), following an explicit diagnostic protocol (Jin et al. 2026b), and unifying understanding with generation in one model (Jin et al. 2026a); a parallel thread pretrained ECG foundation models at scale (Li et al. 2025b,a). Most recently the reading has been wrapped in agents. Zodiac (Zhou et al. 2024) assembles a patient-specific report from a monitoring study under cardiologist validation, but its arrhythmias arrive already detected—a representative strip and a burden table selected upstream—and it stops at interpretation, leaving the management decision to the reviewing cardiologist; its validation scores the report—accuracy against the patient’s own data, completeness, organization, succinctness—not the diagnosis against a label. CardAIC-Agents (Zhang et al. 2025) instead ingests the raw signal and coordinates tools over it, but its task is disease-level diagnosis (heart failure, myocardial infarction) scored by classification accuracy, not the rhythm-level identification of arrhythmias, and its own signal-level wave detection is reported as suboptimal; and while its intermediate visual outputs are reviewed by cardiologists, the management it plans is not among what is scored. A tool-calling dialogue agent (Chung et al. 2026) reads a single ECG over a short exchange. Across these, the object is one recording or one disease label, returned as a reading or a report—none assembles a monitored record, window by window, into a patient-level arrhythmia profile and carries it through to an audited management decision. Our specialist sits at the end of this line, but the system around it does what the line has not: it reconciles a second and weaker signal against the ECG, and it carries the per-window reading up into a patient-level assessment.

LLM agents. Beyond a single model, agentic systems let a language model plan, call tools, and revise its answer, and multi-agent debate improves reasoning by having several

Table 5: Nine-class window counts over the full pool and the patient-disjoint train/val/test split.

class	total	train	val	test
normal	75,039	51,827	12,069	11,143
atrial fibrillation	44,583	36,291	3,868	4,424
PVC	10,940	9,586	789	565
PAC	9,997	8,491	990	516
MAT	9,683	7,801	900	982
SND	6,493	5,312	648	533
SVT	5,215	4,383	515	317
AVB	3,528	2,819	361	348
VT	548	449	54	45
total	166,026	126,959	20,194	18,873

Table 6: Class distribution of the 3,000-window evaluation sample (ECG nine-class).

class	<i>n</i>	%
normal	1,771	59.0
atrial fibrillation	703	23.4
MAT	156	5.2
PVC	90	3.0
SND	85	2.8
PAC	82	2.7
AVB	55	1.8
SVT	51	1.7
VT	7	0.2
total	3,000	

instances argue a question to an adjudicated verdict (Du et al. 2023; Liang et al. 2024; Kim, Kim, and Yoon 2024). Our fusion is an instance of this pattern, but with a difference the setting forces: the debaters are not interchangeable and consensus is not the goal. Because the ECG and the PPG see physically different things, authority is granted feature by feature to the signal that can see a feature, and a referee adjudicates rather than a majority—so a correct minority reading is not voted down.

Medical multi-agent systems. As agents matured, the pattern was carried into medicine, first by wrapping individual clinical tools in a single tool-calling agent (Li et al. 2024), then by coordinating several specialized tools under a reasoning loop for a given modality (Fallahpour et al. 2025), and more recently by adapting the collaboration itself—varying structure with case difficulty (Kim et al. 2024) or specializing agents to a signal such as the ECG (Chung et al. 2026). These share our use of a language-model agent over clinical tools, but each stays within a single modality and a single interaction—one image or one ECG, answered in a turn or a short dialogue. What is specific here is the substrate and the target: a continuously monitored pair of signals rather than one image or a single ECG, and an assessment assembled window by window across a whole record into a diagnosis, its significance, and a management decision.

Dataset Details

A window spans twelve consecutive R-peaks rather than a fixed duration—one leading, ten central, one trailing—which at 60–100 bpm is roughly seven to twelve seconds at 100 Hz; windows are cut on the PPG peaks and the ECG is aligned to them by the pulse transit time, the same offset that carries the labels across (main text). The sliding stride is set per class to counteract imbalance: the rare rhythms (SVT, VT, MAT, AVB) advance one peak at a time ($\sim 92\%$ overlap) so nearly every occurrence is kept, while common windows and *normal* advance three ($\sim 75\%$) and are thereby subsampled. An episode is never split across a boundary—a window either contains it in full or lies wholly within it, and partial straddles are discarded—and a window is dropped if any sample fails quality control (a PPG segment flagged by the signal-quality index, an ECG beat marked bad, or a discard label in the mask) or if it holds fewer than twelve annotated beats.

The resulting distribution is dominated by *normal* and atrial fibrillation, with the ventricular and conduction rhythms rare—the imbalance the class-dependent stride is meant to soften. Table 5 gives the nine-class counts over the full pool and the patient-disjoint split, and Table 6 the 3,000-window evaluation sample. The split is by patient and stratified to preserve the overall class proportions across train, validation, and test, so no patient appears on two sides and each split reflects the same rhythm mix. The 3,000-window evaluation sample is drawn from the test set at random under two constraints—every test patient is represented and the class proportions are preserved—so that it stays representative while keeping the LLM-judge and API cost of scoring reasoning tractable.

Specialist QA. The specialist’s Stage-2 targets—the two-to-three-sentence readings it is trained to produce—are written once, offline, before any training. For each window we render its plot with the detected R-peaks marked and pair it with two text inputs: the window’s arrhythmia annotation (the benchmark label, with per-beat type, amplitude, QRS width, and R–R timing) as a ground-truth *anchor*, and a set of signal measurements summarizing rate and rhythm, the rate profile, and per-interval morphology. A multimodal model (GPT-4o) reads the plot and these inputs and writes an observation-to-evidence-to-diagnosis answer whose conclusion matches the anchor and whose numbers are drawn from the supplied measurements. The anchor guides accuracy but is instructed not to be copied verbatim, so the target reads as a reading of the trace rather than a restatement of the label.

Reasoning benchmark (T1/T2/T3). The reasoning probes are evaluated against a gold reference answer built per window in a generate–review–regenerate pipeline. Each window carries up to three fixed-template questions by its ground-truth class: T1 (diagnosis with supporting features) on every window, T2 (differential, excluding each alternative) on every arrhythmia, and T3 (the one change that would alter the call) on the confusable classes (AF, MAT, PVC, PAC, VT, SVT); under PPG the P-wave and QRS-morphology cues are replaced by pulse-wave morphology and the differential is told to keep any ECG-only-separable pair merged. Golds are generated by GPT-4o from four inputs—the ground-truth

diagnosis as anchor, a short expert description, the measured signal features, and the question—under the constraint that the answer conclude with the anchor diagnosis and cite only numbers present in those features, so the reference is grounded in the label rather than in the generator’s free recall. Every gold is then audited by a *different* model (Gemini-2.5-Pro) against the same anchor and evidence and graded CORRECT/MINOR/INCORRECT; the INCORRECT cases are regenerated with the audit notes attached. Using one model to write and a different one to check, with the ground-truth label—not either model’s opinion—as the standard, keeps the benchmark from being anchored to a single model’s judgement.


Specialist Details

Each specialist is trained in two stages (Figure 5). A *time-series encoder* is first pretrained on the raw waveform by masked reconstruction, adapting the masked-autoencoder scheme (He et al. 2022) to the signal: it is split into patches, a fraction are masked, and a Transformer encoder–decoder reconstructs them, giving a representation of local morphology without any labels. Once pretrained, the decoder is discarded and only the encoder is kept, supplying the per-patch features used in Stage 2. The *reasoning model* is then a LLaVA-7B model (Liu et al. 2023a)—a Llama-7B backbone with a CLIP vision encoder—that reads a plot of the window, encoded by CLIP, alongside a text prompt carrying the measured features of the signal and the per-patch features from the Stage-1 encoder. Presenting a raw signal together with text so that a language model can reason over both (Alnegheimish et al. 2024; Zhou and Yu 2024; Quinlan, Li, and Zhu 2026) is the setting this follows, here for the ECG and the pulse rather than a generic series. It is fine-tuned to answer, in two or three sentences, whether the window contains an arrhythmia and, if so, its class and the evidence.

The text prompt is where the grounding enters: every quantity in it is a measured value or an encoder feature, not a number the model is left to read off the trace. Its fields, named as in Figure 5, are:

- **{arrhythmia taxonomy}** — the nine (ECG) or seven (PPG) classes with their defining criteria, including the QRS-width and P-wave distinctions the pulse cannot resolve.
- **{det feat}** — the measured features: beat count, heart rate, median R–R, R–R coefficient of variation, maximum $|\Delta R-R|$, median QRS width, and premature-beat count.
- **{peak coordinates}** — the detected R-peaks as (index, amplitude).
- **{raw values}** — the min–max normalized waveform, subsampled.
- **{timeseries patch}** — per-patch embeddings from the Stage-1 encoder (the morphological hints).
- **{clinical reasoning}** — a scaffold laying out the R–R, QRS-width, and amplitude sequences and walking four perspectives (rhythm regularity, rate profile, QRS morphology, premature-beat morphology) toward a conclusion.

Stage 1: Time series encoder

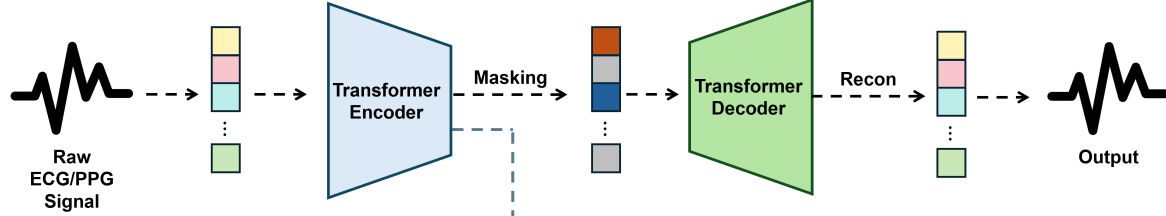

Stage 2: Arrhythmia reasoning model

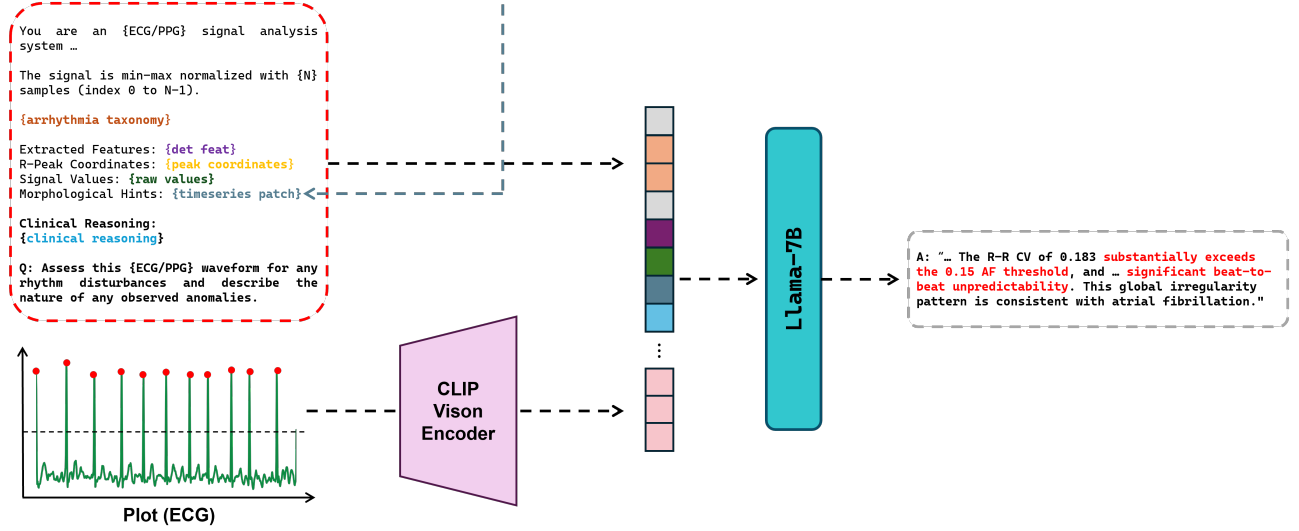


Figure 5: The specialist. Stage 1 pretrains a time-series encoder by masked reconstruction of the raw waveform. Stage 2 is the LLaVA-7B reasoning model (Llama-7B with a CLIP vision encoder) reading a plot of the window together with a prompt built from the measured features and the Stage-1 encoder, whose per-patch features enter that prompt as **{timeseries patch}**.

The prompt closes with the question—assess the window and describe any anomaly. On the PPG, which cannot resolve QRS width or premature-beat morphology, the QRS and premature perspectives drop out and only rhythm regularity, rate, and perfusion remain—the reason its taxonomy merges to seven classes.

Training objectives. Each stage optimizes a single loss term. Stage 1 is a peak-weighted masked reconstruction: the loss is a mean-squared error computed only at masked positions, upweighted around peaks. With M_i the masked positions of sample i and P_i the peaks detected on its target trace,

$$\mathcal{L}_1 = \frac{1}{B} \sum_{i=1}^B \text{mean}_{t \in M_i} \left[w_i(t) (\hat{x}_i(t) - x_i(t))^2 \right],$$

$$w_i(t) = \begin{cases} 5 & \text{if } \exists p \in P_i : |t - p| \leq 10, \\ 1 & \text{otherwise.} \end{cases}$$

Only masked positions contribute; predictions are interpolated back to each sample’s original resolution before the

error is taken; and the peaks are found on the target signal itself (local maxima above the mean, minimum spacing), not from ground-truth annotations, so the objective needs no labels. Stage 2 is the standard answer-only autoregressive cross-entropy,

$$\mathcal{L}_2 = - \sum_{t \in A} \log p_{\theta}(y_t \mid y_{<t}, \text{prompt}, z_{\text{enc}}),$$

summed over the answer tokens A (with the question and padding tokens masked out). Neither stage uses an auxiliary term.

Training. Stage 1 pretrains the encoder as a TimeMAE on each window interpolated to 1024 samples, with patch size 16 (64 patches), mask ratio 0.75, and rotary position embeddings, over a 6-layer 8-head Transformer ($d_{\text{model}}=512$); it is trained for 100 epochs (batch 4096, learning rate 10^{-4}). Stage 2 fine-tunes LLaVA-7B with a 64-query projector; only the Stage-1 signal encoder is frozen, and everything else—the projector, the CLIP vision encoder, and the language backbone—is trained, for 200k steps (learning rate

Table 7: **ECG and PPG specialist** classification (tool level). ECG is 9-class; PPG is 7-class.

	macro-F1	bal. acc	raw acc
<i>ECG (9-class)</i>			
Specialist (Ours)	0.393	0.423	0.762
Gemini-2.5-Pro	0.131	0.148	0.314
Qwen2.5-VL-72B	0.093	0.117	0.587
MedGemma-27B	0.090	0.127	0.578
Llama-3.2-90B-Vision	0.089	0.094	0.347
GPT-4o	0.080	0.126	0.335
<i>PPG (7-class)</i>			
Specialist (Ours)	0.367	0.380	0.783
Qwen2.5-VL-72B	0.198	0.214	0.514
Gemini-2.5-Pro	0.194	0.212	0.560
Llama-3.2-90B-Vision	0.169	0.171	0.518
GPT-4o	0.147	0.162	0.549
MedGemma-27B	0.131	0.146	0.470

5×10^{-6} , bf16, FSDP with gradient checkpointing) on two A100-80GB GPUs.

Results. Each specialist is evaluated alone on its own signal—the ECG specialist on the nine-class task, the PPG specialist on the seven-class merged task—against the same five baselines, each given that signal, its plot, and the taxonomy. The written reading is scored by the same judge on four axes (1–4: correctness, faithfulness, validity, completeness; SUM/16). At 7B, several times smaller than every baseline, each specialist leads all of them (Table 7): macro-F1 0.393 against 0.131 on the ECG and 0.367 against 0.198 on the PPG. Raw accuracy separates them far less (0.762 against Qwen’s 0.587), which is why we do not report it alone—Qwen holds that score at a macro-F1 of 0.093. The specialist also leads every reasoning axis in both settings (SUM 12.57 on the ECG and 13.23 on the PPG, against 9.65 and 10.52 for the strongest baseline; Table 8), by the widest margin on completeness (3.47 vs. 2.50) and validity (3.06 vs. 2.33), and the narrowest on faithfulness (2.64 vs. 2.28).

Agent Details

Deterministic tools. The analysts are deterministic tools computed from the raw signal with NeuroKit2 and SciPy, and their role is to re-examine the specialist’s reading: each returns a measured quantity that the class and evidence the specialist proposed must be checked against, so the language model argues over tool outputs rather than the waveform. Three tools form a shared backbone across both signals—peak detection, rhythm analysis over the peak-to-peak series, and an artifact/quality gate—and the rest are modality-specific, which is where sensor physics enters. The ECG agent has seven tools in all, including those the pulse cannot reproduce:

- **QRS morphology** — beat width and shape, separating wide ventricular beats from narrow ones.
- **Atrial activity** — P-wave presence and organization, separating AF from MAT and SVT.

Table 8: **ECG and PPG specialist** reasoning (tool level), LLM-judge on four axes (1–4): correctness, faithfulness, validity, completeness; SUM/16.

	corr	faith	valid	compl	SUM
<i>ECG (9-class)</i>					
Specialist (Ours)	3.40	2.64	3.06	3.47	12.57
Qwen2.5-VL-72B	2.72	2.28	2.33	2.32	9.65
MedGemma-27B	2.65	1.88	2.18	2.34	9.06
GPT-4o	2.01	1.97	1.88	2.50	8.35
Gemini-2.5-Pro	1.97	1.67	1.82	2.08	7.53
Llama-3.2-90B-Vision	1.99	1.56	1.70	1.87	7.11
<i>PPG (7-class)</i>					
Specialist (Ours)	3.35	3.25	3.25	3.38	13.23
GPT-4o	2.65	2.64	2.60	2.63	10.52
Gemini-2.5-Pro	2.64	2.57	2.53	2.61	10.35
Qwen2.5-VL-72B	2.54	2.70	2.52	2.50	10.26
Llama-3.2-90B-Vision	2.53	2.66	2.48	2.46	10.12
MedGemma-27B	2.61	2.28	2.29	2.30	9.48

- **AV conduction** — pause structure, separating sinus-node dysfunction from AV block.

The PPG agent has nine, weighted toward what the pulse does carry:

- **Perfusion** — perfusion index, for the haemodynamic significance of a fast rhythm.
- **Pulse-wave / beat morphology** — pulse-wave shape and beat matching against a template, for ectopy and aberrant beats.

One of these tools localizes the abnormal span independently of the class, for grounding.

Within-agent debate. Inside each agent the window is read by the same analyst–devil’s-advocate debate used for window fusion (main text), with two differences. It has no fixed referee: three analysts issue independent verdicts and revise over two rounds to a majority position (Du et al. 2023). And it is gated by a critic—if the devil’s advocate’s objection is left unresolved, the reading is failed and returned for another round of tool use, so evidence is re-gathered rather than assumed. The physical constraints still hold: the PPG keeps the merged VT/MAT and PVC/PAC classes and defers to the ECG, and a PPG claim about a P wave or QRS width is treated as a fabrication.

Patient-Level Pipeline

Retrieval. Retrieval draws on two sources: an online PubMed search through the NCBI E-utilities, and a curated local base of about twenty guideline snippets, each a short passage with its PMID embedded; local topics span R–R thresholds, AF on the PPG, ectopy and PVC burden, SVT, rate ranges, and the class-specific management references (2023 AF and CHA₂DS₂-VASc, 2017 VT and sudden death, 2018 bradycardia and AV block). Matching is TF-IDF cosine similarity over the local base and a live query against PubMed. An LLM draws two or three queries from the fused patient state and metadata—about significance, management,

and risk, not the diagnosis, which is already settled—and the top snippets are retrieved; the PMID in each is surfaced as an explicit tag that later stages cite verbatim. For the reported numbers we fix retrieval to the curated base so that every citation resolves to a manually verified reference and the evaluation is exactly reproducible; the PubMed path is used when the system runs live.

Critic. The critic checks the draft P1–P3 against a fixed rule set and edits only violations:

- **Grounding** — remove any finding not supported by the ECG or PPG state.
- **Citations** — every risk or management claim must cite a PMID that appears in the retrieved set; mismatches are remapped, fabricated PMIDs deleted, missing ones added.
- **Anticoagulation** — not indicated for ectopy or on CHA₂DS₂-VASc alone; a recommendation without documented AF is corrected.
- **Redundant tests** — the agent already holds the ECG, so any “obtain an ECG” line is replaced with the ECG’s actual finding.
- **Class specificity** — a management line that would read the same for any rhythm is rewritten to the diagnosed class.

If the critic’s own output is malformed, the draft is kept.

Patient-level debate. The three patient stages reuse the same debate primitive, differing only in what they weigh. Aggregation runs analysts to a referee *without* a devil’s advocate and is descriptive: one lens tracks the dominant rhythm and burden, the other enumerates every distinct arrhythmia without ranking, and the PPG records its own AF confidence without asserting an atrial rhythm it cannot confirm. Fusion and reasoning both add a devil’s advocate. Fusion reconciles the two states—ECG for QRS and P-waves, PPG for pulse regularity and perfusion—and then applies severity, its referee suppressing PPG-only AF false positives and taking the dominant rhythm as the one the ECG shows across most windows. Reasoning adds the metadata to produce P1–P3, its devil attacking mis-prioritized claims (benign ectopy raised over a present AF or VT, urgency not matching the evidence) and its referee requiring that P2 and P3 cite only retrieved PMIDs and that P3 be class-specific.

Scoring Rubrics

Why anchored rubrics. We score reasoning with an LLM judge because the quality we care about—whether a reading is correct, grounded in real evidence, and clinically actionable—has no reference string to match and is exactly what reference-free LLM evaluation was introduced to capture (Zheng et al. 2023; Liu et al. 2023b; Gu et al. 2025). A free-form quality score from such a judge is unstable and drifts between calls, so we do not ask for one: following the form-filling paradigm of G-Eval (Liu et al. 2023b), each axis is decomposed into explicit per-grade anchors, and the judge selects the grade whose description the output meets rather than inventing a number. Anchoring the grades this way is

what makes the scores reproducible and auditable—the reason a given output received a given grade is written into the rubric—and it constrains the judge to the clinical criteria below rather than its own priors. To remove the self-preference bias documented for LLM judges (Liu et al. 2023b; Zheng et al. 2023), the judge (DeepSeek-V4-Pro) is drawn from outside the baseline set, so no system grades its own outputs.

Where the anchors come from. The main-text Metrics section defines the four reasoning axes and the T1–T3 and P1–P3 tasks; here we give the per-grade anchors the judge applies. The four window axes (4/3/2/1) are shared across the specialist reading and the window probes, specialized by task (T1 diagnosis, T2 differential, T3 counterfactual) and by setting (both: nine-class; PPG: seven-class merged). The patient anchors (4/3/2/1/0, adding a 0 for output that does not address the axis at all) are written against clinical practice: acuity follows ACLS triage, and management and significance follow the 2023 AF / CHA₂DS₂-VASc (Joglar et al. 2024), 2017 VT/SCD (Al-Khatib et al. 2018), and 2018 bradycardia (Kusumoto et al. 2019) guidelines, so that a top-grade answer is not merely fluent but concordant with the standard of care. The PPG anchors additionally encode signal physics—a claim to a P wave or a QRS width from the pulse is capped as unfaithful, because the pulse cannot show either.

Window-level rubric: four axes, shared by the specialist and T1/T2/T3 (grades 4→1)

T1 — diagnosis

Correctness

- 4 exact GT class.
- 3 adjacent family (PVC↔PAC, VT↔MAT, AVB↔SND).
- 2 right normal-vs-arrhythmia side, wrong family.
- 1 wrong side or nonsensical.

Faithfulness

- 4 engages measured rate / R–R CV / QRS, numbers within ~25%, no hallucination.
- 3 one minor issue (a number off 25–50%).
- 2 a qualitative claim contradicts evidence, or a number off >50%.
- 1 fabricated or 2×-wrong values.

Validity

- 4 cited features correct and sufficient.
- 3 valid but incomplete.
- 2 some features irrelevant or contradictory.
- 1 features do not support the conclusion.

Completeness

- 4 all discriminators addressed (rate, R–R regularity, QRS/ventricular morphology, pauses).
- 3 one missing.
- 2 two or more missing.
- 1 conclusory, no feature support.

T2 — differential

Correctness

- 4 correct primary and the key alternatives, incl. the most dangerous.
- 3 correct primary, differential incomplete.

- 2 primary a near-miss, or a weak set.
- 1 wrong primary or irrelevant differentials.

Faithfulness

- 4 every exclusion cites a real evidence-based reason.
- 3 one exclusion weakly grounded.
- 2 an exclusion contradicts evidence or is empty.
- 1 exclusions largely fabricated.

Validity

- 4 each alternative ruled out by a correct discriminator.
- 3 one weak discriminator.
- 2 some exclusions illogical.
- 1 logic broken or circular.

Completeness

- 4 covers the relevant alternatives and explicitly the most dangerous.
- 3 one notable alternative missing.
- 2 omits the most dangerous.
- 1 no real differential.

T3 — counterfactual

Correctness

- 4 correct pivotal change and resulting diagnosis.
- 3 correct change but adjacent result (or vice versa).
- 2 a plausible change, but not the gold one.
- 1 wrong or irrelevant.

Faithfulness

- 4 the change is evidence-relevant and realistic.
- 3 minor imprecision.
- 2 contradicts current evidence, or vague.
- 1 fabricated or irrelevant.

Validity

- 4 causal link textbook-correct.
- 3 plausible but loose.
- 2 weak or partial link.
- 1 change does not lead to the result.

Completeness

- 4 both the change and the result explicit.
- 3 one slightly vague.
- 2 only one given.
- 1 no real counterfactual.

PPG (seven-class) specialization. Under PPG, correctness is satisfied by naming a “premature/ectopic beat” (→PVC/PAC) or a “broad/fast run” (→VT/MAT) and keeping an ECG-only-distinguishable pair merged, deferring to the ECG; excluding an indistinguishable alternative is penalized. Faithfulness caps at 2 for any claim to see P-waves or QRS, which the pulse cannot show, and completeness rewards flagging that an ECG is required where atrial or ventricular detail is needed.

Patient-level rubric: P1/P2/P3, two axes each (grades 4→0; the two axes sum to 8)

P1 — impression

Correctness — were the arrhythmias present named?

- 4 named every arrhythmia present and asserted none that were absent.
- 3 named the present arrhythmias but missed one, or added an absent one as a secondary mention.

- 2 named a present arrhythmia only as a secondary mention, and asserted an absent one as dominant.
- 1 mentioned a present arrhythmia at least once.
- 0 named none of the arrhythmias present.

Faithfulness — do the cited features fit the named arrhythmia?

- 4 cited only features characteristic of the named arrhythmia (e.g. AF → irregularly irregular R-R).
- 3 cited features of both the present and an absent arrhythmia.
- 2 cited a feature of an absent arrhythmia, though the feature itself is a valid basis for that arrhythmia (right features, wrong diagnosis).
- 1 cited a feature that supports no arrhythmia, or one physically impossible (a P-wave/PR interval on PPG, R-R CV 0.06).
- 0 cited no features (diagnosis named only).

P2 — clinical significance

Prioritization — is the most dangerous finding foregrounded? (severity, not frequency)

- 4 foregrounded the most dangerous finding as the one to act on, with severity stated correctly.
- 3 correct priority, but severity over- or under-stated.
- 2 mentioned the dangerous finding but did not prioritize it, or prioritized a less dangerous one.
- 1 did not mention the most dangerous finding; presented only a less dangerous one as the one to act on.
- 0 prioritized no finding at all.

Context — is the risk tied to why this patient is at risk? (e.g. 70 yr, hypertension → AF → stroke)

- 4 “this patient is at risk *because...*” holds—patient features connect to a specific risk.
- 3 the link holds but omits one relevant feature (e.g. names AF risk without the comorbidity driving it).
- 2 patient data restated but not tied to risk—true of any patient.
- 1 used a feature unrelated to the arrhythmia as the risk basis, or judged risk backwards.
- 0 used no patient information.

P3 — management

Management — is the plan specific to this rhythm and this patient?

- 4 specific to the rhythm and patient, guideline-concordant, with a correct citation.
- 3 appropriate but generic, or the citation is loose.
- 2 a generic plan that would read the same for any rhythm (no rhythm-specific decision).
- 1 an inappropriate plan (anticoagulation without documented AF, discharging a sustained VT).
- 0 gave no management plan.

Urgency & safety — is the urgency right and free of over-/under-treatment?

- 4 urgency correct and management connects to real findings safely.
- 3 urgency slightly off but still safe.
- 2 urgency clearly wrong, or a generic “urgent” label not tied to a specific finding.
- 1 harmful (discharges a VT, wrong anticoagulation, a life-threatening rhythm missed).
- 0 did not address urgency at all.

Baseline Prompts

Every baseline is given the same information the task provides for answering—the raw ECG and PPG, their plots, the taxonomy, and, at the patient level, the metadata and every window—and is asked for identical outputs. What our system adds is not more information but its orchestration over that shared input, so the comparison isolates the agent construction. At the window level a baseline receives both the ECG and the PPG—their raw sample values and their plots—together with the taxonomy, and is asked for the class, the span, and T1/T2/T3 in one call. The difference from our agent is the tool layer it does not receive: the extracted measurements, R-peak coordinates, morphological hints, and the clinical-reasoning scaffold. All baseline calls use temperature 0 and request a single JSON object.

Baseline prompt — window level (both signals, nine-class)

System. “You are an expert cardiologist. Analyze the signal(s) and image(s) and answer ALL items.”

User message. A multimodal message interleaving two images and text:

[ECG plot image] [PPG plot image]

You are given two time-aligned signals of the same cardiac window. The first image is the ECG, the second is the PPG. ECG signal values (index 0...N-1): 87, 42, 11, 17, 3, 9, ... (integer, $\times 100$). PPG signal values: 88, 22, 14, 28, 8, 11, ... Then the nine-class arrhythmia reference, each class with its criteria (AF: irregularly irregular R-R, CV>0.15, absent P waves, narrow QRS; VT: fast run of wide QRS; and so on for SVT, PVC, PAC, SND, MAT, AVB, Normal).

Requested output. A single JSON object { "class", "span", "T1", "T2", "T3" }, where:

- class: the single best-fitting class from the taxonomy; span: the sample range [start, end] of the abnormal segment.
- T1: diagnose the rhythm and the key features (rate, R-R regularity, P-wave, QRS morphology) that support it.
- T2: the plausible differentials and why each is excluded, the most dangerous first.
- T3: the single change that would most alter the diagnosis, and to what.

Decoding. temperature=0, max_tokens=900, response_format=json_object.

Unlike our specialist, the baseline is given no extracted features, R-peak coordinates, morphological hints, or reasoning scaffold; it reads the class directly from the raw signal and plots.

Baseline prompt — patient level (P1/P2/P3)

System. “You are a cardiologist reviewing a patient’s peri-operative continuous cardiac monitoring. For the same patient you are given several time windows; for each window you receive the ECG and PPG raw waveforms (100 Hz, ~ 10 s) and their plots, plus clinical metadata. Produce an integrated

patient-level assessment. Ground every statement in the actual waveform/plot findings and the metadata—do not invent findings that are not supported by the data.”

User message. The patient’s metadata followed by every window in temporal order, laid out as:

PATIENT METADATA: Age=78; Sex=female; BMI=25.0; Hypertension=yes; Diabetes=no. This patient has 5 monitoring windows, in temporal order below.

=== Window 1 === ECG raw (100 Hz): [...];
PPG raw (100 Hz): [...]; [ECG plot image]
[PPG plot image]
=== Window 2 === ... (through window 5)

Every window thus supplies four inputs—its ECG and PPG raw waveforms and its two plot images—and all of the patient’s windows are included.

Requested output. Compact JSON { "P1", "P2", "P3" }, grounded only in the provided windows and metadata:

- P1 (integrated diagnosis): across all windows, which arrhythmia(s) appear in which window(s), and the dominant / most representative rhythm, citing specific waveform features.
- P2 (contextual significance): given the metadata (age, surgery, comorbidities), the clinical significance and the single finding most important to act on, and why.
- P3 (urgency & management): the appropriate response (further monitoring vs. intervention, and urgency level) and the specific evidence supporting it.

Table 9: Per-class F1. Setting both is nine-class (ECG specialist); PPG is seven-class merged (PPG specialist).

class	spec.	CARDIO.	gpt-4o	gem.	qwen	medg.	llama
<i>Setting: both (nine-class)</i>							
normal	0.93	0.95	0.84	0.42	0.77	0.75	0.00
atrial fibrill.	0.68	0.83	0.65	0.42	0.05	0.00	0.18
MAT	0.42	0.43	0.01	0.00	0.00	0.00	0.00
SVT	0.11	0.08	0.04	0.04	0.00	0.00	0.00
VT	0.00	0.00	0.25	0.25	0.00	0.00	0.00
PVC	0.50	0.61	0.21	0.48	0.43	0.00	0.00
PAC	0.17	0.45	0.30	0.04	0.20	0.00	0.00
AVB	0.18	0.12	0.00	0.00	0.00	0.00	0.00
SND	0.55	0.57	0.00	0.02	0.00	0.00	0.00
macro-F1	0.39	0.45	0.26	0.19	0.16	0.08	0.02
<i>Setting: PPG (seven-class merged)</i>							
normal	0.94	0.93	0.72	0.70	0.67	0.66	0.70
atrial fibrill.	0.76	0.75	0.30	0.48	0.38	0.18	0.26
VT/MAT	0.22	0.38	0.00	0.00	0.00	0.00	0.04
SVT	0.00	0.12	0.00	0.00	0.00	0.03	0.06
PVC/PAC	0.53	0.51	0.01	0.18	0.33	0.05	0.13
AVB	0.00	0.04	0.00	0.00	0.00	0.00	0.00
SND	0.13	0.24	0.00	0.00	0.00	0.00	0.00
macro-F1	0.37	0.43	0.15	0.19	0.20	0.13	0.17

Experiments Details

Per-class classification. Table 9 breaks the window-level macro-F1 out by class, and it locates where our margin comes

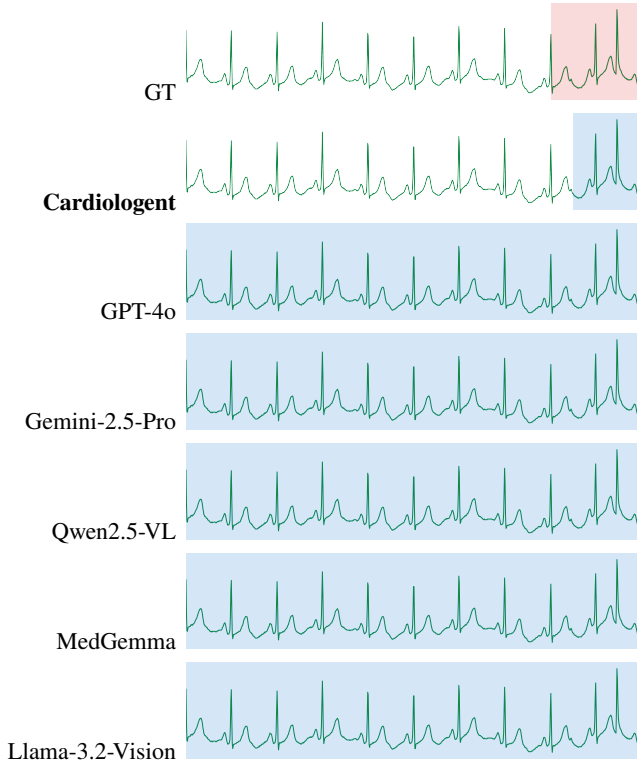


Figure 6: Grounding on ECG window 4925_w0099 (ground-truth PAC). Each row is one model’s localized abnormal span on the same trace; the top row is the ground truth. CARDIOLOGENT localizes a tight span on the PAC, whereas the baselines mark almost the entire window, inflating IoU while collapsing MCC toward zero.

from. It is not a uniform edge on the common rhythms: it is that the general baselines *collapse to zero* on the rare arrhythmias. On *both* signals, every baseline scores 0.00–0.01 on MAT, SVT, AVB, and SND—four of the nine classes they effectively cannot detect at all—while our system holds usable F1 on each (0.43, 0.08, 0.12, 0.57). This is what the 1.8× macro-F1 gap (0.449 vs. 0.255) is made of: not a higher ceiling on normal and AF, where the stronger baselines are already competent, but coverage of the rare classes where they read nothing. Within our own stack, the agent layer helps most on the classes the specialist reads least well—PAC more than doubles (0.17→0.45) once the atrial-activity and premature-beat tools weigh in, and AF (0.68→0.83) and PVC (0.50→0.61) improve likewise—though it is not uniformly additive: on SVT and AVB the extra deliberation slightly lowers an already low score. The help is largest exactly where the raw signal is thinnest: under *PPG*, where the specialist is weak on almost every arrhythmia, the agent lifts nearly all of them—SVT and AVB from 0.00 to 0.12 and 0.04, SND 0.13→0.24, and the merged VT/MAT 0.22→0.38—so the tools and debate matter more, not less, when the pulse alone carries the reading. The one class we do not recover is VT: with 45 test windows on *both* it is too rare for the specialist to learn, and here two baselines edge us (0.25 vs. 0.00), though

under *PPG* the merged VT/MAT class reaches 0.38. Figure 7 walks through a window where the ECG and PPG agents first disagree and window fusion settles them on the correct class.

Grounding. Figure 6 shows where each model localizes the abnormality on ECG window 4925_w0099, whose ground truth is a premature atrial contraction (PAC). A PAC is a focal event confined to a few beats, so a correct localization is a narrow span. CARDIOLOGENT places a tight span on the PAC, while the general baselines mark essentially the whole window as abnormal—failing to localize. Marking everything inflates the overlap-based IoU yet carries no localizing information, which is exactly what the per-sample MCC exposes: a whole-window prediction overlaps the truth by chance (high IoU) but falls to $MCC \approx 0$. The split shows in the window-level grounding results: Llama posts a competitive IoU (0.533, second only to ours on *both*) but an MCC of 0.001—overlap without localization—and across all baselines none exceeds MCC 0.18, with MedGemma going negative, even though several reach moderate IoU. CARDIOLOGENT leads on both (0.584 / 0.305). We therefore report both. Let P be the set of samples a model predicts abnormal and G the true abnormal samples:

$$IoU = \frac{|P \cap G|}{|P \cup G|},$$

$$MCC = \frac{TP \cdot TN - FP \cdot FN}{\sqrt{(TP+FP)(TP+FN)(TN+FP)(TN+FN)}}.$$

IoU rewards overlap, while MCC scores per-sample agreement and cannot be inflated by a blanket abnormal prediction (which drives TN to zero).

Reasoning by axis. Table 10 gives every rubric axis behind the reasoning scores, with the per-probe SUM (window SUM over four axes, max 16; patient SUM over two axes, max 8). Read by axis rather than by SUM, it shows *where* the margin lives and where it does not.

Window level. Our advantage is a reading whose four axes hold *together*: the cited features are real (faithful), they actually support the diagnosis (valid), and the discriminators are covered (complete), not just a correct label. This shows in how the axes move as a block rather than one at a time. On T1 our SUM leads (12.51 vs. GPT-4o’s 11.57 on *both*), and the widest single-axis gaps are in validity and completeness (3.22 vs. Gemini’s 1.99; 3.03 vs. 2.84)—the axes that ask whether the evidence is wired to the conclusion, which is what the tools and the debate enforce. A baseline that lands the right label often does not carry the rest of the reading with it: Gemini’s T1 correctness is 2.24 but its faithfulness falls to 1.59, an answer that is right without being grounded. T2—the differential—is the hardest probe for everyone (correctness the lowest cell for every model, our 2.34 vs. the baselines’ 1.02–2.14), since ruling alternatives out demands more than naming the rhythm, and it is where completeness most separates us (2.75 vs. Gemini’s 2.44). T3 (counterfactual) is where the stronger baselines come closest, as reasoning about a hypothetical change draws on general-purpose ability rather than on reading the trace (12.02 vs. our 12.88 on *both*). The few cells a baseline tops are single axes

Table 10: Reasoning by rubric axis (LLM-judge, 1–4). Window probes have four axes (SUM/16); patient tasks two axes (SUM/8). Best per column in bold.

<i>Window level — Setting both</i>															
model	T1 (diagnosis)					T2 (differential)					T3 (counterfactual)				
	corr	faith	valid	compl	SUM	corr	faith	valid	compl	SUM	corr	faith	valid	compl	SUM
CARDIOLOGENT	3.54	2.72	3.22	3.03	12.51	2.34	2.60	2.81	2.75	10.50	2.44	3.47	3.21	3.76	12.88
gpt-4o	3.17	2.55	2.92	2.94	11.57	1.65	2.40	2.69	2.25	8.99	2.34	3.37	2.79	3.53	12.02
gemini	2.24	1.59	1.99	2.84	8.66	2.14	2.52	2.84	2.44	9.94	2.02	2.98	2.83	3.48	11.30
qwen	2.86	2.32	2.73	2.83	10.74	1.13	1.32	1.54	1.55	5.55	1.41	2.05	2.15	2.61	8.21
medgemma	2.83	1.71	2.46	2.75	9.75	1.02	1.13	1.33	1.39	4.87	1.24	1.95	2.12	1.90	7.21
llama	1.12	1.05	1.09	1.20	4.47	1.06	1.22	1.30	1.20	4.78	1.30	2.26	1.66	2.36	7.58
<i>Window level — Setting PPG</i>															
model	T1 (diagnosis)					T2 (differential)					T3 (counterfactual)				
	corr	faith	valid	compl	SUM	corr	faith	valid	compl	SUM	corr	faith	valid	compl	SUM
CARDIOLOGENT	3.32	3.48	3.35	2.74	12.90	1.96	2.23	2.11	1.95	8.26	2.32	3.36	2.80	3.70	12.19
gpt-4o	2.10	2.13	2.07	1.90	8.21	1.16	2.47	2.07	1.31	7.01	2.12	3.03	2.72	3.43	11.30
gemini	2.82	2.32	2.59	2.77	10.50	1.30	1.92	1.93	1.41	6.56	1.38	2.37	1.89	3.64	9.29
qwen	2.77	2.71	2.71	2.85	11.03	1.09	2.03	2.05	1.48	6.65	1.17	2.28	1.71	2.47	7.63
medgemma	2.75	2.48	2.61	2.77	10.60	1.06	1.76	1.71	1.24	5.78	1.17	2.27	1.74	2.04	7.22
llama	1.70	1.35	1.41	1.33	5.80	1.58	2.05	1.92	1.53	7.07	1.56	2.37	1.78	2.67	8.38
<i>Patient level — judge (65 patients)</i>															
model	P1 (diagnosis)			P2 (significance)			P3 (management)								
	correct	faithful	SUM	priority	context	SUM	mgmt	urgency	SUM						
CARDIOLOGENT	1.97	2.69	4.66	2.34	2.03	4.37	2.14	2.42	4.55						
gemini	0.89	2.31	3.20	1.69	1.62	3.31	1.85	2.09	3.94						
gpt-4o	0.55	2.23	2.78	1.48	1.35	2.83	1.37	1.57	2.94						
qwen	0.12	2.00	2.12	1.06	1.17	2.23	1.29	1.57	2.86						
medgemma	0.00	1.91	1.91	0.85	1.25	2.09	1.32	1.57	2.89						
llama	0.12	1.94	2.06	0.88	0.62	1.49	1.38	1.52	2.91						
<i>Patient level — cardiologists (20 patients)</i>															
model	P1 (diagnosis)			P2 (significance)			P3 (management)								
	correct	faithful	SUM	priority	context	SUM	mgmt	urgency	SUM						
CARDIOLOGENT	2.23	2.35	4.58	2.35	2.35	4.70	2.52	2.50	5.03						
gemini	1.62	1.70	3.33	1.82	1.93	3.75	2.08	2.15	4.22						
gpt-4o	1.40	1.57	2.98	1.65	1.82	3.48	1.88	1.98	3.85						
qwen	0.68	0.80	1.48	1.18	1.30	2.48	1.43	1.50	2.92						
medgemma	0.55	0.82	1.38	1.30	1.32	2.62	1.27	1.52	2.80						
llama	0.40	0.60	1.00	0.88	1.00	1.88	1.05	1.20	2.25						

of a single probe—Gemini’s T2 validity on *both* (2.84 vs. 2.81), GPT-4o’s T2 faithfulness on *PPG*—and none carries the probe it sits in. The consistent read is that our reasoning is not merely more often right but more often *coherent*: the diagnosis, its evidence, and its differential agree, which is exactly the property a downstream patient-level decision has to rest on.

Patient level. Here the axis split is starkest, and it exposes a specific failure. Faithfulness—citing real, on-signal features—is something the baselines do nearly as well as we do: on P1 their faithfulness sits at 2.23–2.31 against our 2.69. But the diagnosis those features are attached to is wrong. P1 correctness separates the systems completely: 1.97 for us against 0.89, 0.55, and 0.00 for Gemini, GPT-4o, and MedGemma. The general models are fluent about the signal and confidently mistaken about the patient—a SUM would

let their faithfulness mask the failed diagnosis, and only the per-axis view makes it visible. The same holds through P2 and P3, where using the metadata to place the risk (context, our 2.03 vs. Llama’s 0.62) and choosing a rhythm-specific management (our 2.14 vs. 1.29–1.85) are what the baselines miss; CARDIOLOGENT leads on every patient axis with no baseline topping any. At the patient level, then, the baselines stop at a plausible-sounding account while the judgment that a clinician would act on—the right diagnosis, tied to the right patient context, with rhythm-specific management—is reached only by our system; on the axis that matters most for a bedside decision, a confidently written wrong answer is the baselines’ typical failure, and the one our pipeline is built to avoid.

Table 11: **Patient-level diagnosis accuracy** ($n=65$): exact class, same management family, and any mention. $\pm$ is a 95% bootstrap CI ($B=10,000$).

	exact	family	mention
CARDIOLOGENT	0.48 $\pm$ 0.12	0.60 $\pm$ 0.12	0.77 $\pm$ 0.10
GPT-4o	0.17 $\pm$ 0.09	0.23 $\pm$ 0.10	0.26 $\pm$ 0.11
Gemini-2.5-Pro	0.12 $\pm$ 0.08	0.15 $\pm$ 0.08	0.34 $\pm$ 0.12
Qwen2.5-VL-72B	0.08 $\pm$ 0.07	0.12 $\pm$ 0.08	0.25 $\pm$ 0.11
MedGemma-27B	0.00 $\pm$ 0.00	0.00 $\pm$ 0.00	0.22 $\pm$ 0.10
Llama-3.2-90B-Vision	0.00 $\pm$ 0.00	0.00 $\pm$ 0.00	0.00 $\pm$ 0.00

Table 12: Patient-level score by the patient’s dominant class: SUM over P1–P3, max 24. The *all* row is the patient-level assessment summed over its three tasks. Best per row in bold, second best underlined.

class	n	CARDIOLOGENT	gem.	gpt-4o	qwen	medg.	llama
<i>LLM judge</i>							
AF	13	21.4	<u>14.2</u>	6.5	6.8	6.7	7.2
SVT	13	6.8	7.9	6.2	6.9	<u>7.0</u>	6.2
VT	9	8.3	13.7	<u>11.7</u>	6.1	5.4	6.0
PVC	7	21.3	13.3	<u>14.3</u>	10.0	7.7	5.4
PAC	8	12.1	8.9	<u>10.4</u>	7.4	7.4	7.1
MAT	5	10.6	8.0	7.0	7.6	6.8	7.0
SND	7	15.9	6.7	7.0	7.1	<u>7.7</u>	6.7
AVB	3	10.7	6.0	6.0	6.3	<u>6.7</u>	4.7
all	65	13.6	<u>10.4</u>	8.6	7.2	<u>6.9</u>	6.5
<i>Cardiologists (20-patient subset)</i>							
AF	4	20.4	<u>19.4</u>	6.1	6.1	6.5	4.8
SVT	4	13.5	<u>10.0</u>	5.8	4.8	7.9	4.9
VT	3	6.0	<u>14.5</u>	16.3	5.7	7.3	6.0
PVC	2	<u>22.0</u>	12.0	22.5	14.0	7.5	5.2
PAC	2	<u>12.8</u>	10.0	17.2	11.8	4.8	5.0
MAT	2	10.2	4.0	<u>6.0</u>	5.8	4.2	4.8
SND	2	12.8	3.5	3.8	5.0	<u>9.5</u>	4.0
AVB	1	17.0	6.0	<u>10.5</u>	4.0	4.5	8.0
all	20	14.3	<u>11.3</u>	10.3	6.9	6.8	5.1

Patient-level diagnosis accuracy. Table 11 scores the dominant rhythm at three grades: *exact*, the ground-truth dominant rhythm named (or a single-window dangerous finding correctly flagged as suspected); *family*, the named rhythm sharing a management tier with the truth; and *mention*, the true arrhythmia named anywhere in P1. The families group the rhythms that share a management tier: {VT, MAT}, {PVC, PAC}, and {SND, AVB}, with AF, SVT, and normal standing alone, following the ventricular, supraventricular, and bradycardia management guidelines (Al-Khatib et al. 2018; Page et al. 2016; Kusumoto et al. 2019). Unlike the rubric scores, this is a hard label match, so it checks the P1 result against something a graded judgement cannot inflate.

CARDIOLOGENT names the exact dominant rhythm for 0.48 of patients against 0.17 for the strongest baseline, and the ordering is the same at every grade. More telling is how the three grades move. Relaxing *exact* to *family* gains us 0.12 (0.48 $\rightarrow$ 0.60), so a good part of our residual error stays inside the management tier and would not change what is done; for

Table 13: Rater agreement per task, ICC(2,1) with absolute agreement ($\pm$ its 95% CI), on the 20-patient subset. R1, R2 are the two cardiologists.

	R1 $\leftrightarrow$ R2	Judge $\leftrightarrow$ R1	Judge $\leftrightarrow$ R2
P1 (diagnosis)	0.711 $\pm$ 0.105	0.651 $\pm$ 0.117	0.669 $\pm$ 0.140
P2 (significance)	0.612 $\pm$ 0.146	0.714 $\pm$ 0.127	0.549 $\pm$ 0.153
P3 (management)	0.480 $\pm$ 0.170	0.518 $\pm$ 0.172	0.488 $\pm$ 0.166

Case study: Window Fusion

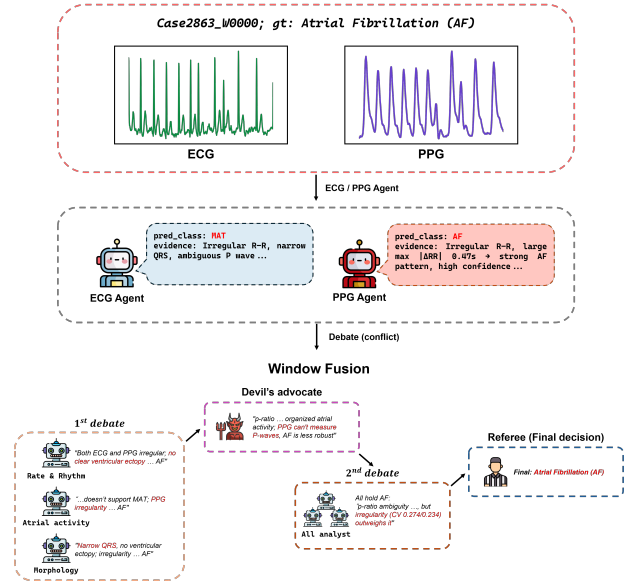


Figure 7: **Case study: Window fusion** (Case2863_W0000; GT: atrial fibrillation). The ECG agent misreads an ambiguous P wave and calls MAT; the PPG agent calls AF from the R–R irregularity. Because the two conflict, fusion runs a two-round analyst–devil’s-advocate debate: the devil’s advocate notes the PPG cannot confirm P waves, the analysts weigh this against the measured irregularity (R–R CV 0.274/0.234), and the referee settles on the correct class, AF. The right answer is reached by weighing evidence, not by deferring to the more confident agent.

the baselines the same relaxation gains almost nothing (GPT-4o 0.17 $\rightarrow$ 0.23, Gemini 0.12 $\rightarrow$ 0.15), meaning that when they are wrong they are wrong outside the tier. Relaxing again to *mention* separates two failure modes: Gemini (0.15 $\rightarrow$ 0.34) and MedGemma (0.00 $\rightarrow$ 0.22) do often name the true rhythm somewhere in P1 but not as the dominant one—the finding is present and mis-ranked—whereas Llama never names it at all. Our own gap between *exact* and *mention* (0.48 vs. 0.77) is the same failure in milder form: for roughly a third of patients the true arrhythmia appears in the impression without being made dominant.

Patient-level score by class. Table 12 breaks the patient-level total (P1–P3, max 24) out by each patient’s dominant rhythm. CARDIOLOGENT leads six of the eight classes, by the widest margins where the window level reads the rhythm

Case study: Patient Level

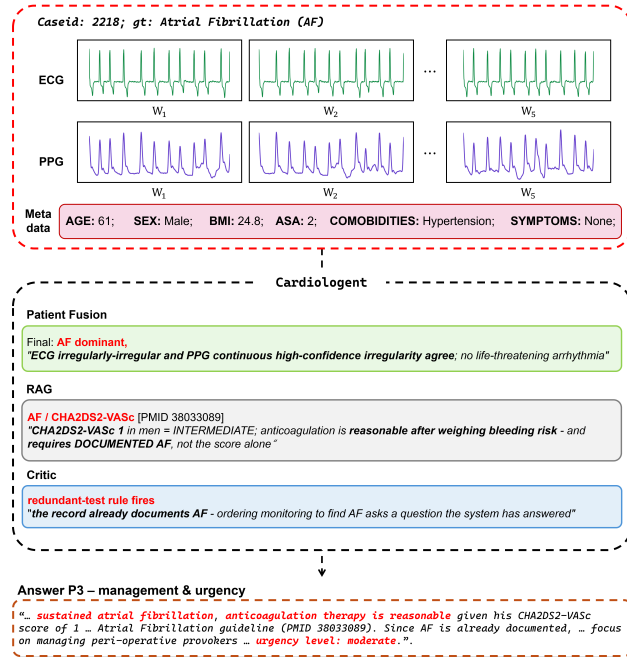


Figure 8: **Case study: Patient-level** (case 2218; GT: atrial fibrillation). The ECG and PPG readings agree, so patient fusion returns AF as the dominant rhythm without a conflict to settle. Retrieval supplies the guideline that sets how strongly the recommendation may be put—CHA₂DS₂-VASc 1 in a man is intermediate, and anticoagulation is licensed by documented AF rather than by the score—and the critic removes the draft’s request for monitoring to find an AF the record already documents. Excerpts are verbatim; ellipses mark omitted text.

well: AF 21.4 against 14.2, PVC 21.3 against 14.3, SND 15.9 against 7.7. The two it does not lead are VT (8.3 against Gemini’s 13.7) and SVT (6.8 against 7.9)—the two rhythms the window level detects least reliably (0.00 and 0.08 F1 on both). A patient-level decision inherits the window readings it is assembled from, and where those fail the decision does too. Per-class counts are small (3–13 patients), so the ordering carries more than the individual values. The cardiologist block, on the 20-patient subset, follows the same pattern—CARDIOLOGENT first overall and on most classes, VT the clearest loss under both raters—with the counts too small (1–4) to read individual classes closely.

Rater agreement by task. Table 13 breaks the overall judge–cardiologist ICC(2,1) out by task. Agreement falls from P1 to P3 for the judge and cardiologists alike—the integrated diagnosis is scored most consistently, and management least, where even the two cardiologists agree only moderately (0.48). Management is where clinical judgement legitimately diverges—the same rhythm admits more than one defensible plan—so the lower agreement there reflects the task rather than the judge: where the judge is noisier is exactly where the two cardiologists disagree more as well.

CARDIOLOGENT is still scored highest on P3 by both raters, so its lead holds even on the task clinicians agree on least. At every task the judge–cardiologist agreement overlaps the human–human agreement within its confidence interval.

Case study: window fusion. Figure 7 traces a single window (GT: atrial fibrillation) through the fusion, and shows why the two-signal read is more than a vote. The ECG agent calls MAT—it sees irregular R–R and a narrow QRS but reads the P wave as present—while the PPG agent, which cannot be misled by an ambiguous P wave, calls AF from the irregularity alone (max $|\Delta RR|$ 0.47 s). The two conflict, so fusion opens a debate rather than averaging. In the first round all three analyst lenses (rate/rhythm, atrial activity, morphology) favor AF: no ventricular ectopy, narrow QRS, and marked irregularity. The devil’s advocate raises the one real objection—an apparent P-ratio could mean organized atrial activity, and the PPG cannot see P waves, so its AF is the less anatomically grounded of the two claims. In the second round the analysts weigh that against the measured irregularity (R–R CV 0.274/0.234) and hold AF, and the referee returns AF. The correct label here is not the more confident agent’s by default; it survives because the objection was voiced and then outweighed by evidence, which is the property the per-sample grounding and the per-axis reasoning scores reward throughout this section.

Case study: patient level. Figure 8 traces one record (case 2218; GT: atrial fibrillation) through the three patient-level stages. *Patient fusion* has no conflict to settle: the ECG is irregularly irregular and the PPG reports continuous irregularity with high confidence, so the two agree and the profile comes back as AF dominant with no life-threatening rhythm—the work falls to the stages after it. *Retrieval* then supplies not the diagnosis but the strength at which a recommendation may be put: for a man with a CHA₂DS₂-VASc of 1 the score is intermediate, so anticoagulation is *reasonable* rather than indicated, and it is licensed by documented AF rather than by the score alone. The *critic* deletes a step the draft had added—an ambulatory monitor to screen for occult AF—because the record already documents AF, and ordering a test to settle a question the system has answered only defers the decision that question was meant to inform. What reaches P3 is the therapy at that strength, without the confirmatory monitoring, at moderate urgency.